


Unison UK-C1
 Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

Donepezil (Agents for alzheimer-type dementia)

ドネペジル (アルツハイマー型認知症治療薬)

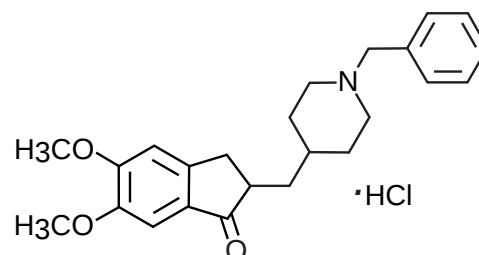

donepezil hydrochloride

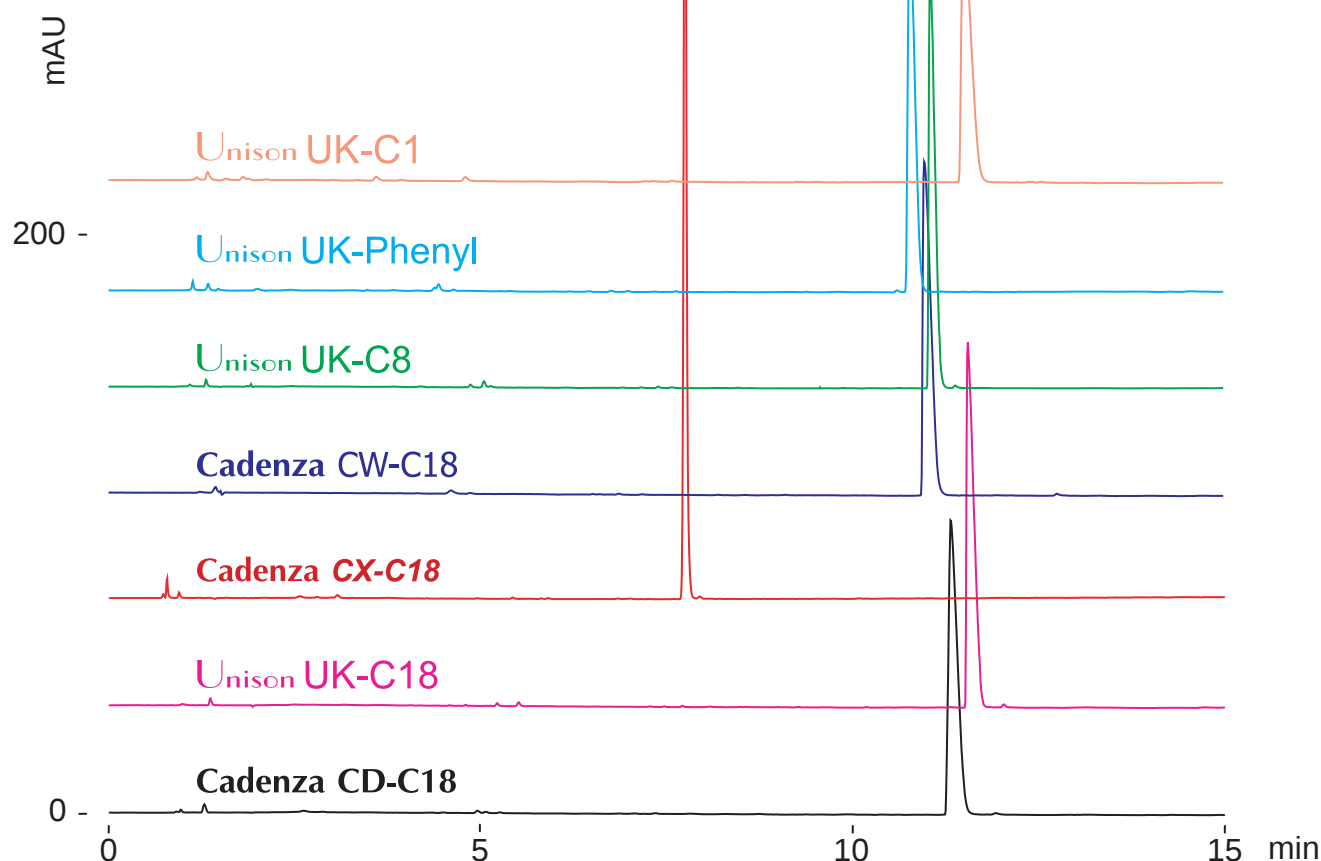

100 x 4.6 mm

A: water / acetonitrile / formic acid = 90 / 10 / 0.5

B: water / acetonitrile / formic acid = 10 / 90 / 0.5

0-50%B (0-15min)

1.0mL/min (11MPa), 37°C, 270nm

10μL (0.1mg/mL, water / formic acid = 100 / 1)