




Unison UK-C8

100 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Loratadine and desloratadine (Antiallergic agent)

ロラタジン及びデスロラタジン（抗アレルギー薬）のLC-MS分析

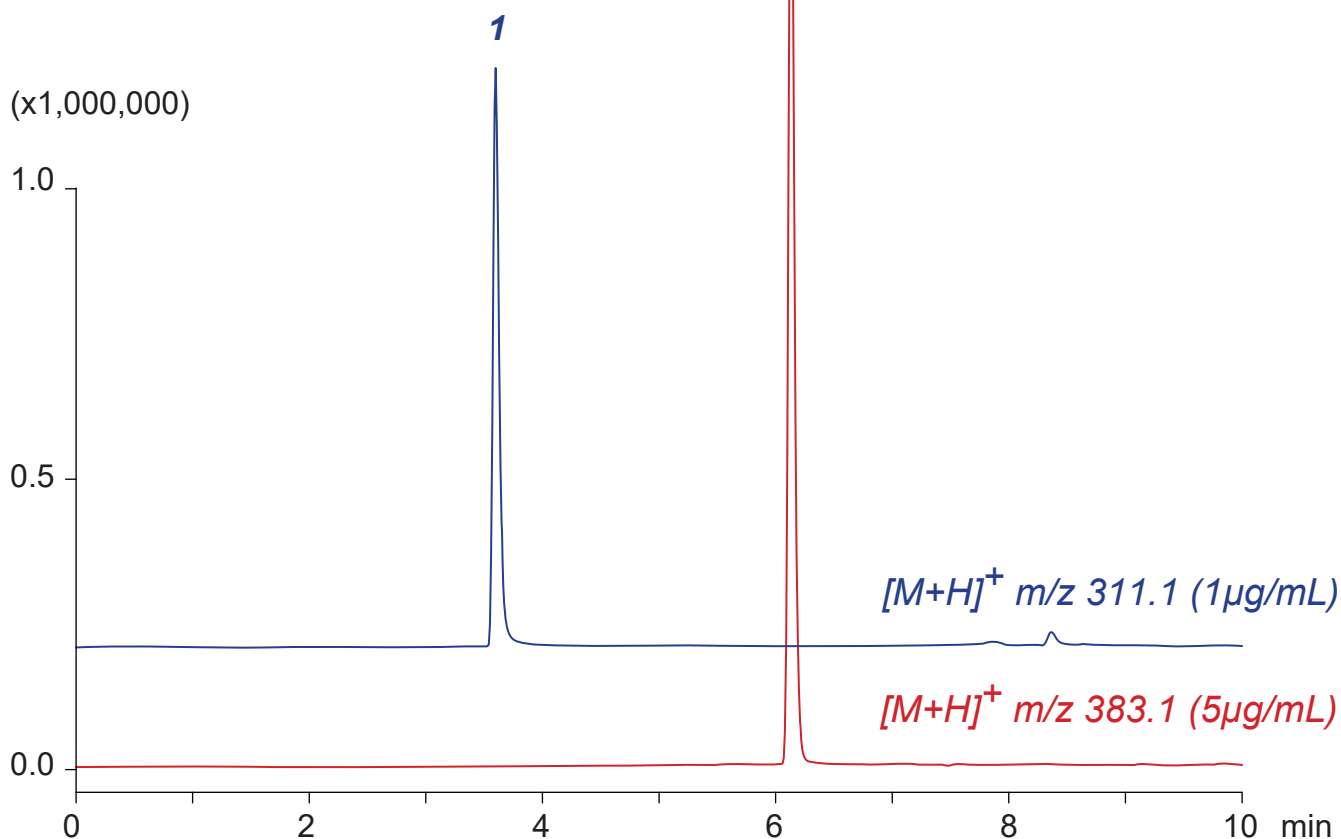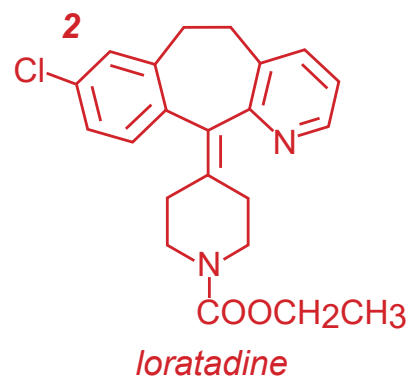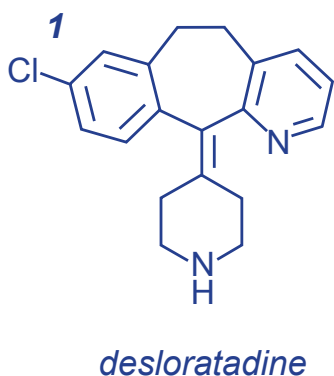

Unison UK-C8, 100 x 3 mm

A: acetonitrile / water / formic acid = 10 / 90 / 0.1

B: acetonitrile / water / formic acid = 90 / 10 / 0.1

0-100%B (0-10min)

0.4mL/min (8MPa), 37°C, Single Quad. MS (ESI)

5 μ L (1-5 μ g/mL, ACN / water / HCOOH = 25/ 75 / 1)