


Unison UK-C1
 Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Dimesitylammonium pentafluorobenzenesulfonate (Condensation organocatalysts)

ジメシチルアンモニウム ペンタフルオロベンゼンスルホナート
(縮合有機触媒)

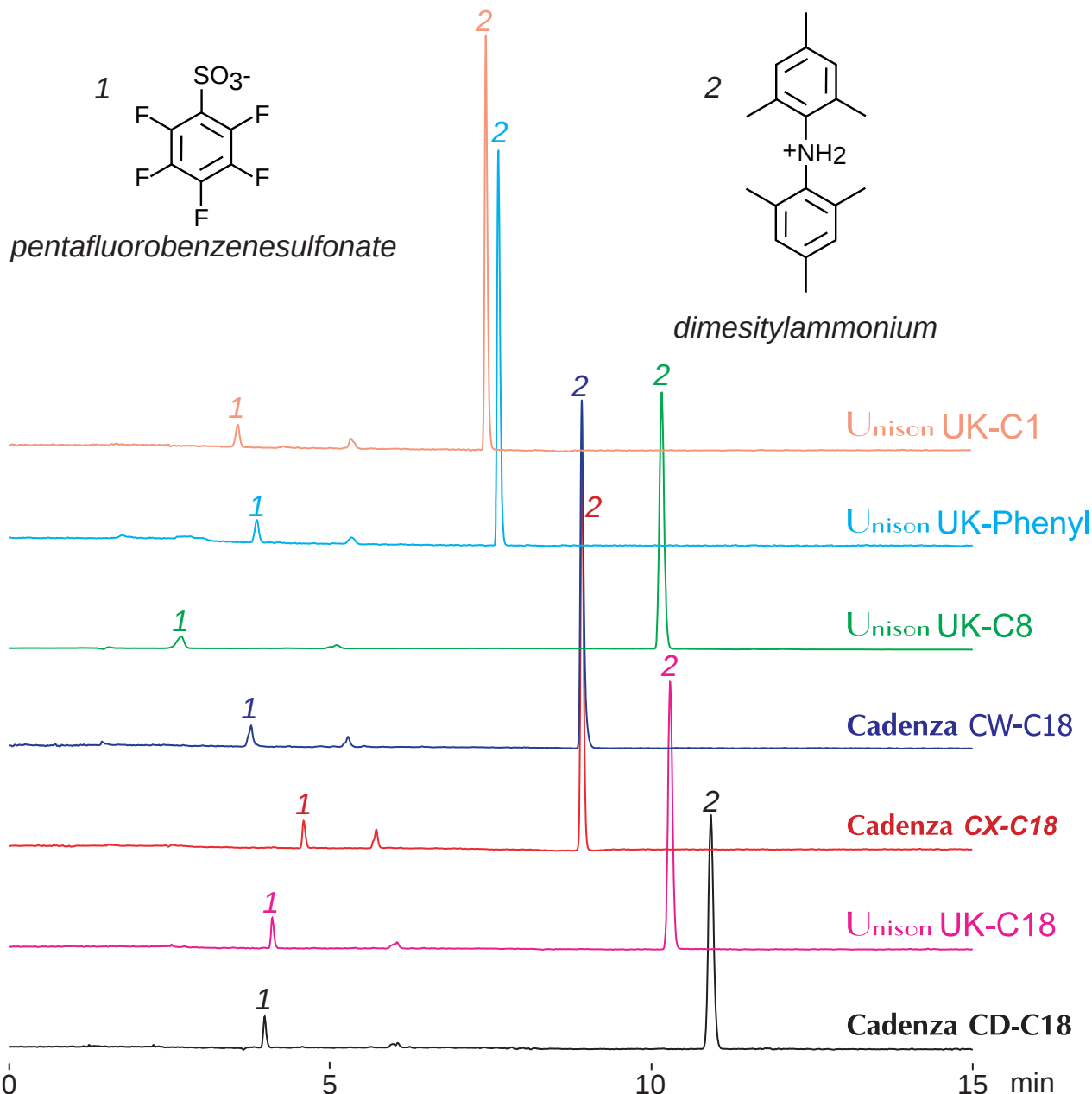

100 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 10 / 90

0-100%B (0-7min), 100%B (7-15min)

1.0mL/min (10MPa), 37°C, 260nm

10μL (0.1mg/mL, 100mM ammonium formate / acetonitrile = 30 / 70)