

Unison UK-Amino


Nardis ND-NX

150 x 4.6 mm

HPLC COLUMN APPLICATION

Entecavir (HBV(hepatitis B virus) therapeutic drug)

エンテカビル (HBV(B型肝炎ウイルス) 治療薬)

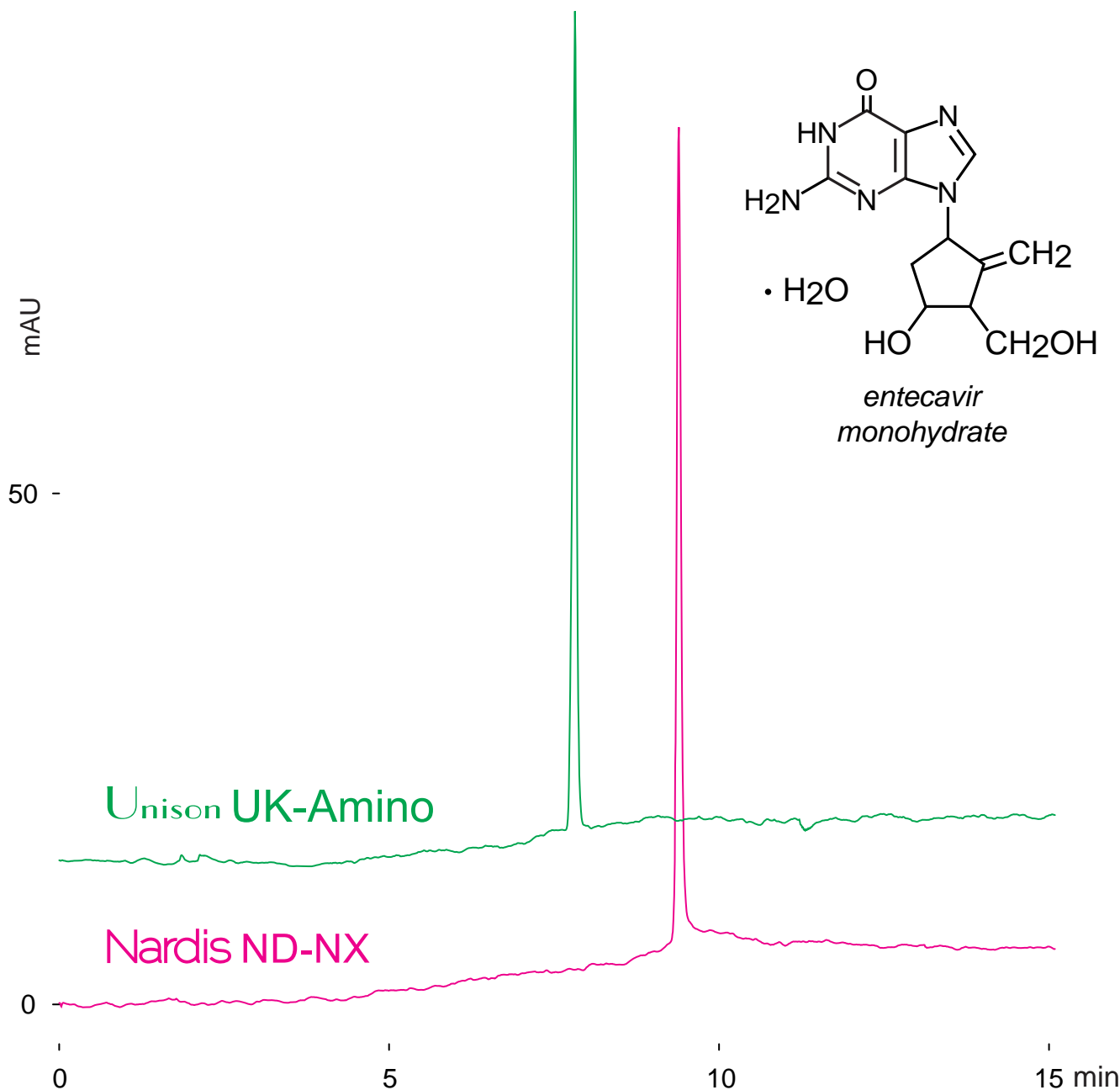

150 x 4.6 mm

A: acetonitrile

B: acetonitrile / 50mM ammonium formate = 50 / 50

0-100%B (0-15min)

1.0mL/min (6MPa), 37°C, 260nm

10µL (50µg/mL, 100mM HCOONH₄/EtOH=50/50)