


Unison UK-C18
 Cadenza CD-C18
 Unison UK-C8
 Unison UK-Phenyl

100 x 4.6 mm

HPLC COLUMN APPLICATION

beta-Carboline alkaloids

カルボリンアルカロイド類

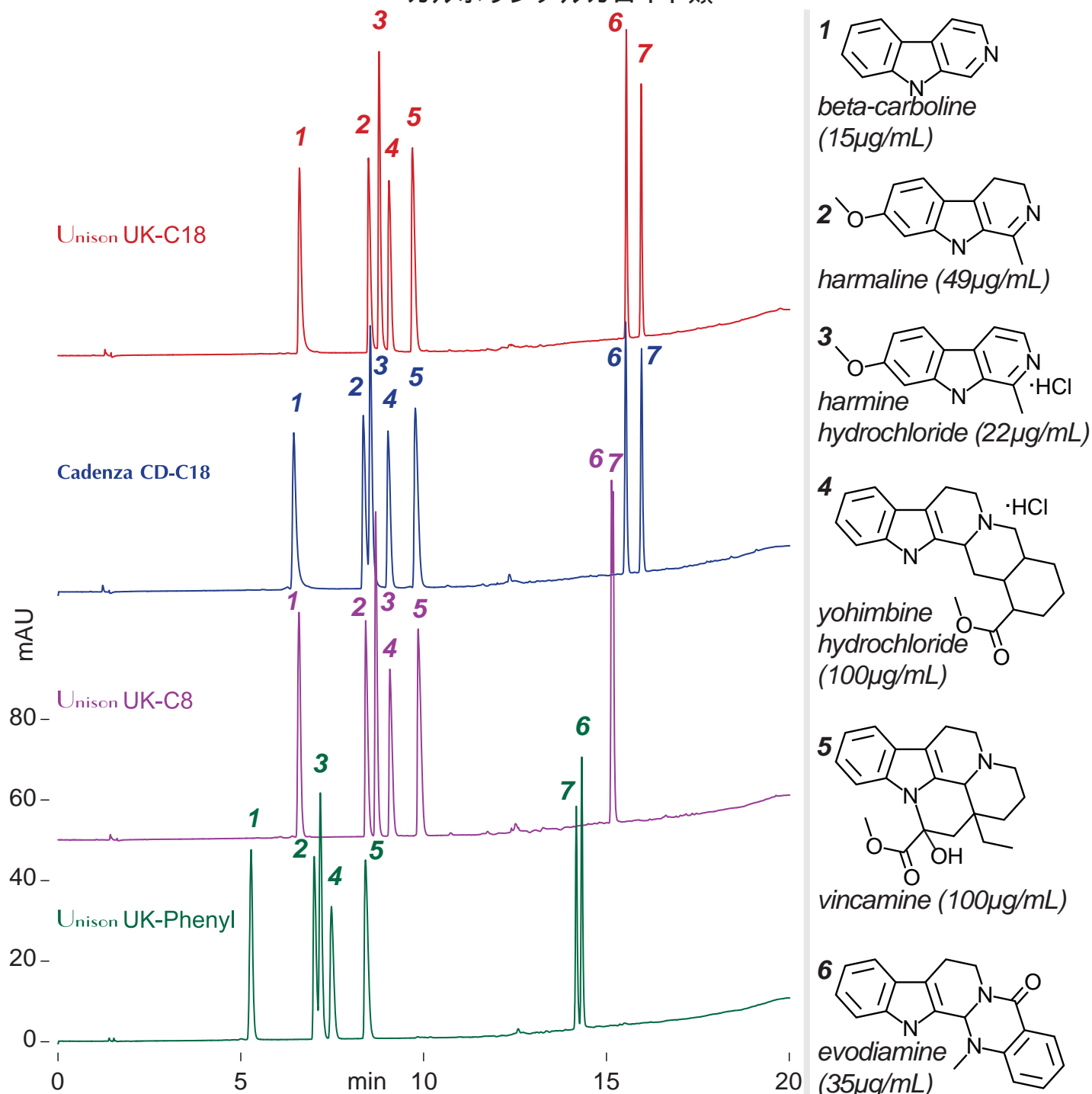

100 x 4.6 mm

A: water / formic acid = 100 / 0.5

B: acetonitrile / water / formic acid = 95 / 5 / 0.5

5 -30%B (0-10min), 30-100%B% (10-18min), 100%B (18-20min)

1mL/min (9-12MPa), 37°C, 260nm

 5 μ L (15-100 μ g/mL, ACN / water / HCOOH = 50 / 50 / 1)