



Nardis ND-NX
Unison UK-Amino

150 x 4.6 mm

HPLC COLUMN APPLICATION

Gastrodin (Phenolic glycosides)

ガストロジン（フェノール配糖体）

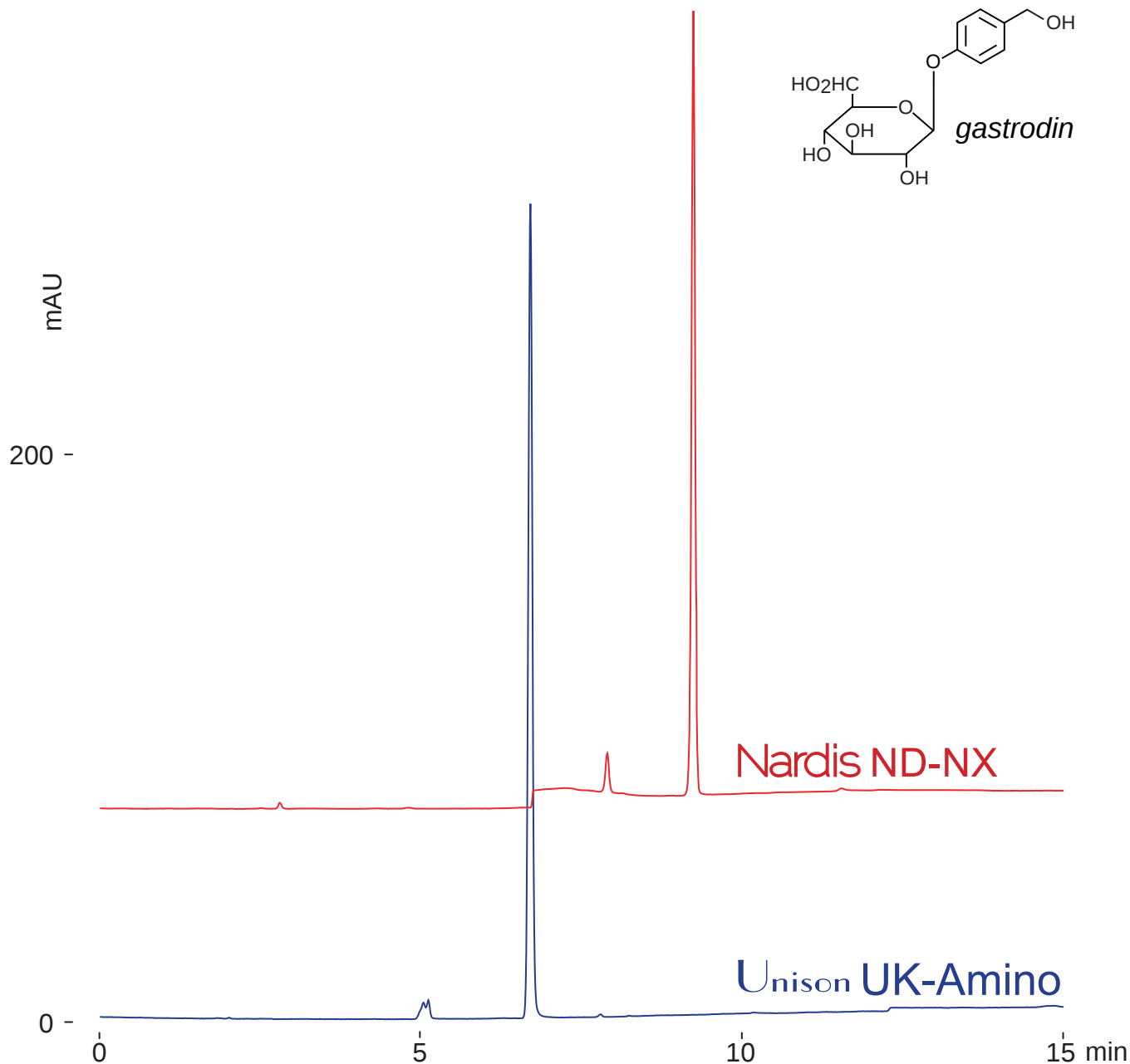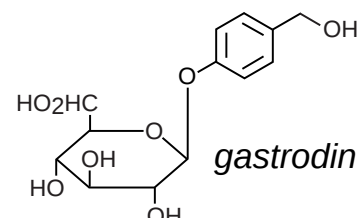

150 x 4.6 mm

A: acetonitrile

B: acetonitrile / 100mM ammonium formate = 50 / 50

5-75%B (0-15min)

1mL/min (5MPa), 37°C, 270nm

5μL (2.0mg/mL, 100mM ammonium acetate / acetonitrile = 50 / 50)