


Presto FF-C18

150 x 3 mm

HPLC COLUMN APPLICATION

Non-Porous ODS-MS: Bovine serum albumin (BSA)

ウシ血清アルブミン (BSA) のノンポーラス ODS カラム -MS 分析

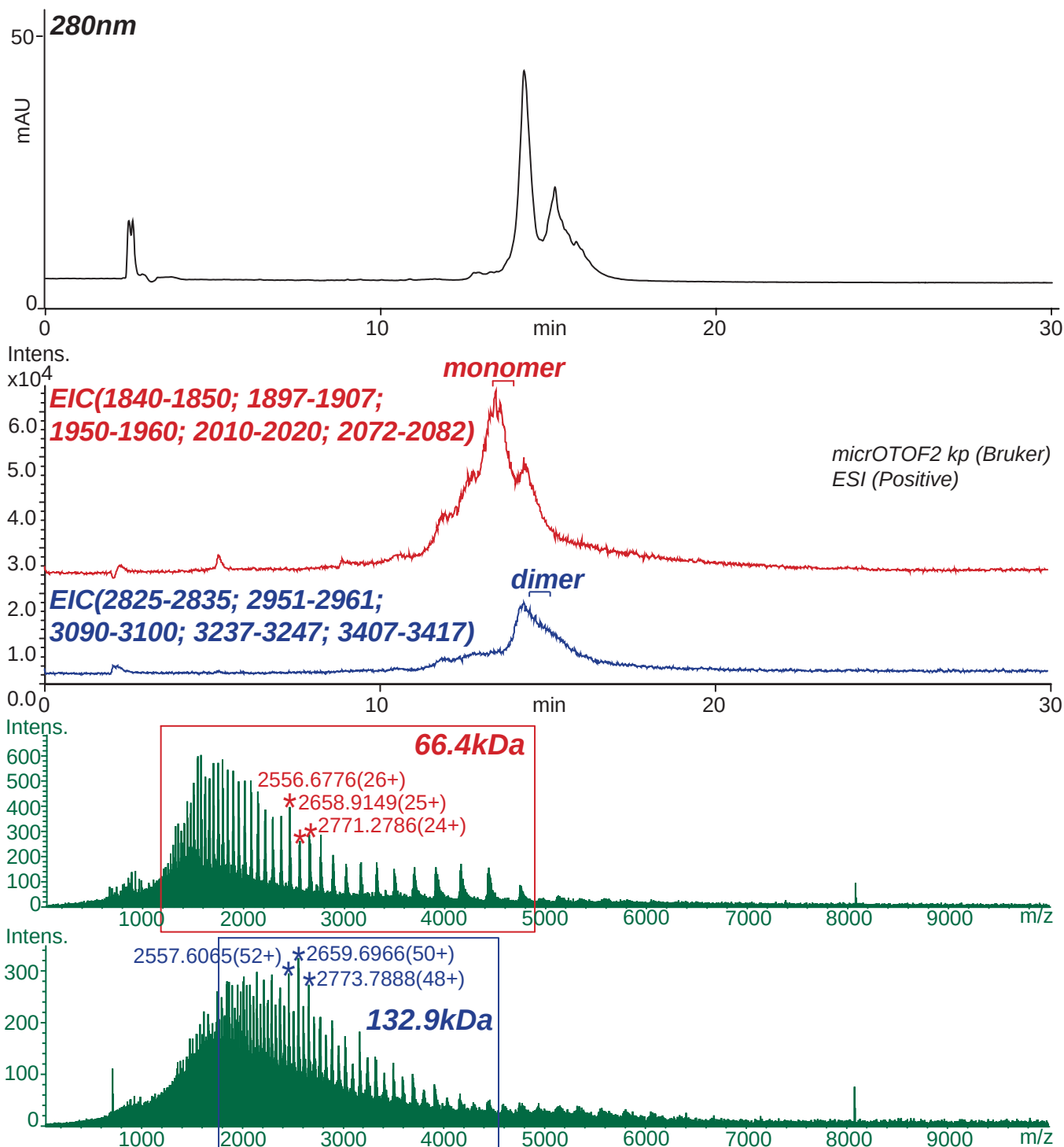

Presto FF-C18, 150 x 3mm

A: 100mM ammonium formate / acetonitrile = 95 / 5

B: 100mM ammonium formate / acetonitrile = 20 / 80

25-70%B (0-30min)

0.2mL/min (26MPa), 37°C

5uL (10mg/mL, 100mM ammonium formate)