


Presto FF-C18

150 x 3 mm

HPLC COLUMN APPLICATION

Non-Porous ODS-MS:Transferrin (holo)
トランスフェリン (ホロ) のノンポーラスODSカラム・MS分析

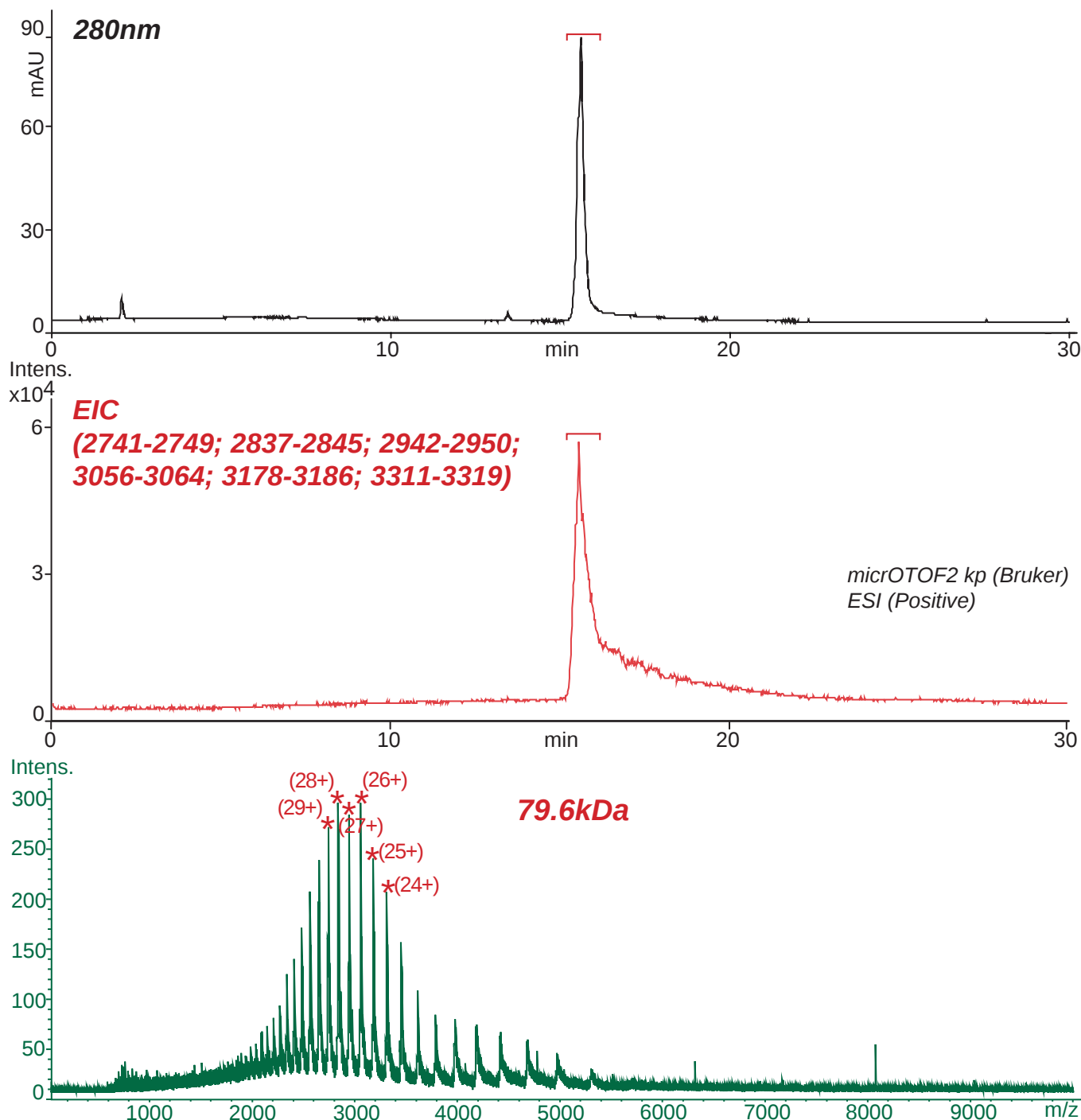

Presto FF-C18, 150 x 3mm

A: 100mM ammonium formate / acetonitrile = 95 / 5

B: 100mM ammonium formate / acetonitrile = 20 / 80

0-100%B (0-30min)

0.2mL/min (24MPa), 37°C

5uL (5mg/mL, 100mM ammonium formate)