


Intrada Organic Acid 150 x 2 mm

HPLC COLUMN APPLICATION

LC-MS: trimellitic acid isomers and relative compounds

トリメリット酸異性体及び関連化合物のLC-MS分析

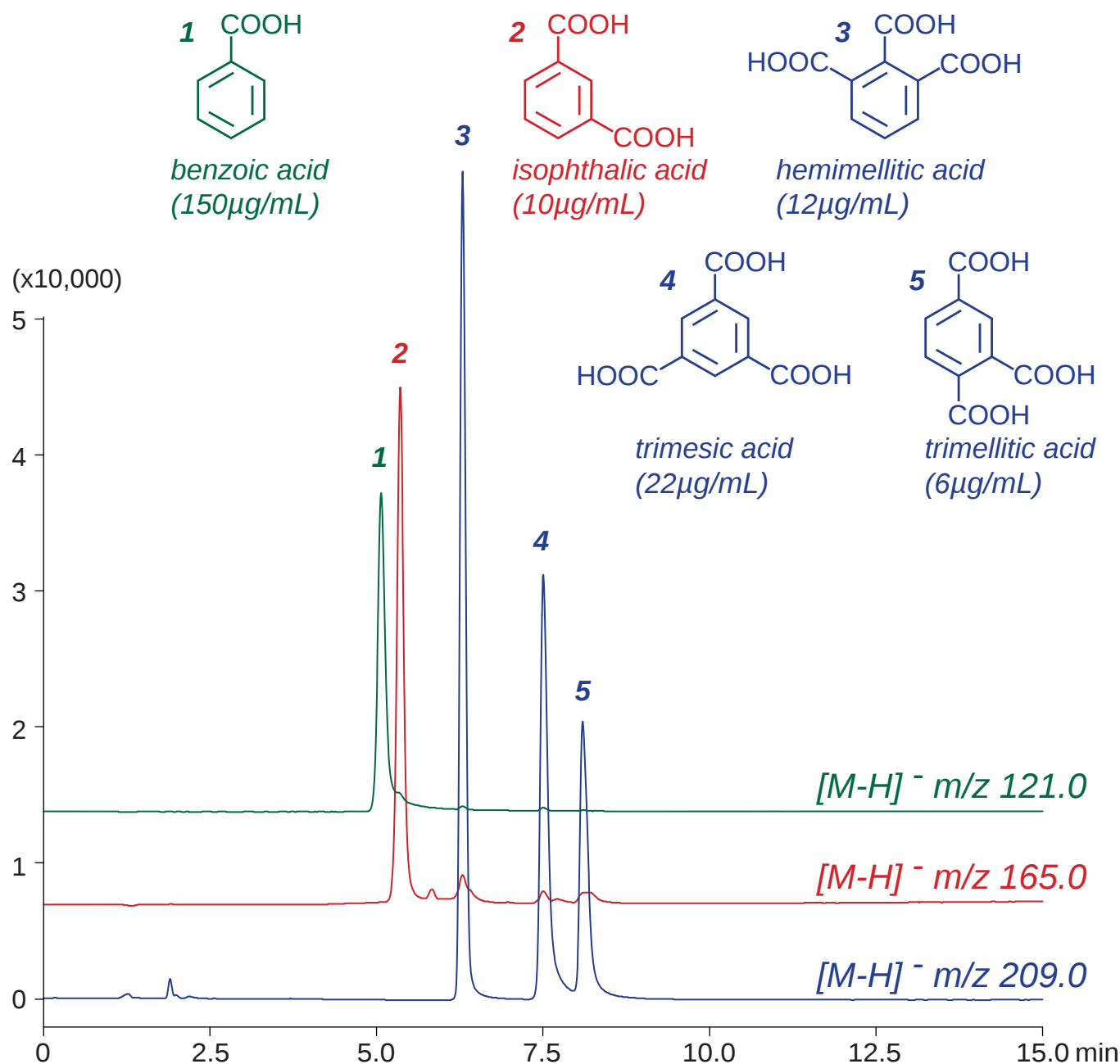

Intrada Organic Acid, 150 x 2 mm

A: water / acetonitrile / formic acid = 90 / 10 / 0.1

B: 100mM ammonium formate / acetonitrile = 80 / 20

0%B (0-1min), 0-70%B (1-12min), 100%B (12-15min)

0.2mL/min (11MPa), 37°C, Single Quad. MS (ESI)

2 μ L (6-150 μ g/mL, 2% formic acid)