

Unison UK-C8

Unison UK-C18


Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

Lilial (Synthetic fragrance)

リリアール (合成香料)

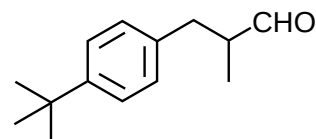

3-(4-tert-butylphenyl)
-2-methylpropanal

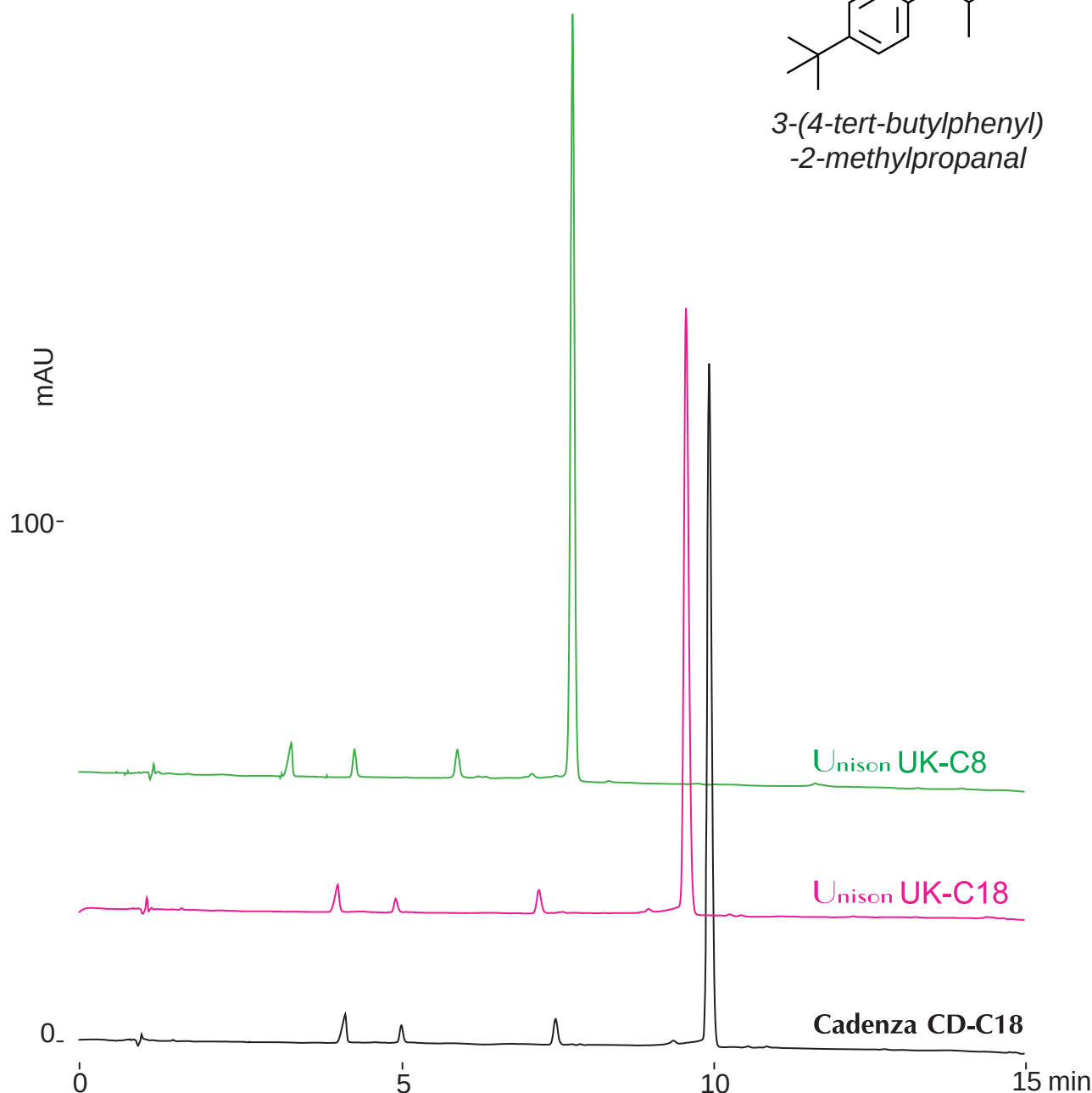

100 x 4.6 mm

A: 100mM ammonium formate / acetonitrile = 90 / 10

B: acetonitrile

35-100%B (0-15min)

1.0mL/min (13-15MPa), 37°C, 260nm

10μL (2mg/mL, 100mM ammonium formate / acetonitrile = 40 / 60)