



Scherzo SS-C18

150 x 3 mm
250 x 3 mm



Application

Neurotransmitters

神経伝達物質

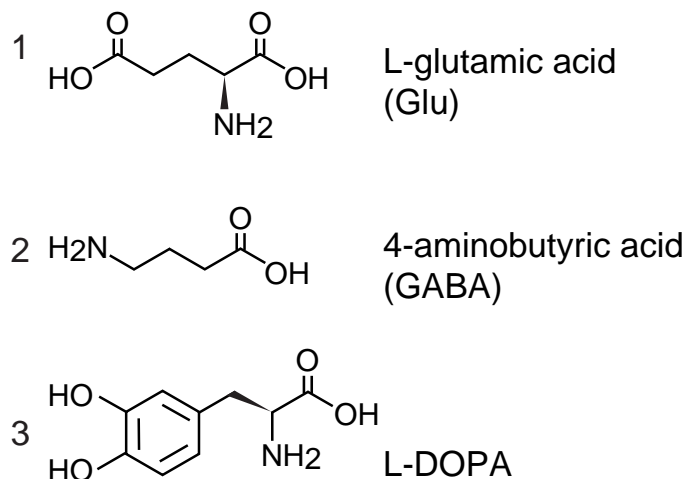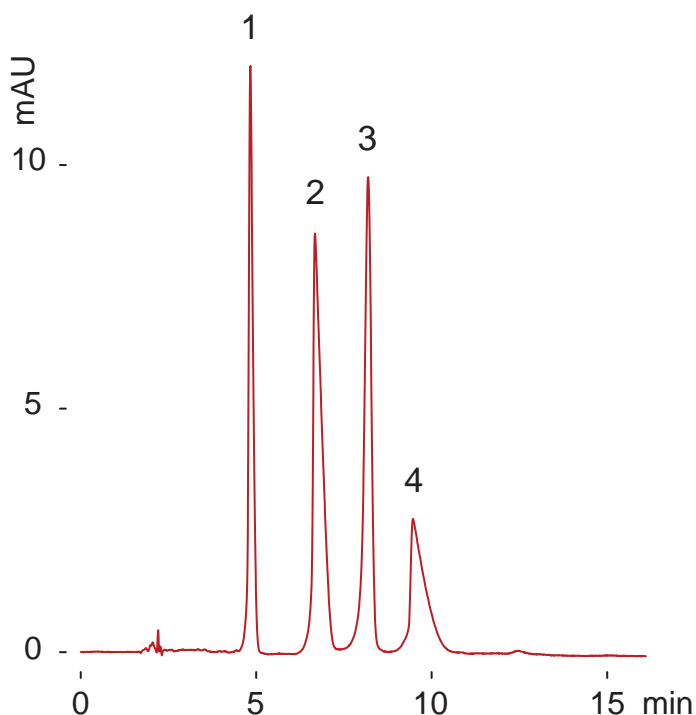

Scherzo SS-C18, 150 x 3 mm
20 mM H₃PO₄ / acetonitrile = 80 / 20
0.4 mL/min (11 MPa), 37 deg.C, 220 nm
5 uL (0.06-5 ug, 40 mM H₃PO₄)

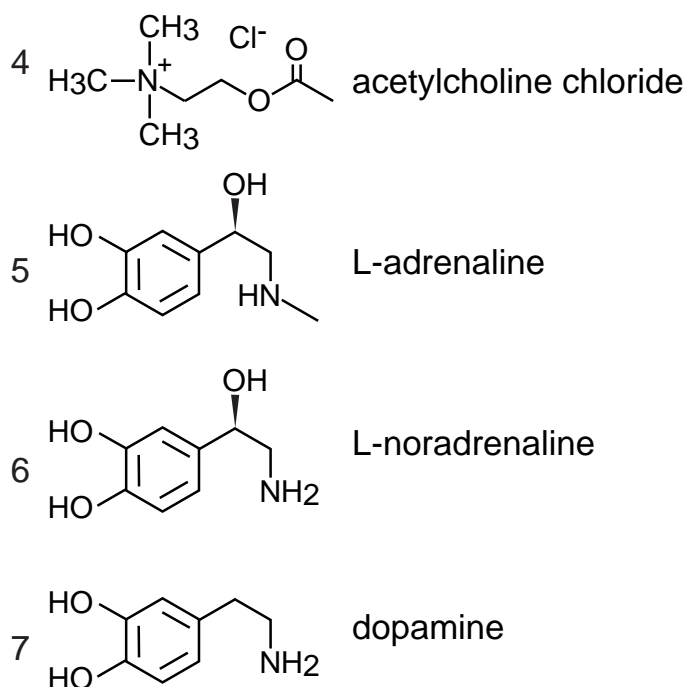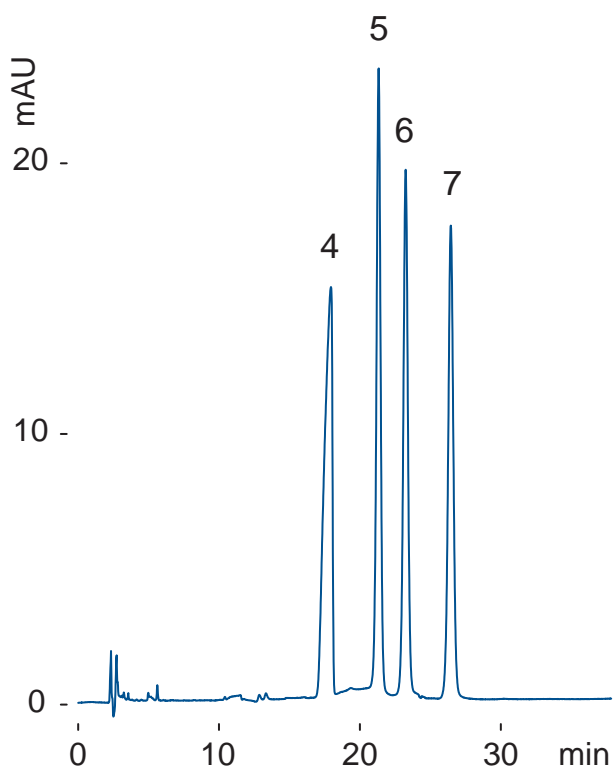

Scherzo SS-C18, 250 x 3 mm
10 mM NH₄H₂PO₄ / acetonitrile = 40 / 60
0.4 mL/min (12 MPa), 37 deg.C, 220 nm,
4 uL (0.09-40 ug)