




Intrada Amino Acid

150 x 3 mm

Application

LC-MS: Valine isomers

バリン異性体のLC-MS分析

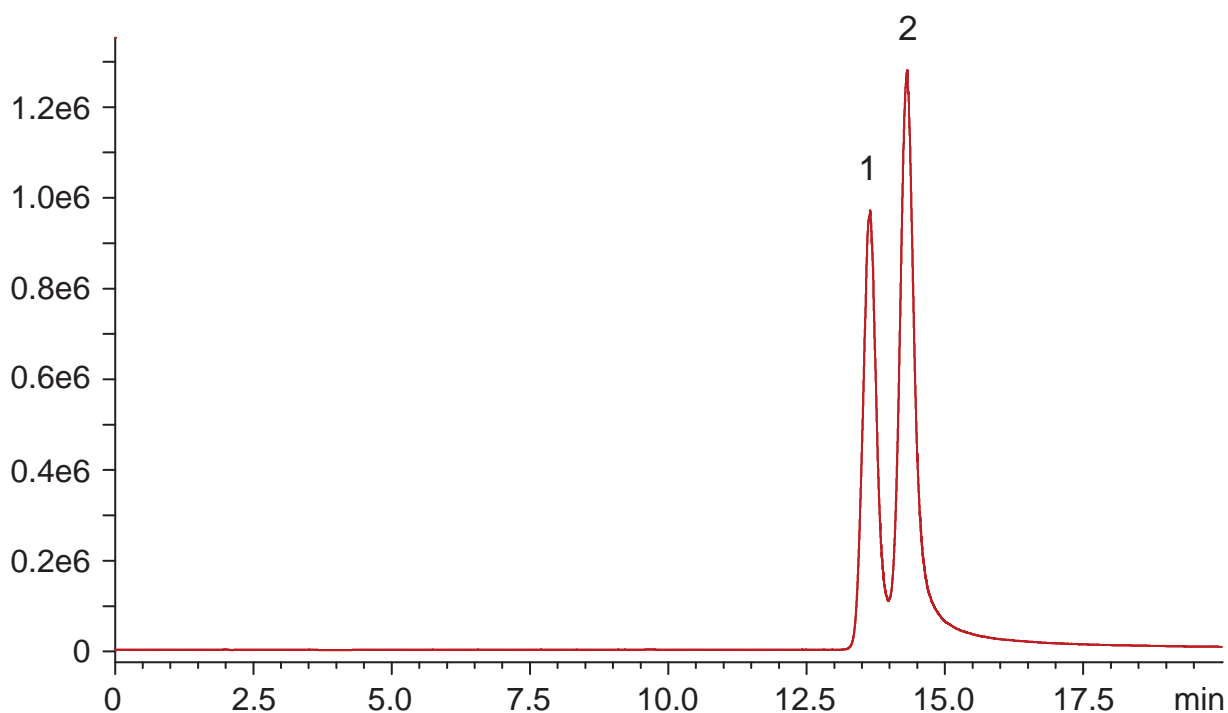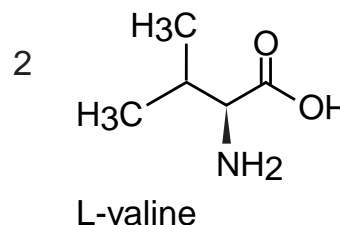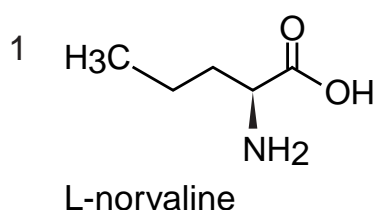

Intrada Amino Acid, 150 x 3 mm

A: acetonitrile / water / formic acid = 85 / 15 / 0.3

B: 100mM ammonium formate

5-8.5 %B (0-18 min), 100 %B (18-20 min)

0.4 mL/min (6 MPa), 25°C, 5uL (1mM, 0.1N HCl)

Single Quad.MS (ESI, positive, m/z 118.2)