


Scherzo SM-C18
Scherzo SS-C18


150 x 3 mm

Application

Sodium Halides
 ハロゲン化ナトリウム類

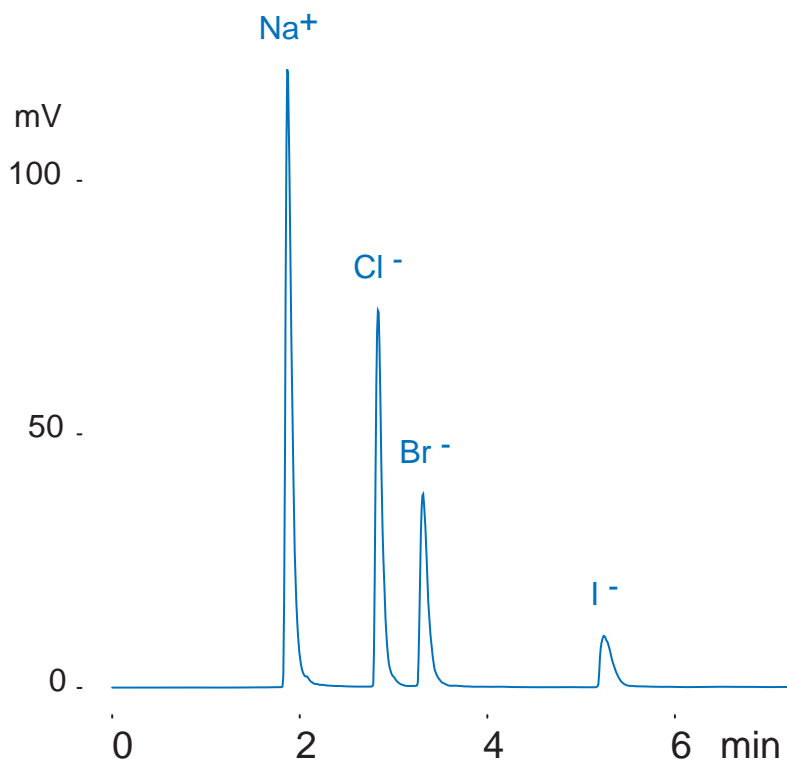

Scherzo SM-C18, 150 x 3 mm
 20 mM ammonium formate
 / 20 mM formic acid = 95 / 5
 0.4 mL / min (12 MPa), 37 deg.C
 ELSD (spray chamber 20 deg.C
 drift tube 45 deg.C)
 2.0 uL (0.7 ug as sodium salts)

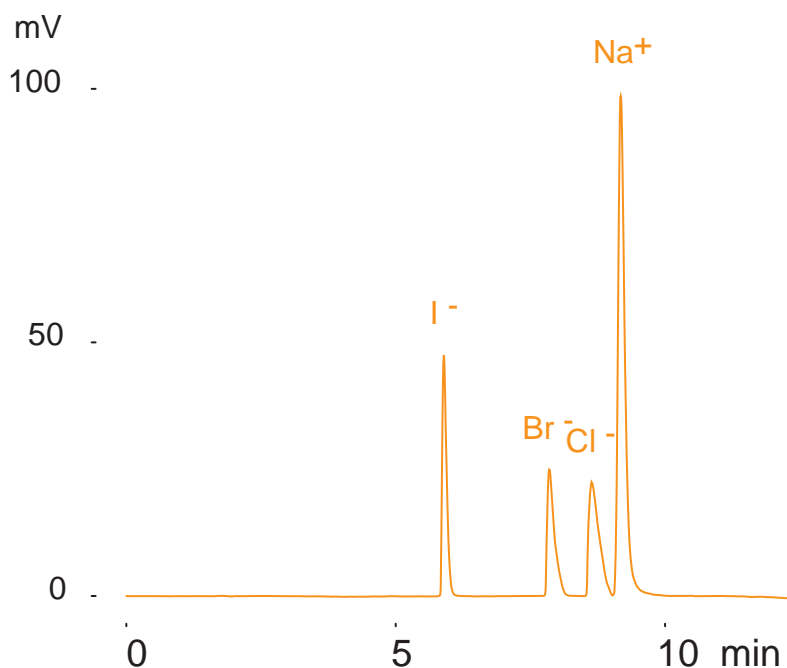

Scherzo SS-C18, 150 x 3 mm
 A: 50 mM ammonium acetate
 / acetonitrile = 15 / 85
 B: 100 mM ammonium acetate
 / acetonitrile = 20 / 80
 0 - 100 %B (0 - 10 min)
 0.4 mL / min (6 MPa), 37 deg.C
 ELSD (spray chamber 20 deg.C
 drift tube 45 deg.C)
 4.0 uL (1.4 ug as sodium salts)