

Unison UK-Phenyl Cadenza CX-C18



Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Clotiazepam (Antianxiety, GABA-A receptor agonist)

クロチアゼパム (抗不安薬, GABA-A 受容体作動薬)

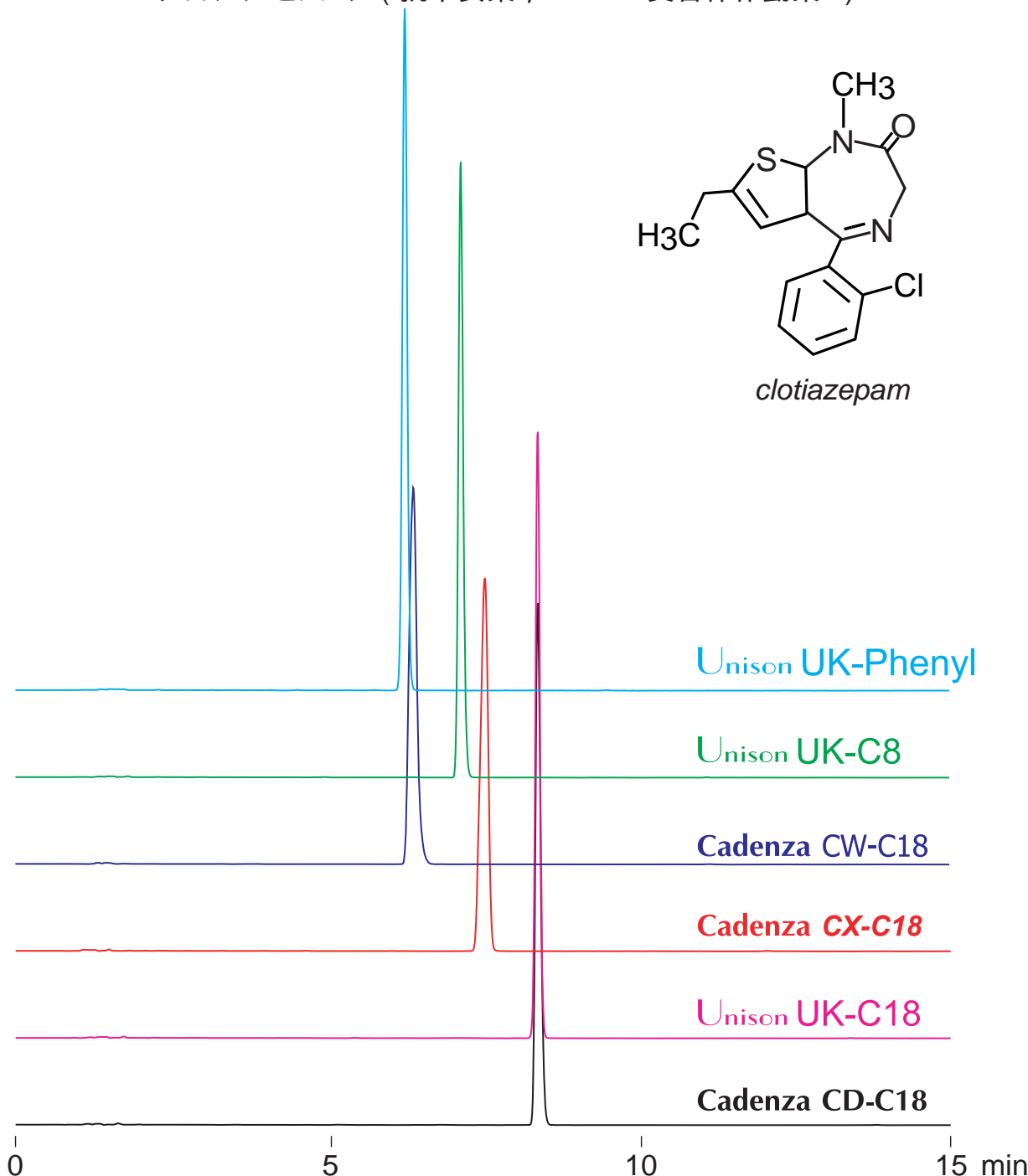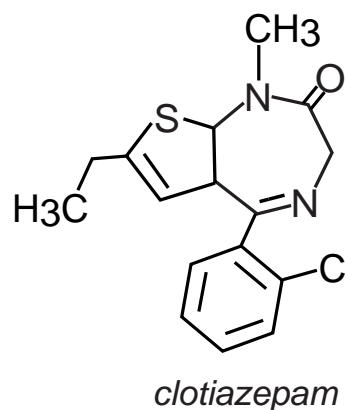

100 x 4.6 mm

A: 100mM ammonium formate

B: acetonitrile / 100mM ammonium formate = 90 / 10

35-100%B (0-15min)

1.0 mL/min (7MPa), 37°C, 260nm

10µL (0.25mg/mL, methanol / formic acid = 100 / 1)