


**DACAPO DX-C18**  
**Cadenza CX-C18**

100 x 4.6 mm

**HPLC COLUMN APPLICATION**

**Basic Blue ( Dye )**  
 ベーシックブルー 26 ( 染料 )

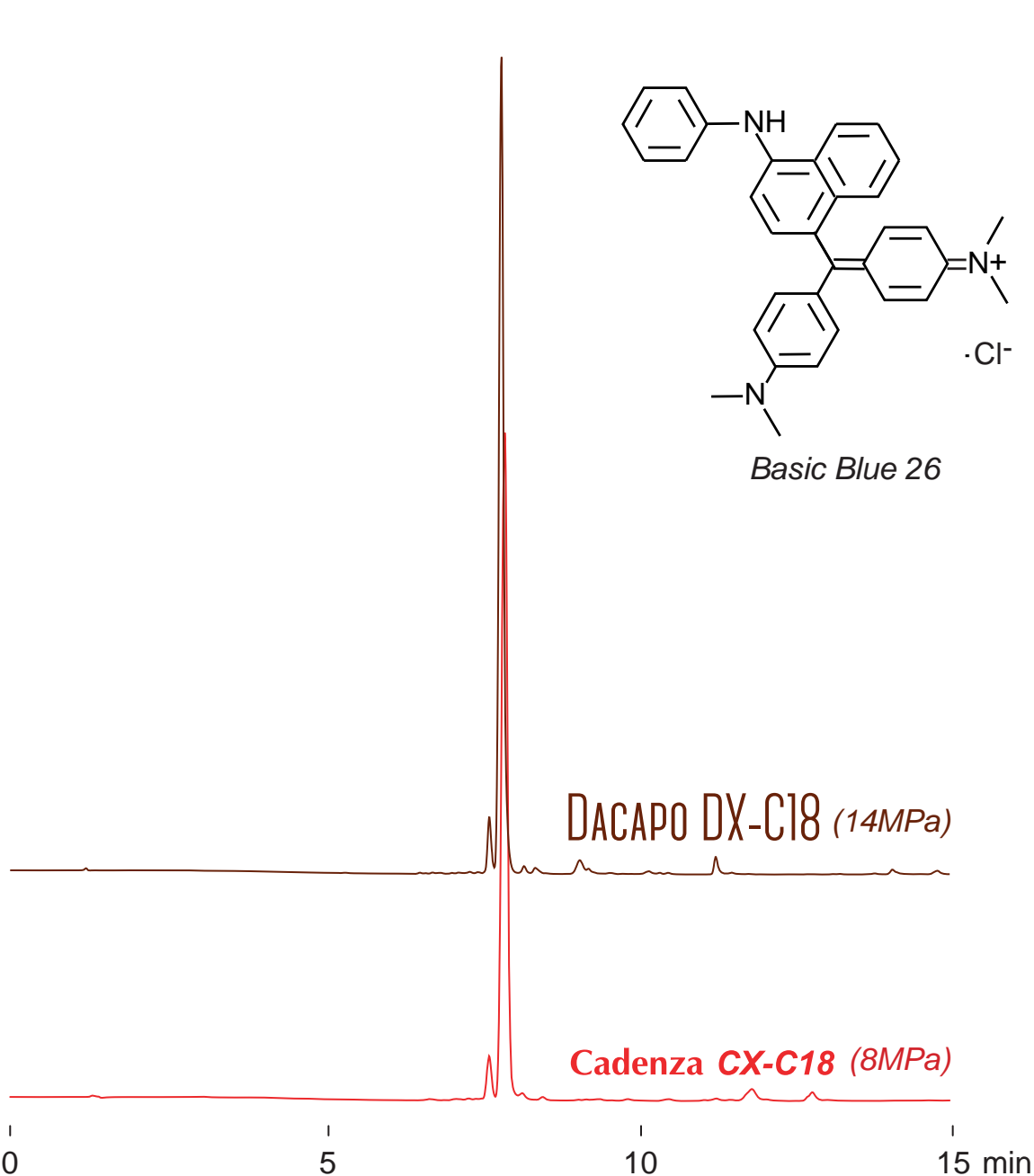

100 x 4.6 mm  
 A: water / formic acid = 100 / 0.1  
 B: acetonitrile / formic acid = 100 / 0.1  
 0-100%B (0-15min)  
 1.0mL/min, 37°C, 260nm  
 10μL (1mg/mL, 100mM ammonium formate)