

Unison UK-Phenyl Cadenza CW-C18

Unison UK-C8

Cadenza CX-C18

Unison UK-C18



Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

## Sunset Yellow ( Food additive )

黄色 5号 ( 食品添加物 )

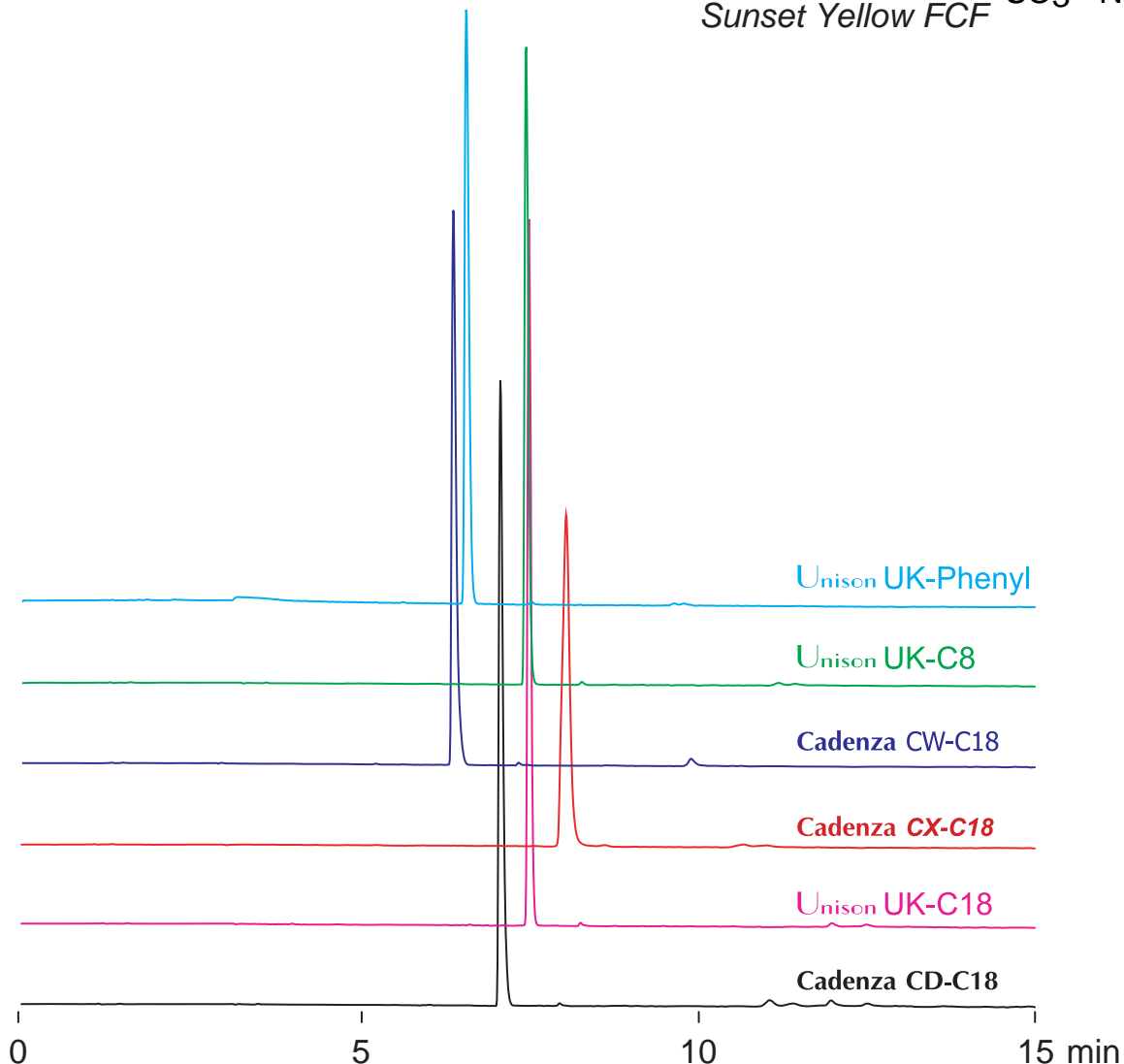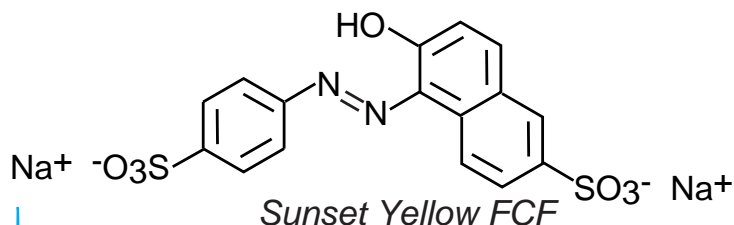

100 x 4.6 mm

A: 100mM ammonium acetate

B: acetonitrile / 100mM ammonium acetate = 90 / 10

0-50%B (0-15min)

1mL/min (9MPa), 37°C, 260nm

10µL (0.25mg/mL, 100mM ammonium formate)