



Cadenza CW-C18 Unison UK-C1
 Cadenza CX-C18 Unison UK-Phenyl
 Cadenza CD-C18 Unison UK-C8
 Unison UK-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Paclitaxel (Antineoplastic)

パクリタキセル (抗悪性腫瘍薬)

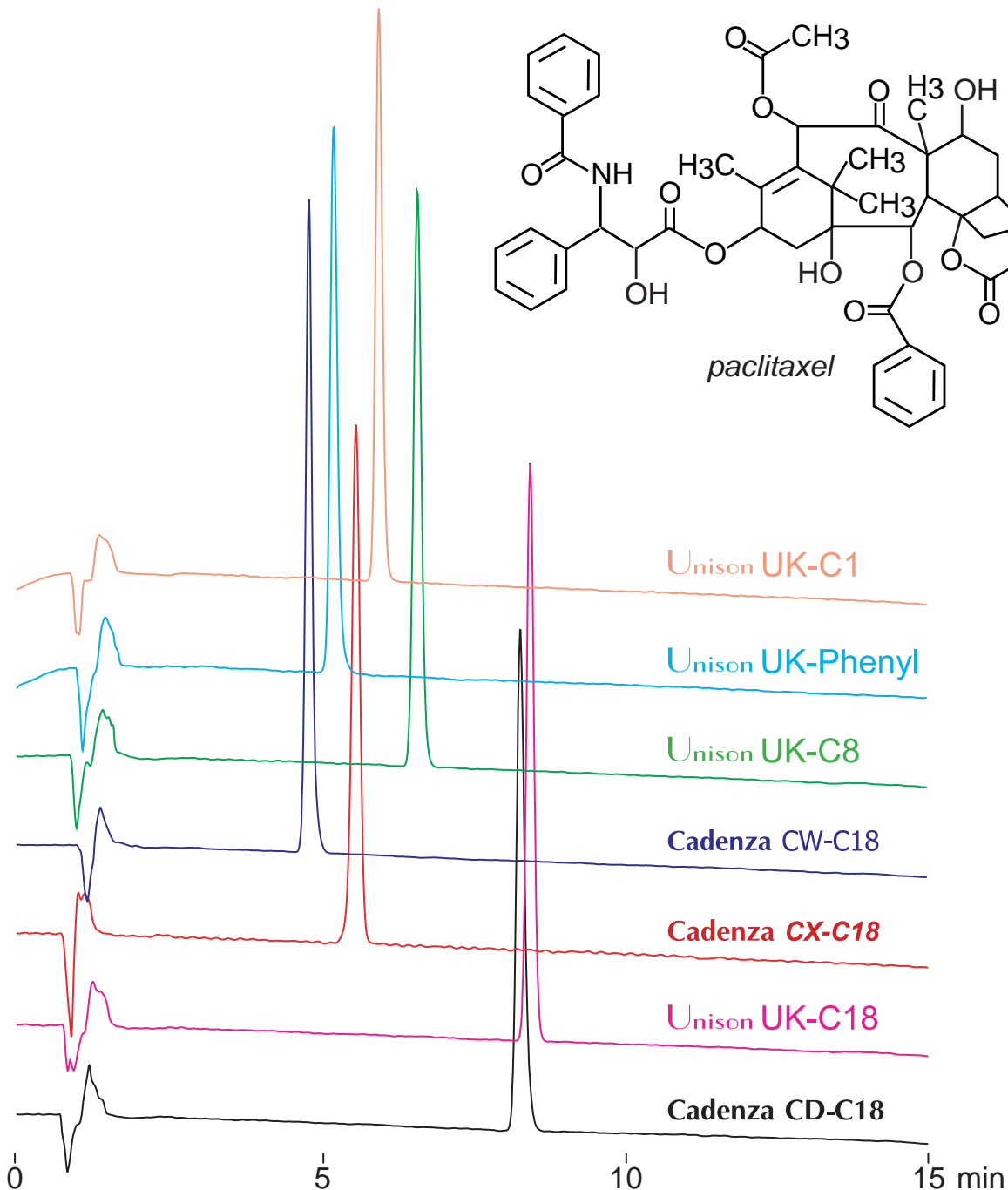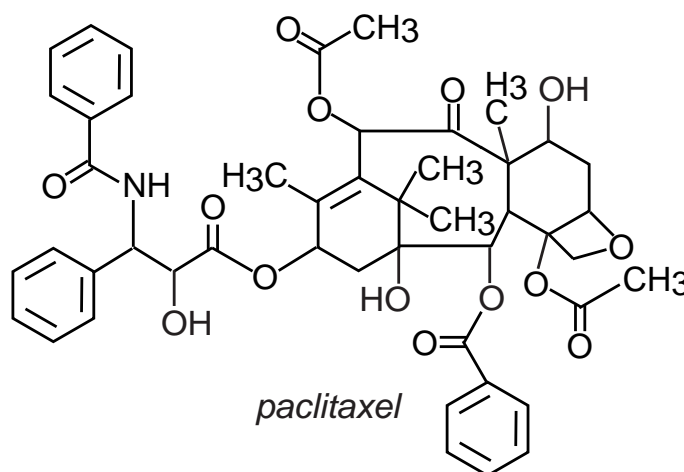

100 x 4.6 mm

A: 100mM ammonium formate

B: acetonitrile / 100mM ammonium formate = 90 / 10

50-65%B (0-15min)

1.0mL/min (8MPa), 37°C, 230nm

10µL (0.1mg/mL, methanol)