


Cadenza CX-C18

150 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Organophosphorus pesticides (3)

有機リン系殺虫剤 (3) の LC-MS 分析

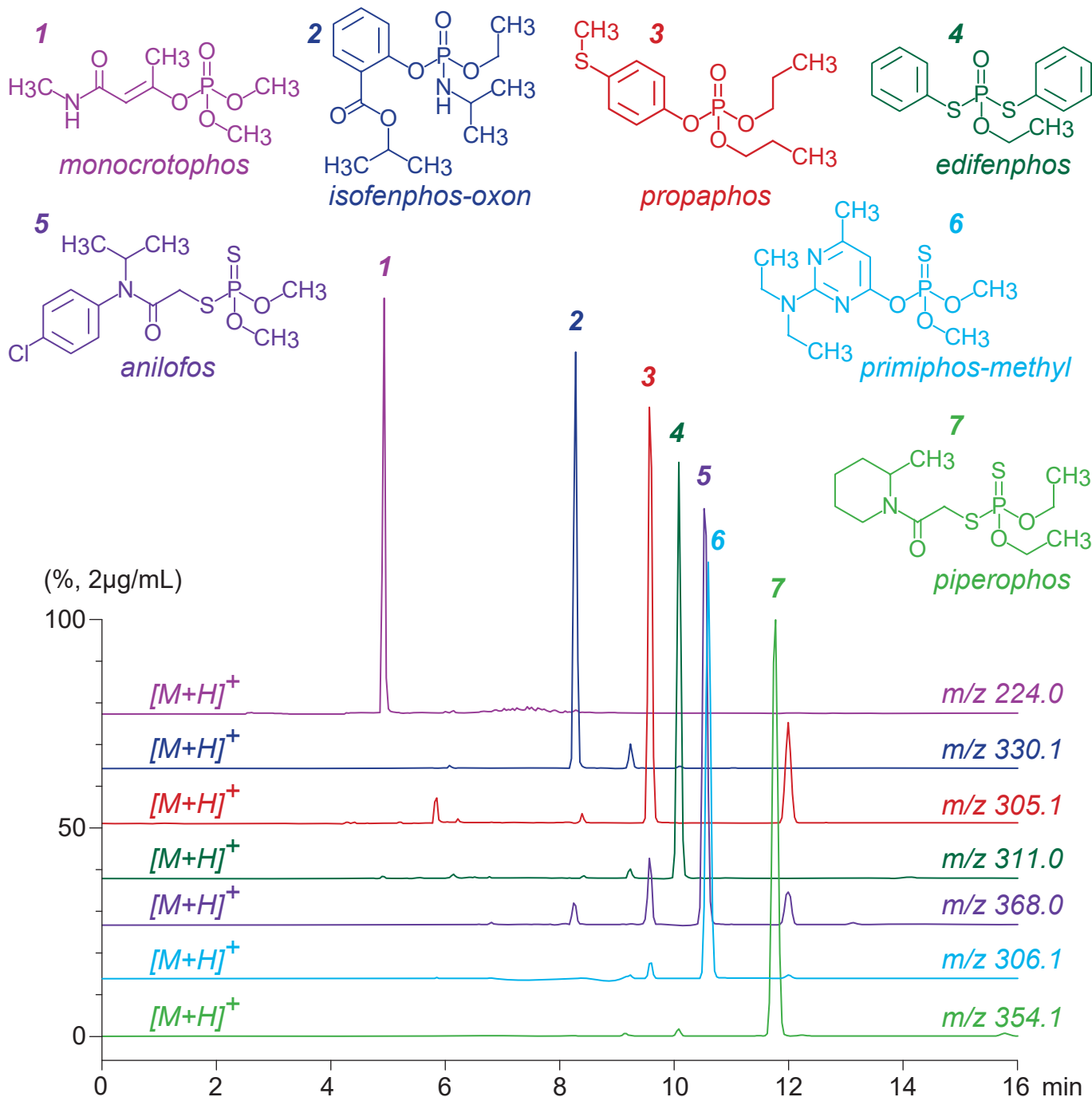

Cadenza CX-C18, 150 x 3 mm
 A: water / formic acid = 100 / 0.1
 B: acetonitrile / formic acid = 100 / 0.1
 0-65%B (0-5min), 65%B (5-16min)
 0.4mL/min (9MPa), 37°C, Single Quad. MS (ESI)
 5µL (2µg/mL, ACN / water / HCOOH = 50 / 50 / 0.1)