


Intrada WP-RP

150 x 4.6 mm

HPLC COLUMN APPLICATION

Polyoxethylene lauryl ether (Nonionic surfactant)

ポリオキシエチレンラウリルエーテル (非イオン性界面活性剤)

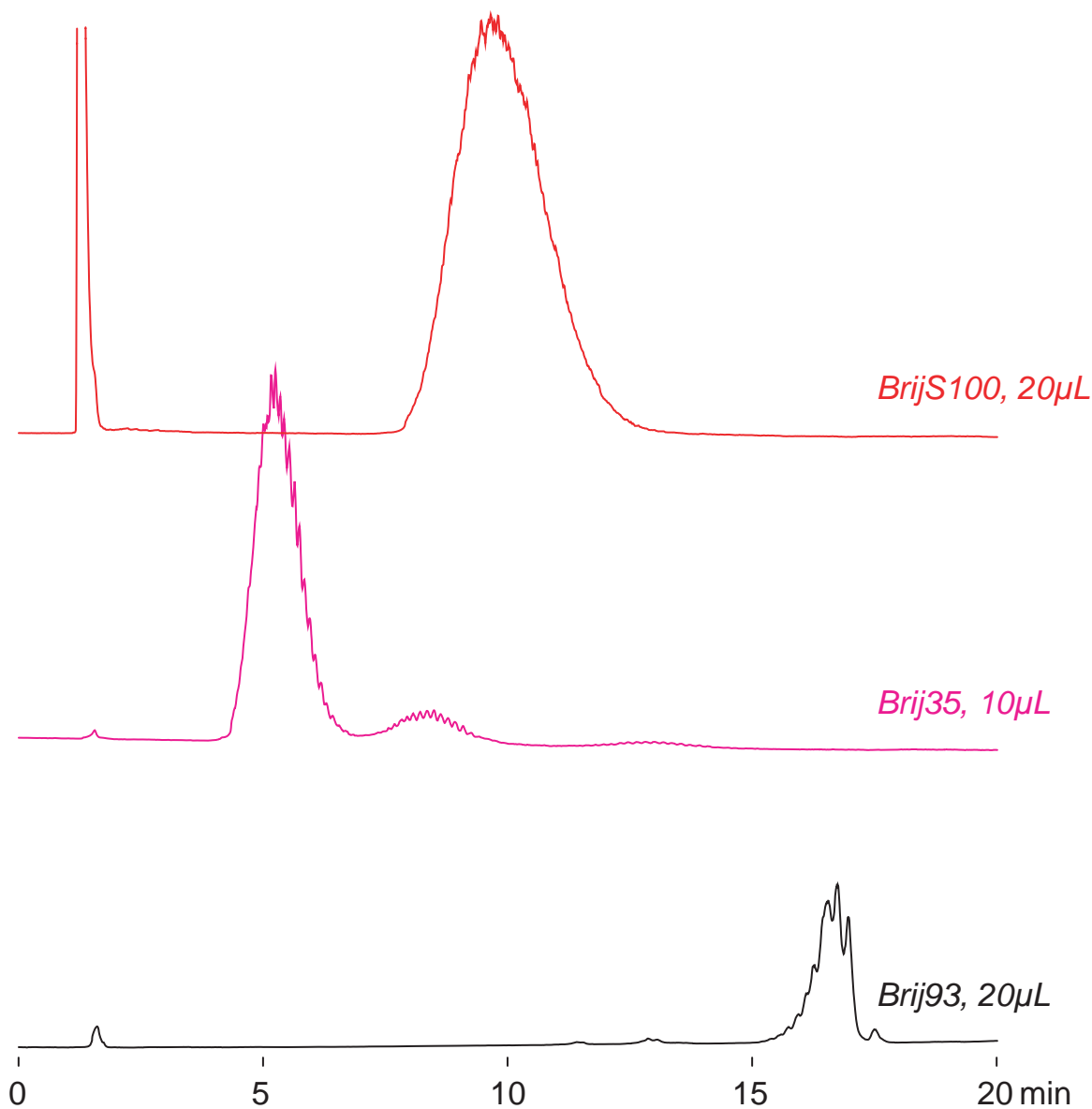

Intrada WP-RP, 150 x 4.6 mm
 A: water / acetic acid = 100 / 0.5
 B: acetonitrile / acetic acid = 100 / 0.5
 80-100%B (0-20min)
 1.0mL/min (7MPa), 37°C
 ELSD (spray chamber 30°C, drift chamber 60°C)
 10-20µL (5mg/mL, ACN / 100mM HCOONH4 = 50 / 50)