

Scherzo SW-C18

Scherzo SM-C18


Scherzo SS-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

Carbazochrome (Hemostatic)

カルバゾクロム (止血薬)

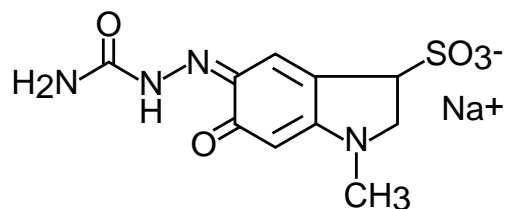

carbazochrome sodium sulfonate

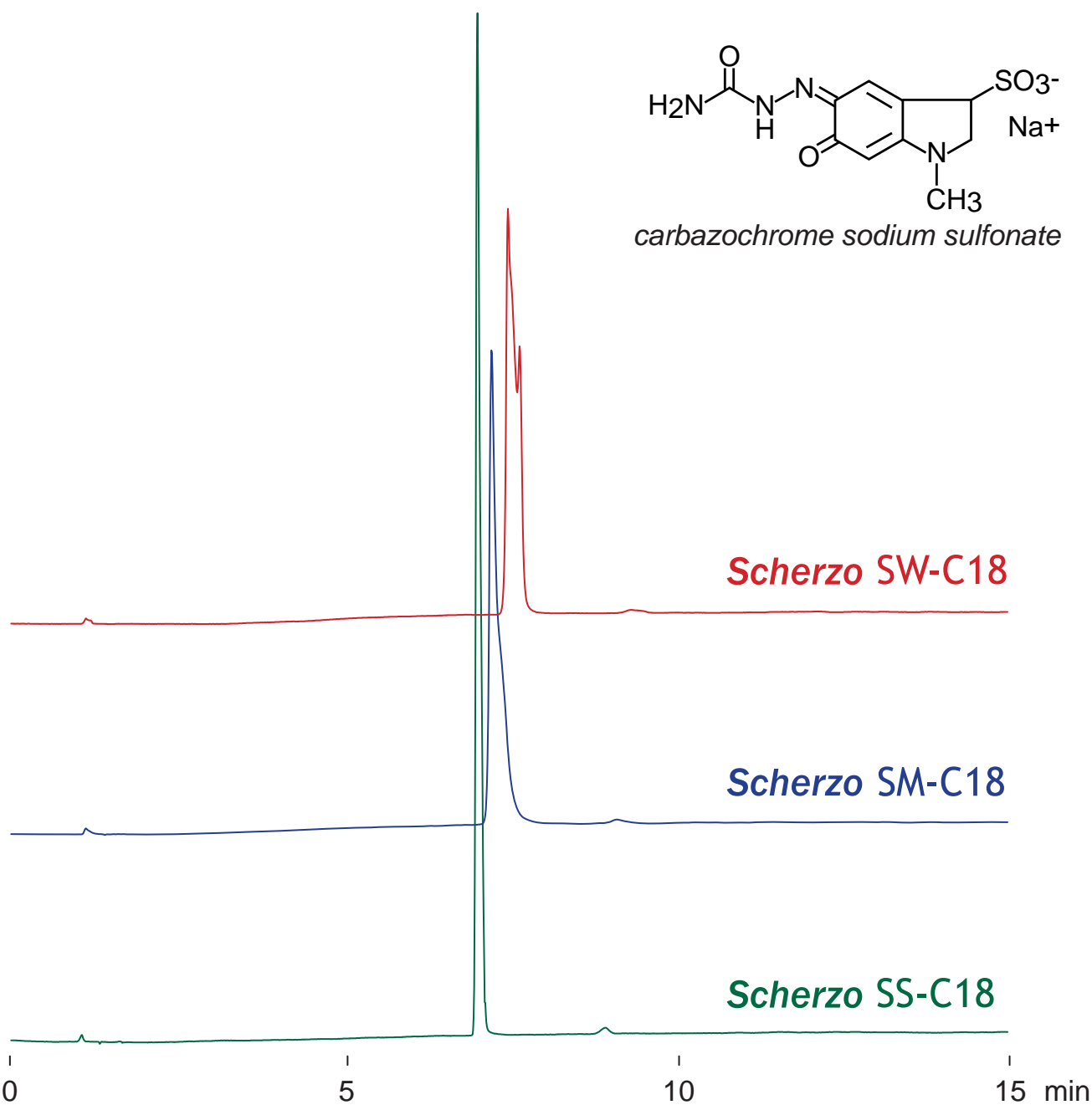

100 x 4.6 mm

A: 10mM ammonium acetate

B: 100mM ammonium acetate / acetonitrile = 75 / 25

0-100%B (0-10min), 100%B(10-15min)

1.0 mL/min (12MPa), 37°C, 360nm

10µL (0.05mg/mL, 100mM ammonium acetate)