


Intrada WP-RP

250 x 4.6 mm

HPLC COLUMN APPLICATION

Transferrin

トランスフェリン

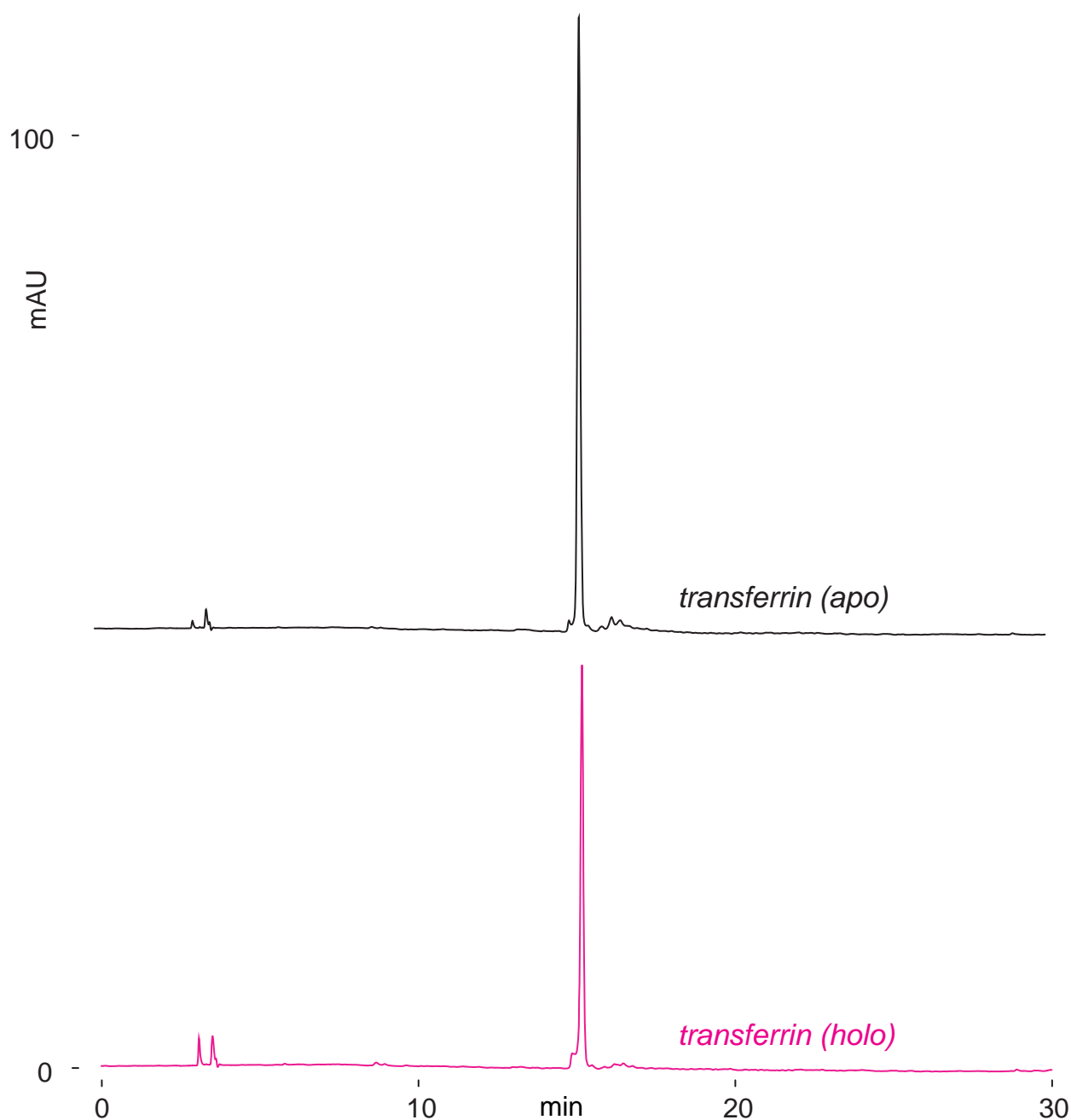

Intrada WP-RP, 250 x 4.6 mm

A: water / acetonitrile / trifluoroacetic acid = 90 / 10 / 0.1

B: water / acetonitrile / trifluoroacetic acid = 10 / 90 / 0.1

0-75%B (0-30min)

1mL/min, 37°C, 280nm

10µL (1mg/mL, water / trifluoroacetic acid = 100 / 0.1)