

Unison UK-Phenyl Cadenza CX-C18



Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Adrafinil (Synthetic stimulant)

アドラフィニル (合成精神刺激薬)

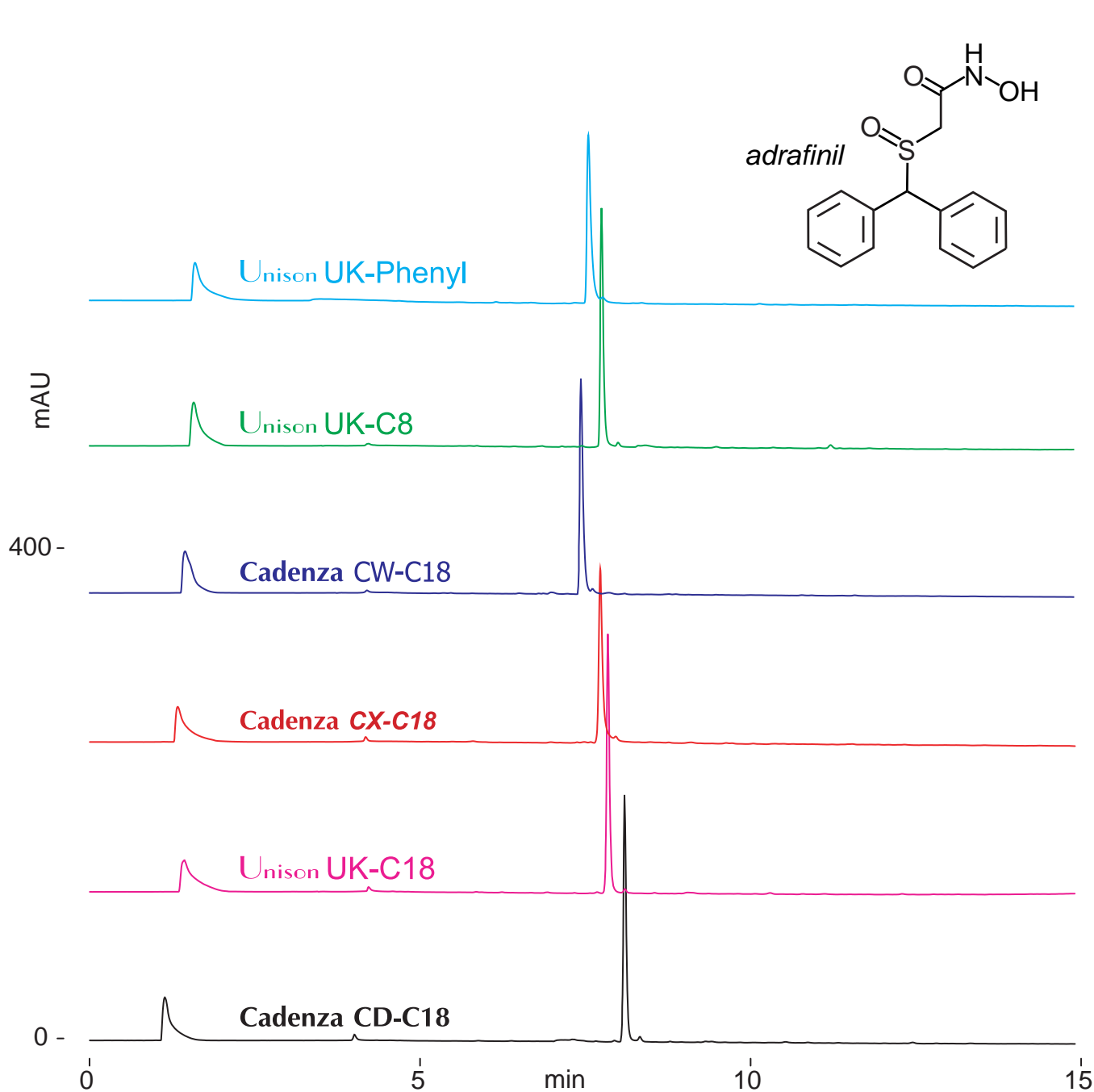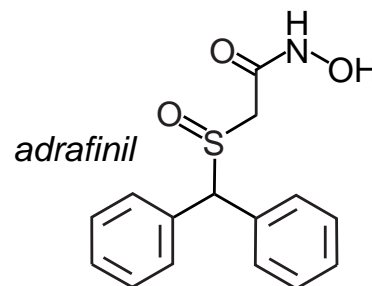

100 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 10 / 90

0-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 250nm

10µL (0.5mg/mL, dimethyl sulfoxide / 100mM ammonium formate = 50 / 50)