

Unison UK-Phenyl

Unison UK-C8

Unison UK-C18


Cadenza CW-C18

Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Quinidine (Antiarrhythmic)

キニジン (抗不整脈薬)

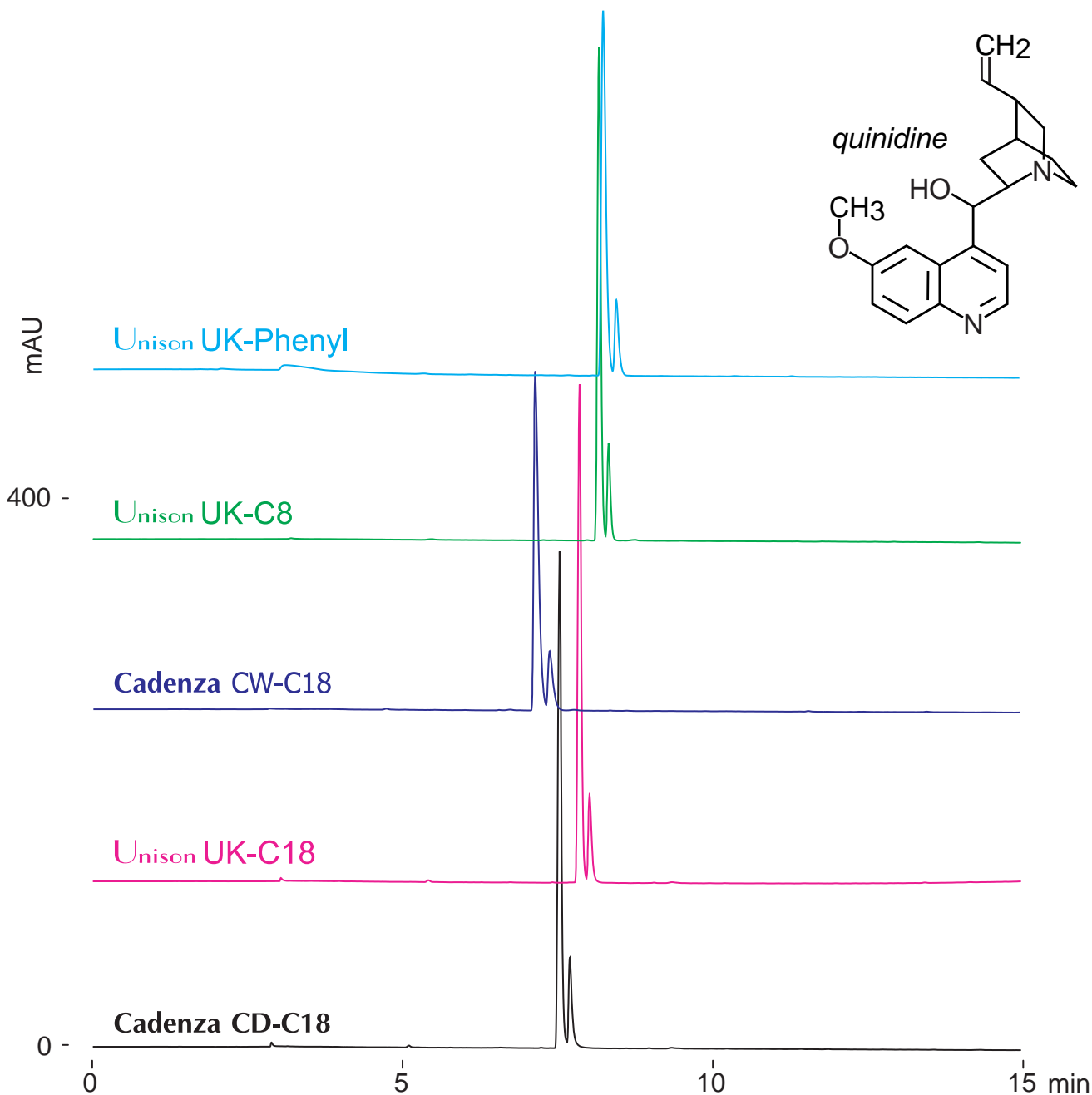

100 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 10 / 90

0-90%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (1mg/mL, 100mM ammonium formate)