
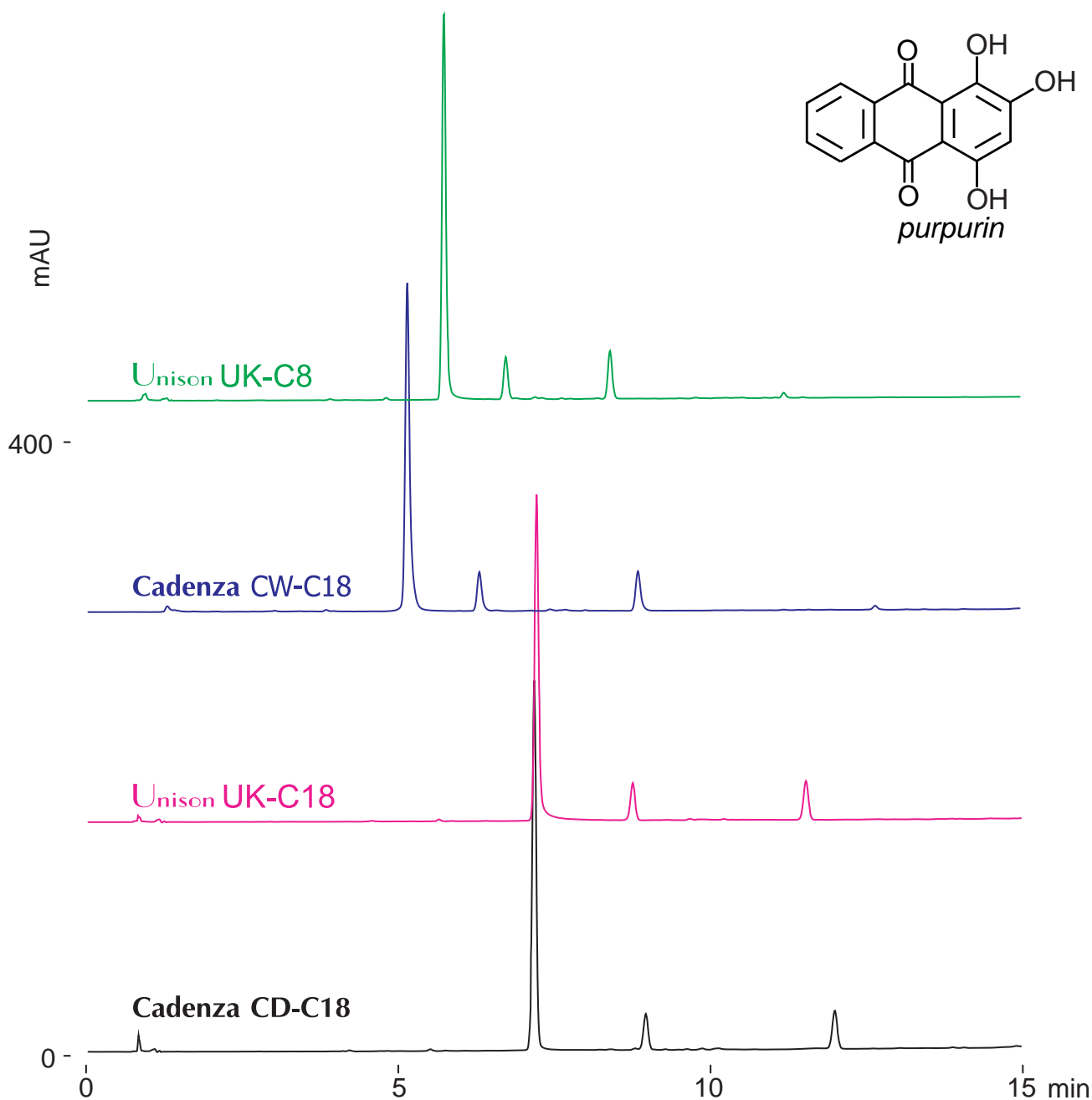


Unison UK-C8
Unison UK-C18

Cadenza CW-C18
Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Purpurin (Dye)
プルプリン (染料)

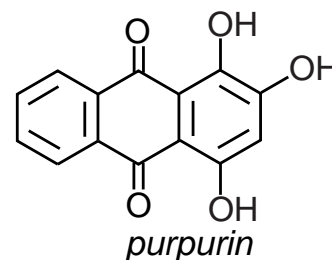

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

40-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.056mg/mL, water / acetonitrile / formic acid = 33 / 67 / 1)