

Cadenza CX-C18


Cadenza CD-C18

100 x 3 mm

HPLC COLUMN APPLICATION


LC-MS: Chlorhexidine (Antibacterial)

クロルヘキシジン (抗菌薬) の LC-MS 分析

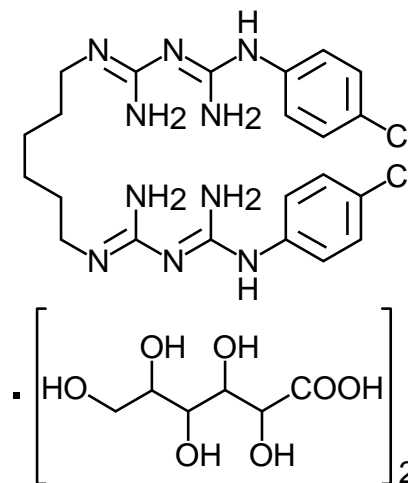

chlorhexidine digluconate

$[M+H]^+$   $m/z$  505.2.0

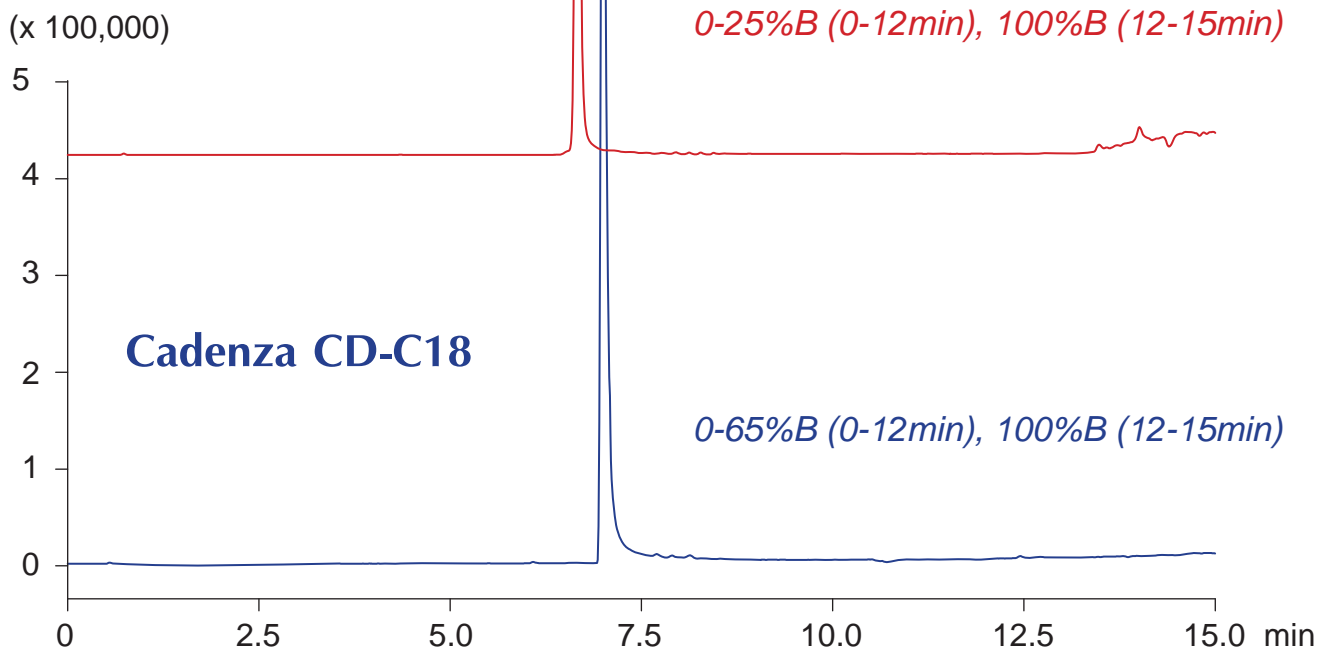

100 x 3 mm

A: acetonitrile / water / formic acid = 10 / 90 / 0.1

B: acetonitrile / water / formic acid = 90 / 10 / 0.1

0.4mL/min (7MPa), 37°C, Single Quad. MS (ESI)

5μL (6μg/mL, water / formic acid = 100 / 1)