

Unison UK-Phenyl Cadenza CX-C18


Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

## Folinic acid (Antifolate)

ホリン酸 (葉酸代謝拮抗薬)

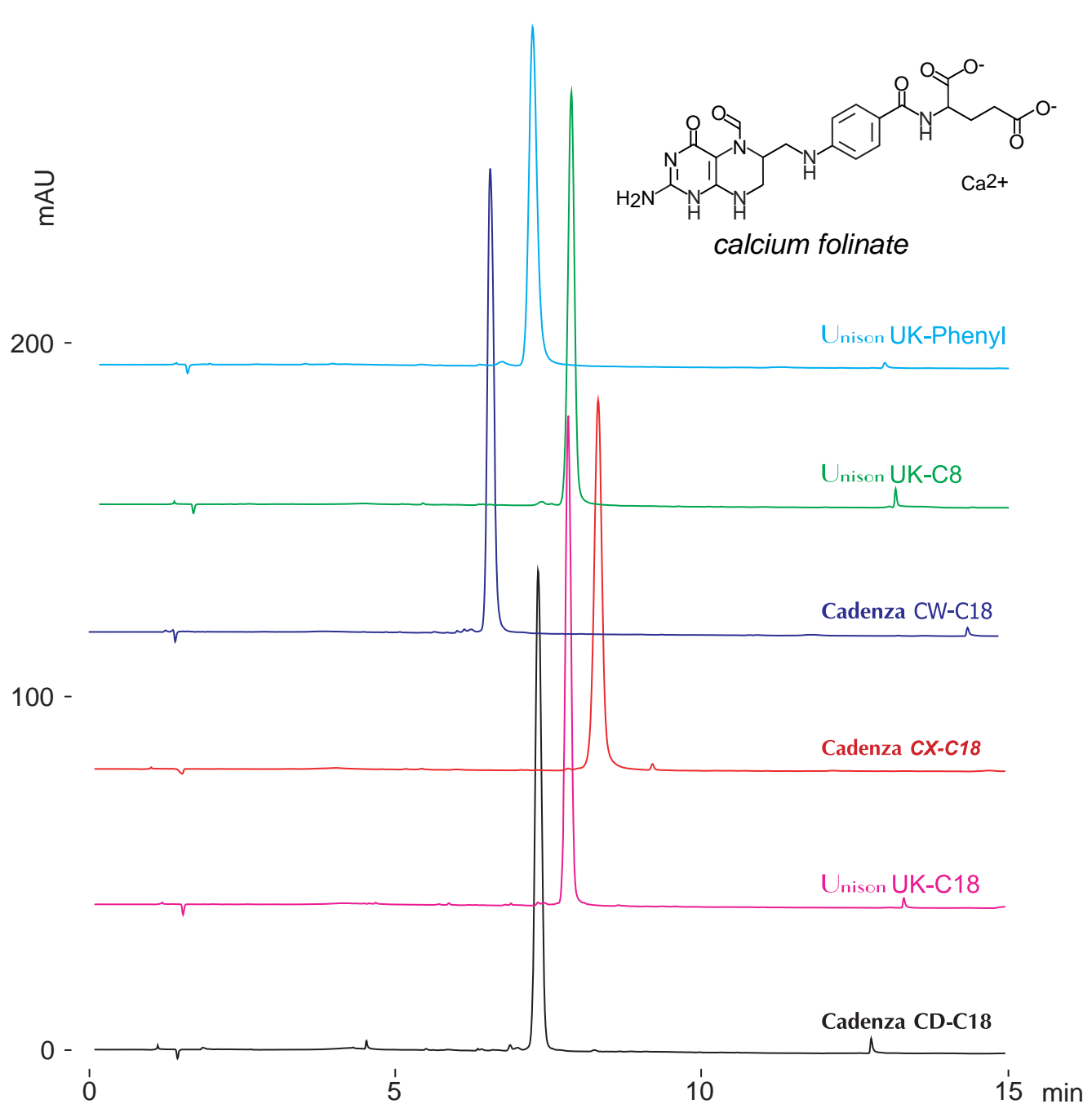

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

0%B (0-1min), 0-45%B (1-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.15mg/mL, 100mM NH<sub>4</sub>OH)