


DACAPO DX-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Noscapine (Antitussive, Isoquinoline Alkaloids)

ノスカピン (鎮咳薬, イソキノリンアルカロイド)

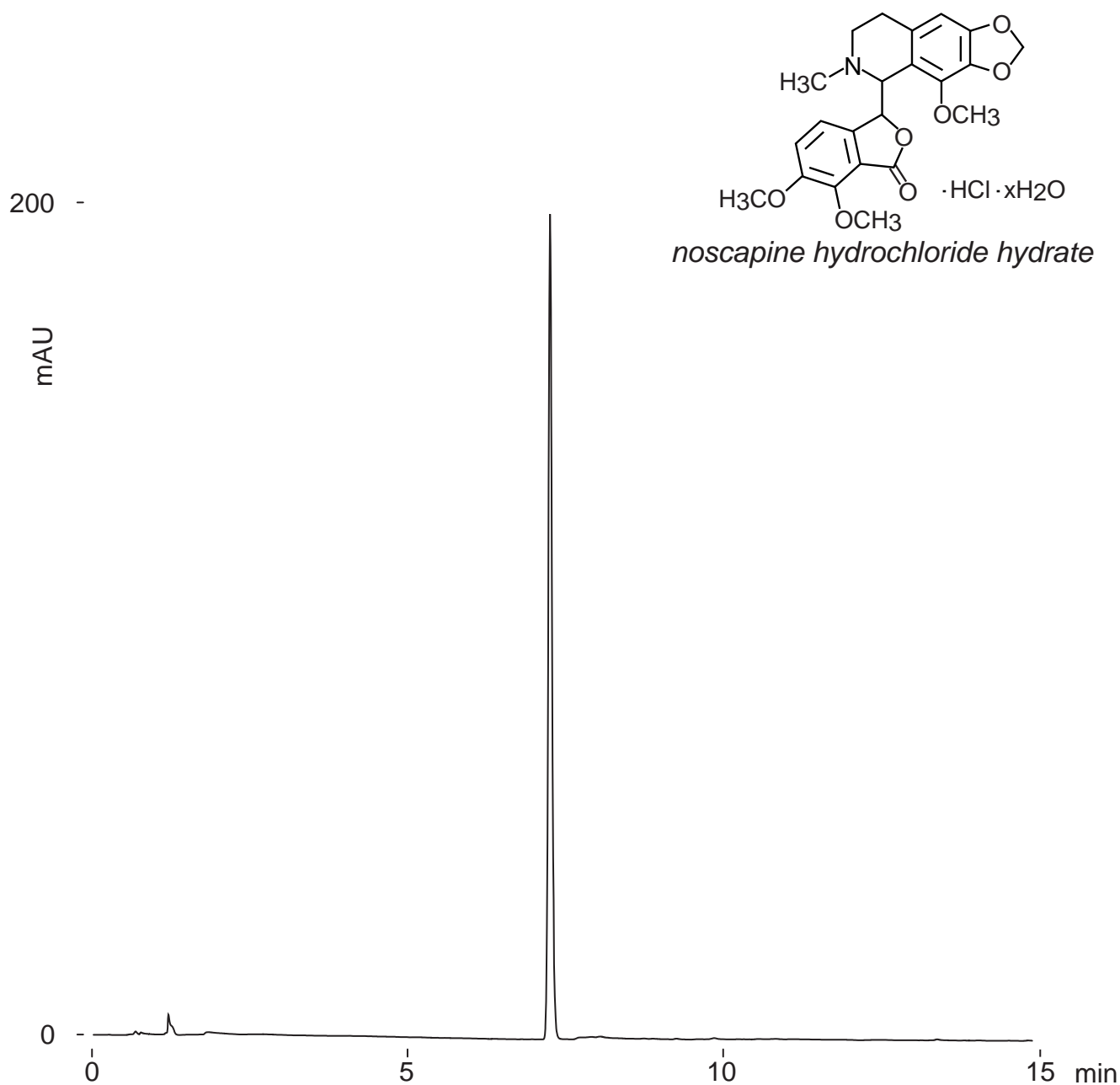

Dacapo DX-C18, 100 x 4.6 mm

A: acetonitrile / 10mM ammonium hydroxide= 10 / 90

B: acetonitrile / 10mM ammonium hydroxide = 90 / 10

0-100%B (0-15min)

1mL/min (11MPa), 37°C, 260nm

10µL (0.45mg/mL, acetonitrile / 1M ammonium hydroxide = 6/5)