

Unison UK-Phenyl Cadenza CX-C18



Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Khellin (Furanochromone, Vasodilator)

ケリン (フラノクロモン, 血管拡張薬)

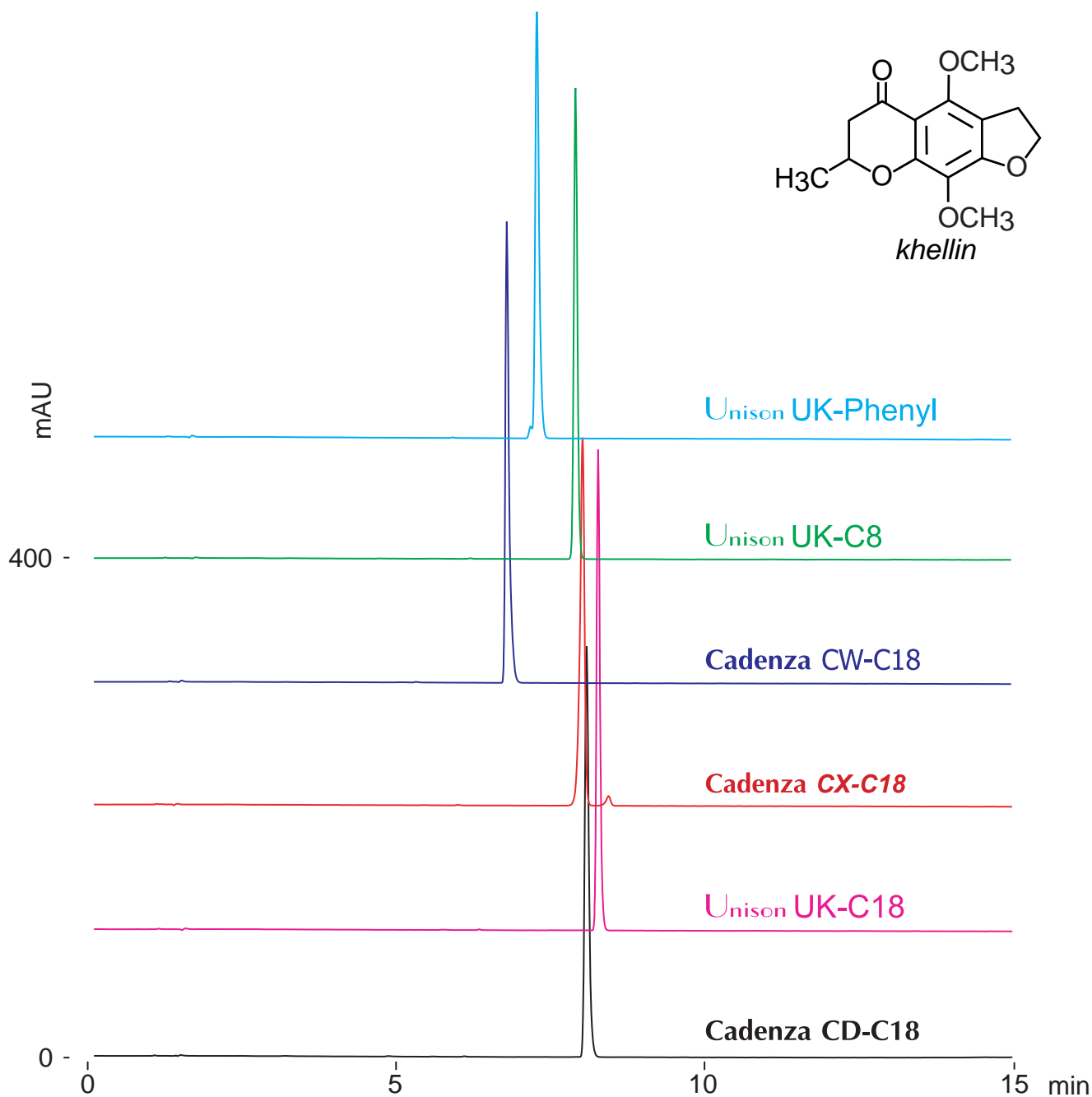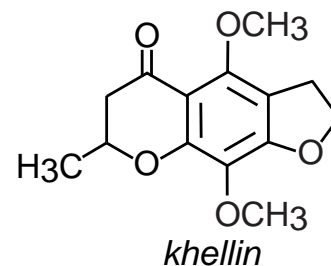

100 x 4.6 mm

A: 100mM ammonium formate / acetonitrile = 90 / 10

B: 100mM ammonium formate / acetonitrile = 10 / 90

0-90%B (0-15min)

1.0mL/min (11MPa), 37°C, 260nm

10µL (0.04mg/mL, 100mM ammonium formate / acetonitrile = 60 / 40)