

Unison UK-Phenyl

Unison UK-C8

Unison UK-C18





Cadenza CW-C18

Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Monocyclic pyrans (Koji fermentation metabolites)

単環式ピラン類 (麹発酵代謝副産物)

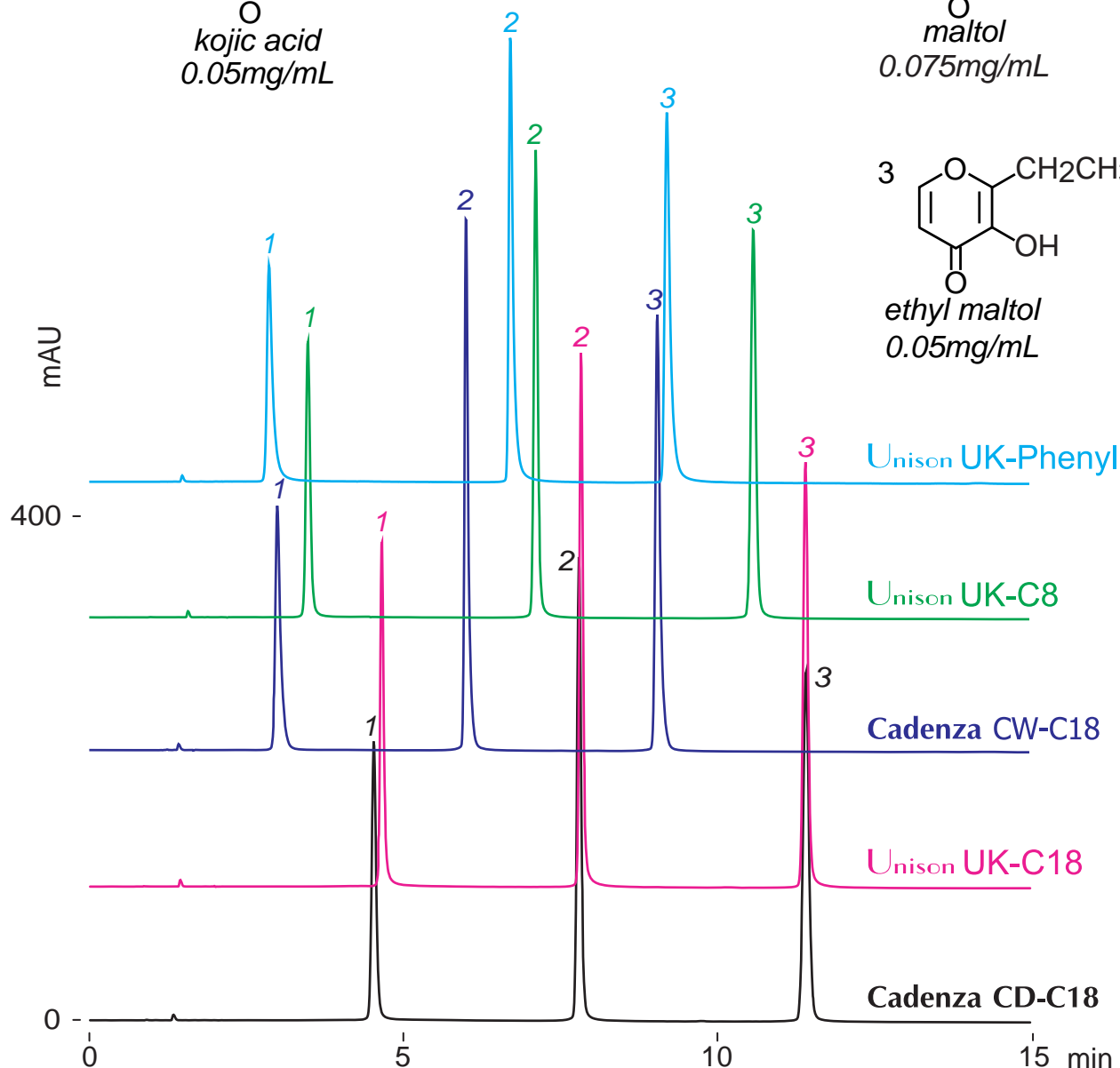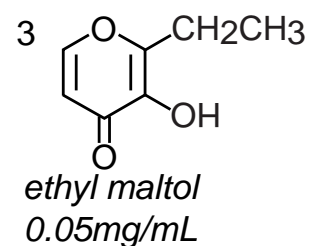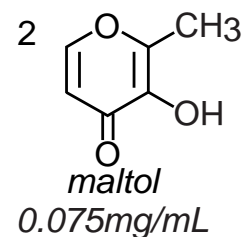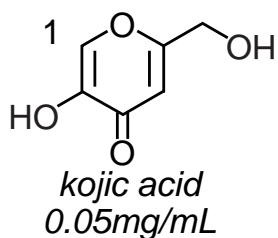

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 70 / 30 / 0.1

0%B (0-1min), 0-100%B (1-15min)

1.0mL/min (10MPa), 37°C, 260nm

10µL (0.05-0.075mg/mL, water / formic acid = 100 / 1)