

Unison UK-Phenyl Cadenza CX-C18


Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

## Hematoxylin (Histology Stains)

ヘマトキシリン (組織学用染色剤)

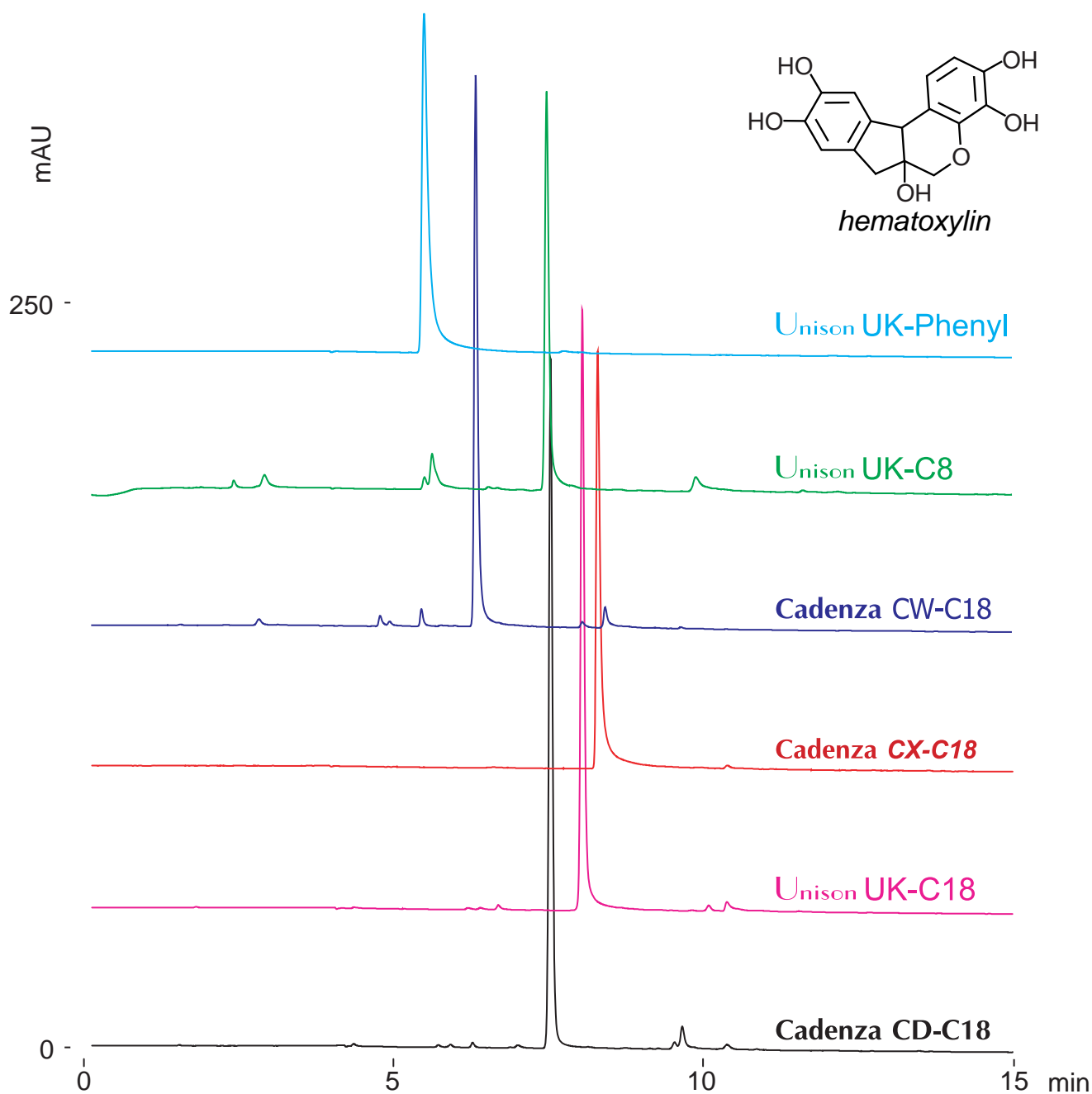

100 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 10 / 90

0%B (0-1min), 0-35%B (1-15min)

1.0mL/min (11MPa), 37°C, 290nm

10µL (0.25mg/mL, 100mM ammonium formate)