


Scherzo SM-C18 Nardis ND-RX
 Scherzo SW-C18 Scherzo SS-C18 100 x 4.6 mm


HPLC COLUMN APPLICATION

5-Formyl-2-furansulfonic acid (Furfural)

5-ホルミル-2-フランスルホン酸 (フルフラール類)

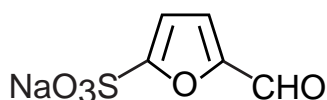

5-formyl-2-furansulfonic acid

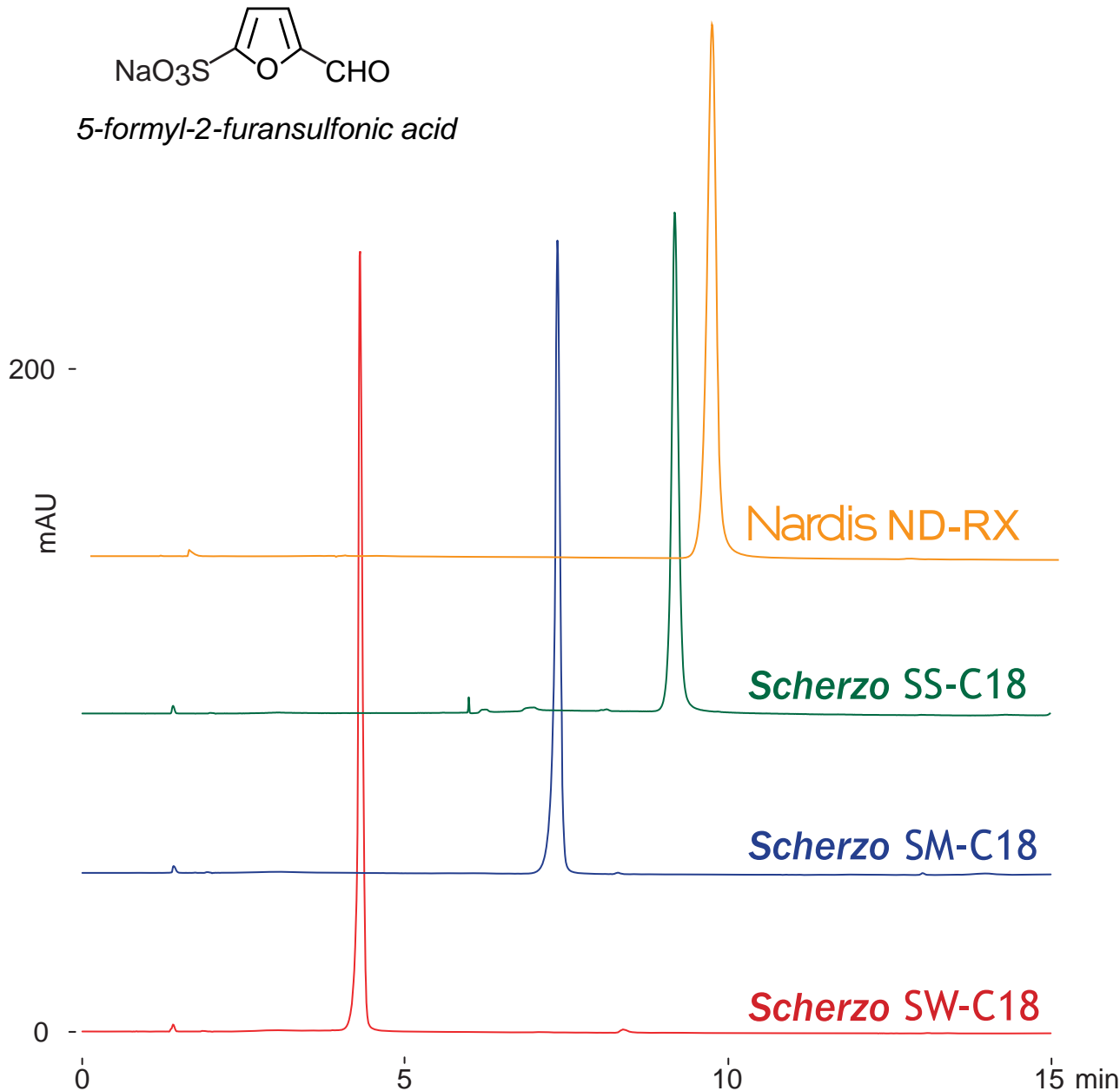

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: 100mM ammonium formate / acetonitrile = 50 / 50

0-100%B (0-15min)

1.0 mL/min (9MPa), 37°C, 280nm

10µL (0.05mg/mL, water / formic acid = 100 / 1)