
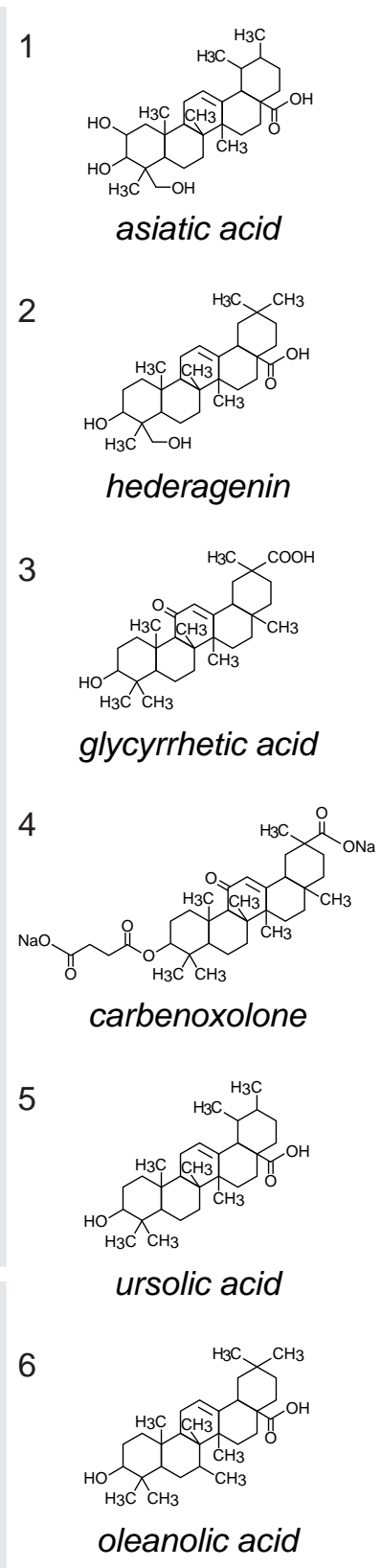


Cadenza CD-C18 Unison UK-C1
 Unison UK-C18 Unison UK-Phenyl
 Unison UK-C8

150 x 4.6 mm

HPLC COLUMN APPLICATION

Triterpenic acids
 トリテルペン酸

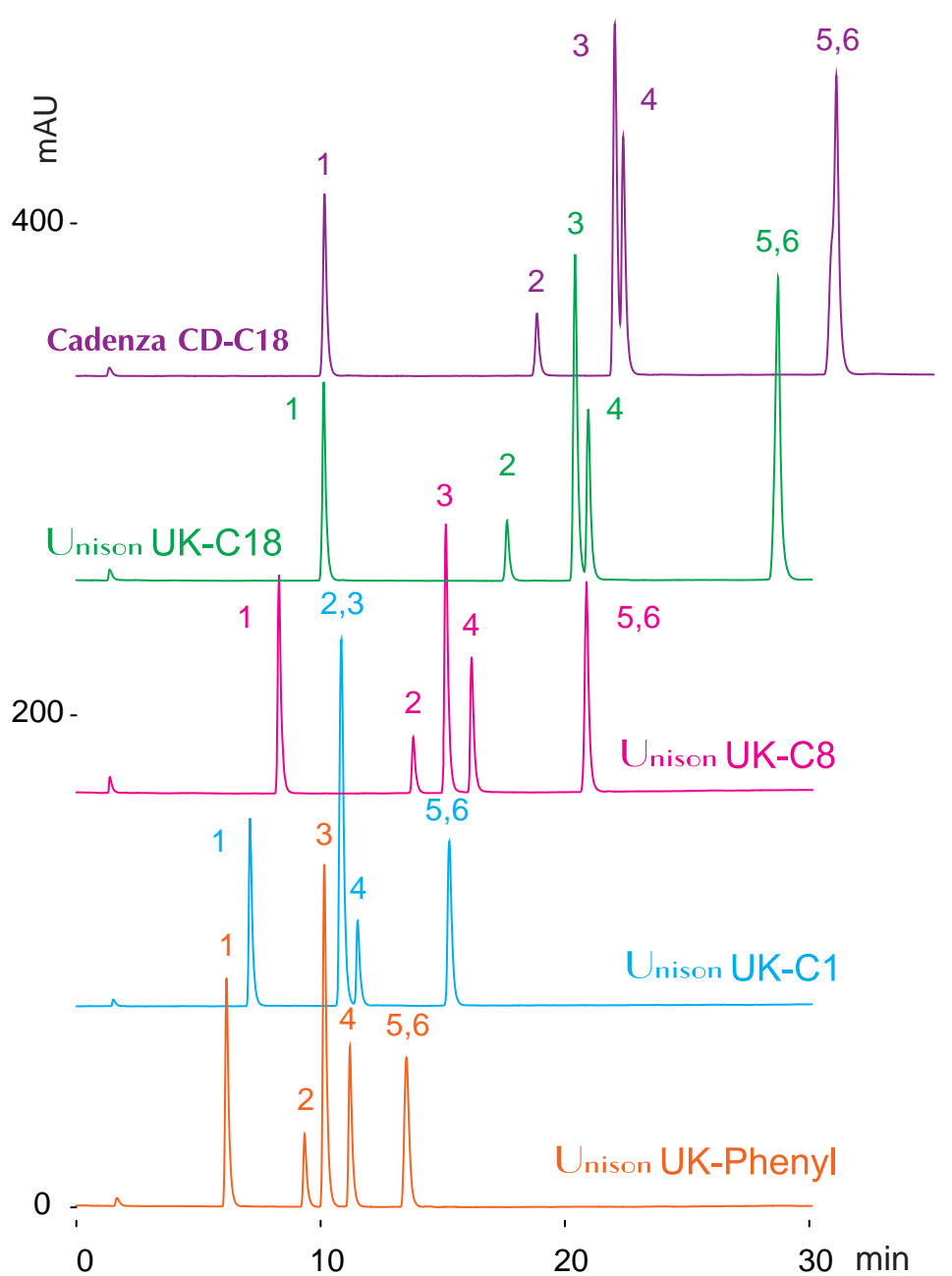

150 x 4.6 mm

A: water / formic acid = 100 / 0.5

B: acetonitrile / formic acid = 100 / 0.5

40-85%B (0-30min)

1.0mL/min (11-12MPa), 37°C

ELSD (spray chamber 70°C, drift chamber 100°C)

10µL (0.1-0.05mg/mL, DMSO/MeOH/EtOH=45/45/10)