

Unison UK-Phenyl Cadenza CX-C18


Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm


HPLC COLUMN APPLICATION

Dicoumarol (Coumarin anticoagulant)

ジクマロール (クマリン系抗凝固薬)

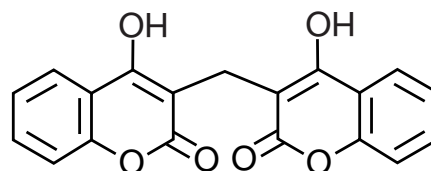

dicoumarol

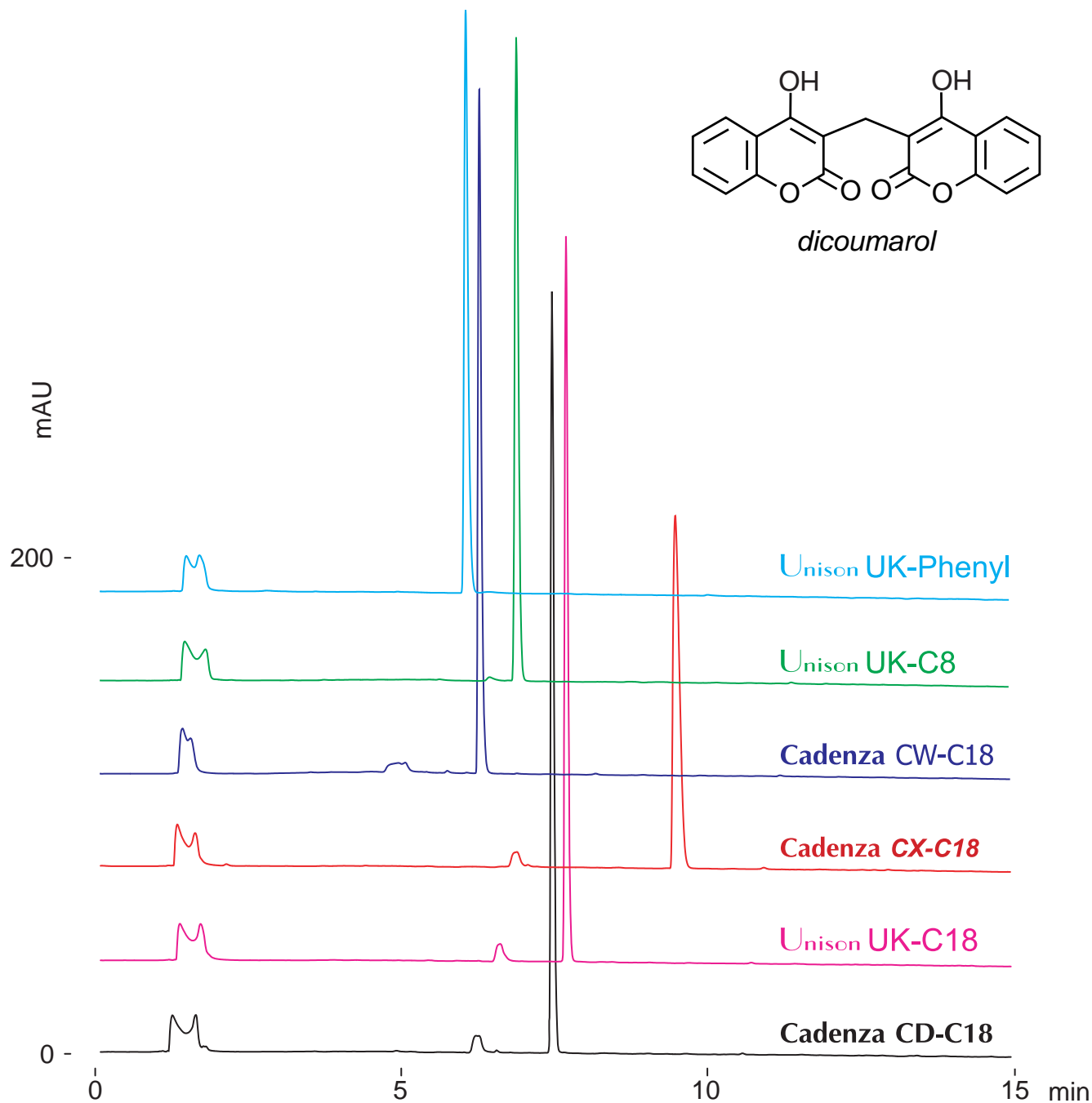

100 x 4.6 mm

A: 100mM ammonium formate / acetonitrile = 90 / 10

B: 100mM ammonium formate / acetonitrile = 10 / 90

0-100%B (0-15min)

1.0mL/min (10MPa), 37°C, 260nm

10µL (0.2mg/mL, 100mM ammonium formate / DMF = 50 / 50)