

Unison UK-Phenyl Cadenza CX-C18



Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Chrysazin (Contact laxatives)

クリサジン (接触刺激性下剤)

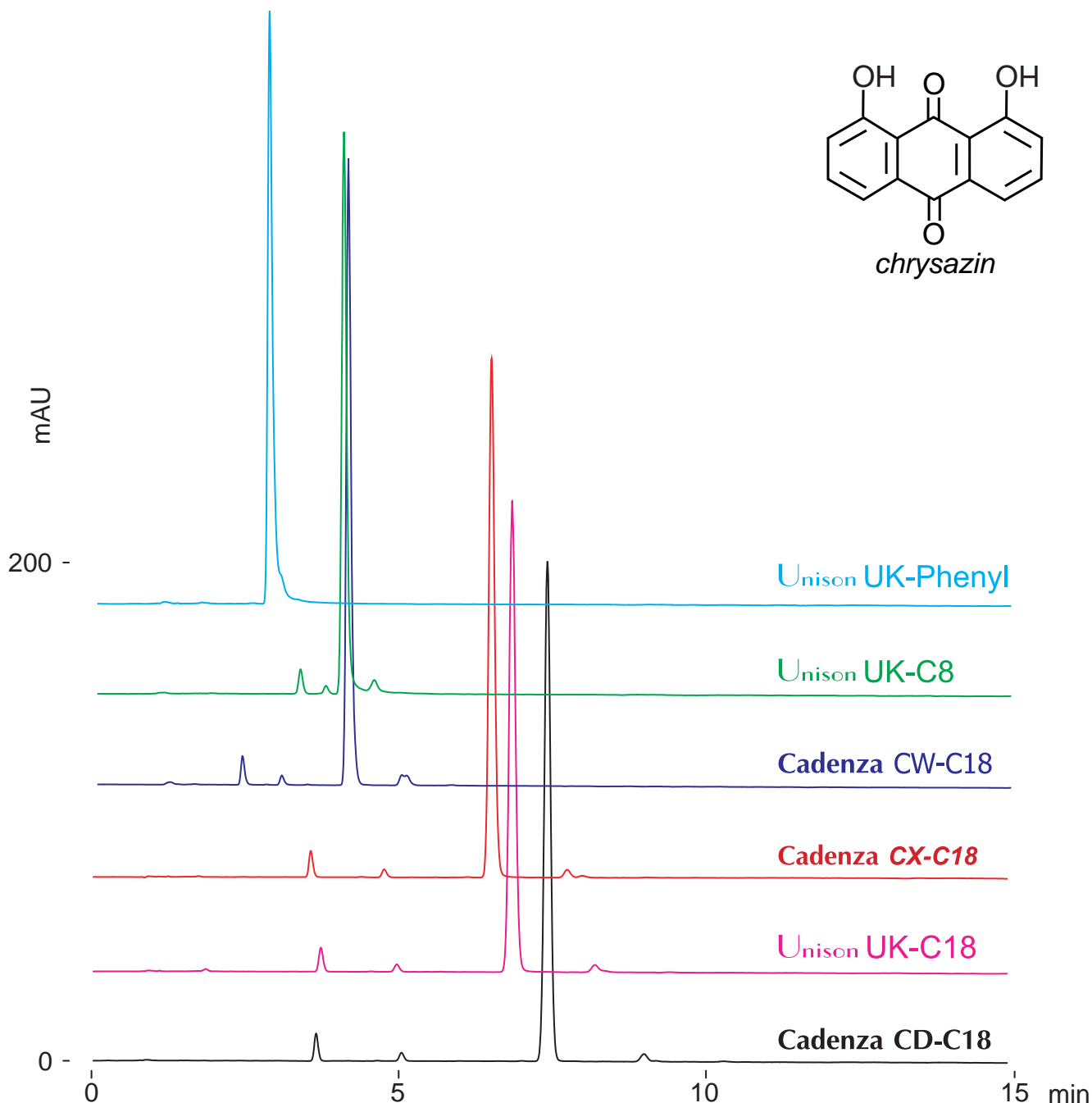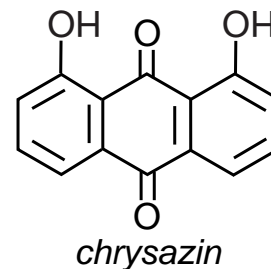

100 x 4.6 mm

A: 100mM ammonium formate / acetonitrile = 50 / 50

B: 100mM ammonium formate / acetonitrile = 10 / 90

10-100%B (0-15min)

1.0mL/min (10MPa), 37°C, 260nm

10µL (0.075mg/mL, 100mM ammonium formate / acetonitrile = 25 / 75)