


**Cadenza CD-C18**

150 x 3 mm

**HPLC COLUMN APPLICATION**

**LC-MS: Cholane derivatives (Cholanic acids)**

コラン誘導体 (コラン酸) の LC-MS 分析

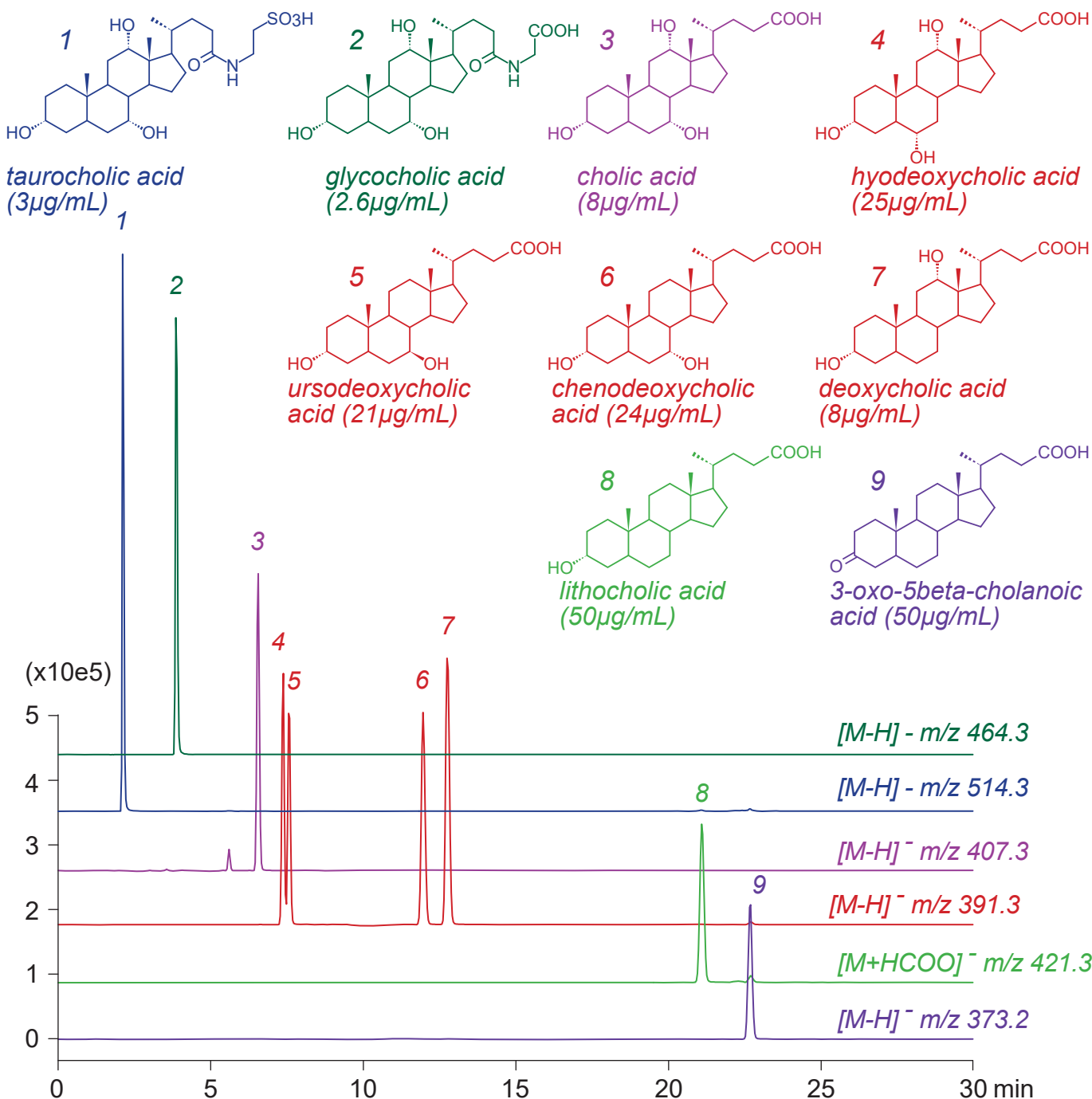

Cadenza CD-C18, 150 x 3 mm

A: acetonitrile / water / formic acid = 10 / 90 / 0.1

B: acetonitrile / water / formic acid = 90 / 10 / 0.1

40-100%B (0-30min)

0.4mL/min (9MPa), 37°C, Single Quad. MS (ESI)

5µL (2.6-50µg/mL, MeOH / water / HCOOH = 50 / 50 / 0.1)