

Unison UK-Phenyl Cadenza CX-C18


Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm


HPLC COLUMN APPLICATION

Aloe emodin (Stimulant laxative)

アロエエモジン (刺激性下剤)

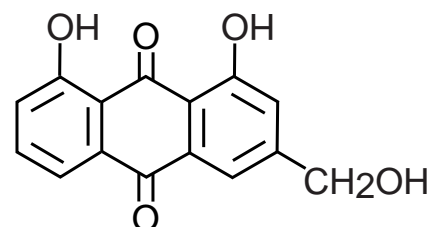

aloe emodin

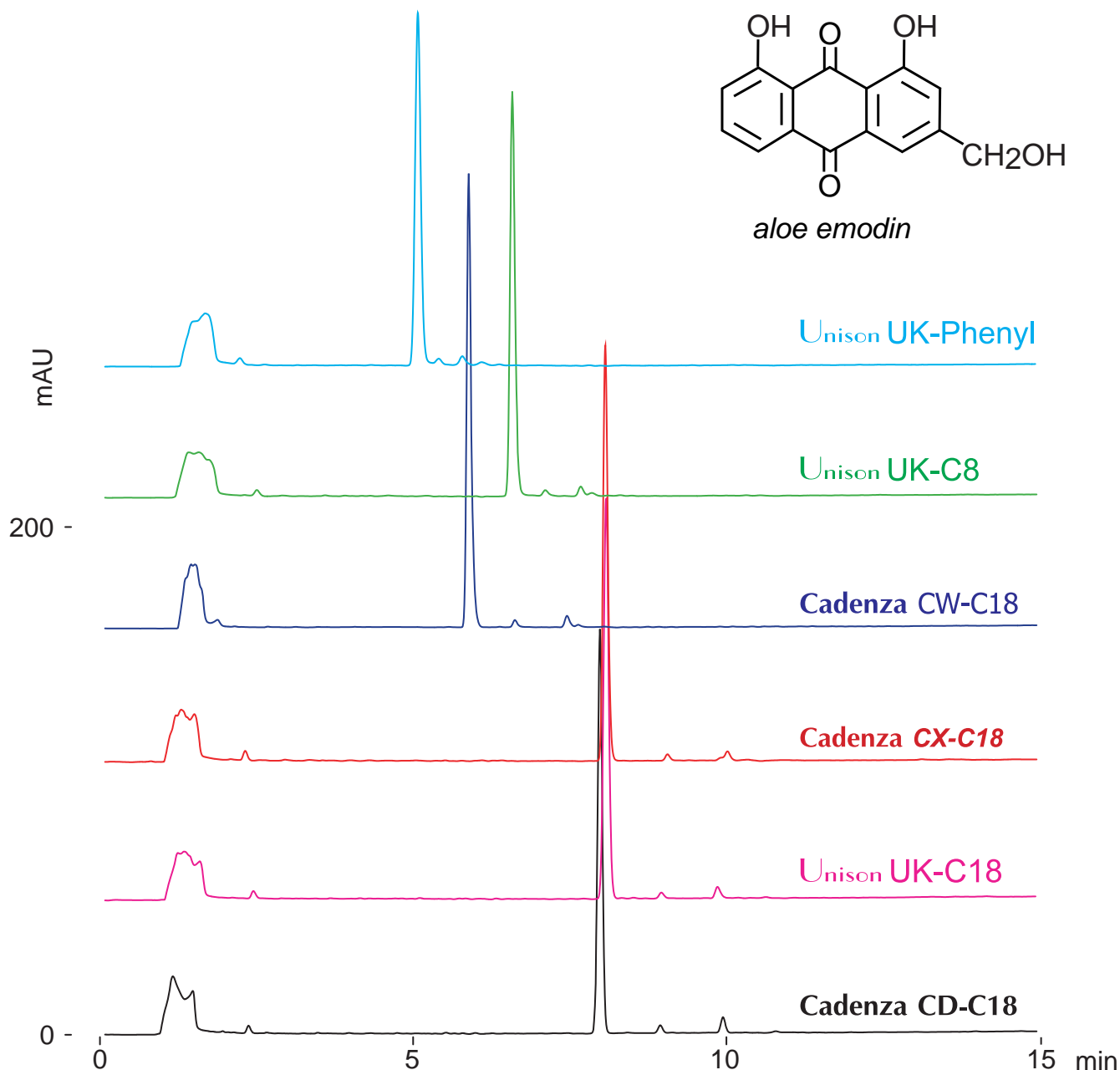

100 x 4.6 mm

A: water / acetonitrile / formic acid = 90 / 10 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

20-85%B (0-15min)

1.0mL/min (10MPa), 37°C, 260nm

10µL (0.05mg/mL, water / DMF / formic acid = 10 / 90 / 1)