



Unison UK-Phenyl Cadenza CX-C18

Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

**Brucine (Indole alkaloid)**  
 ブルシン (インドールアルカロイド)

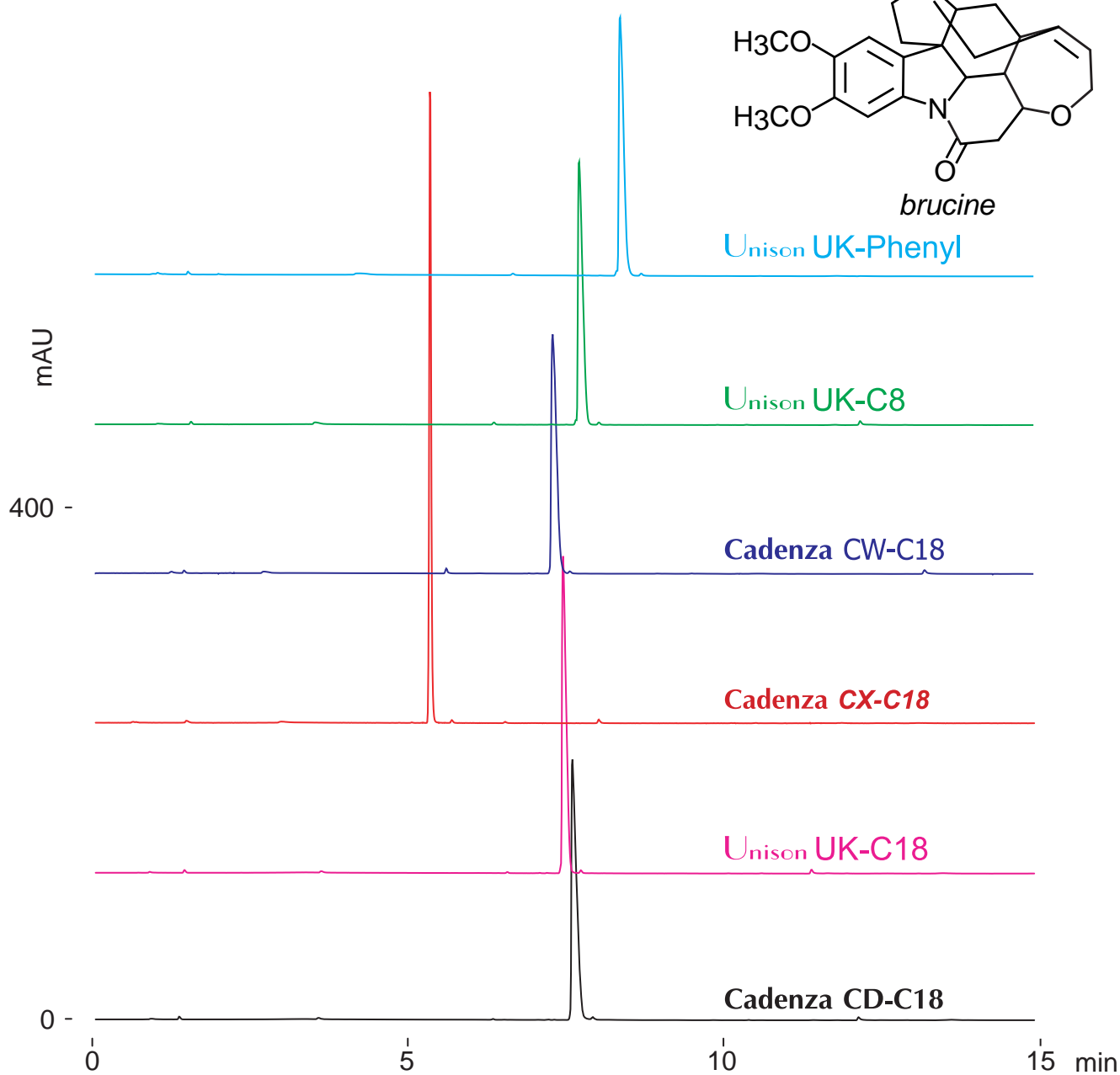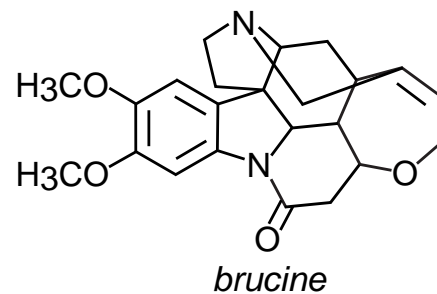

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

0-50%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.1mg/mL, water / formic acid = 100 / 1)