

Unison UK-Phenyl Cadenza CX-C18



Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

## Bergenin (Isocoumarin derivatives)

ベルゲニン (イソクマリン誘導体)

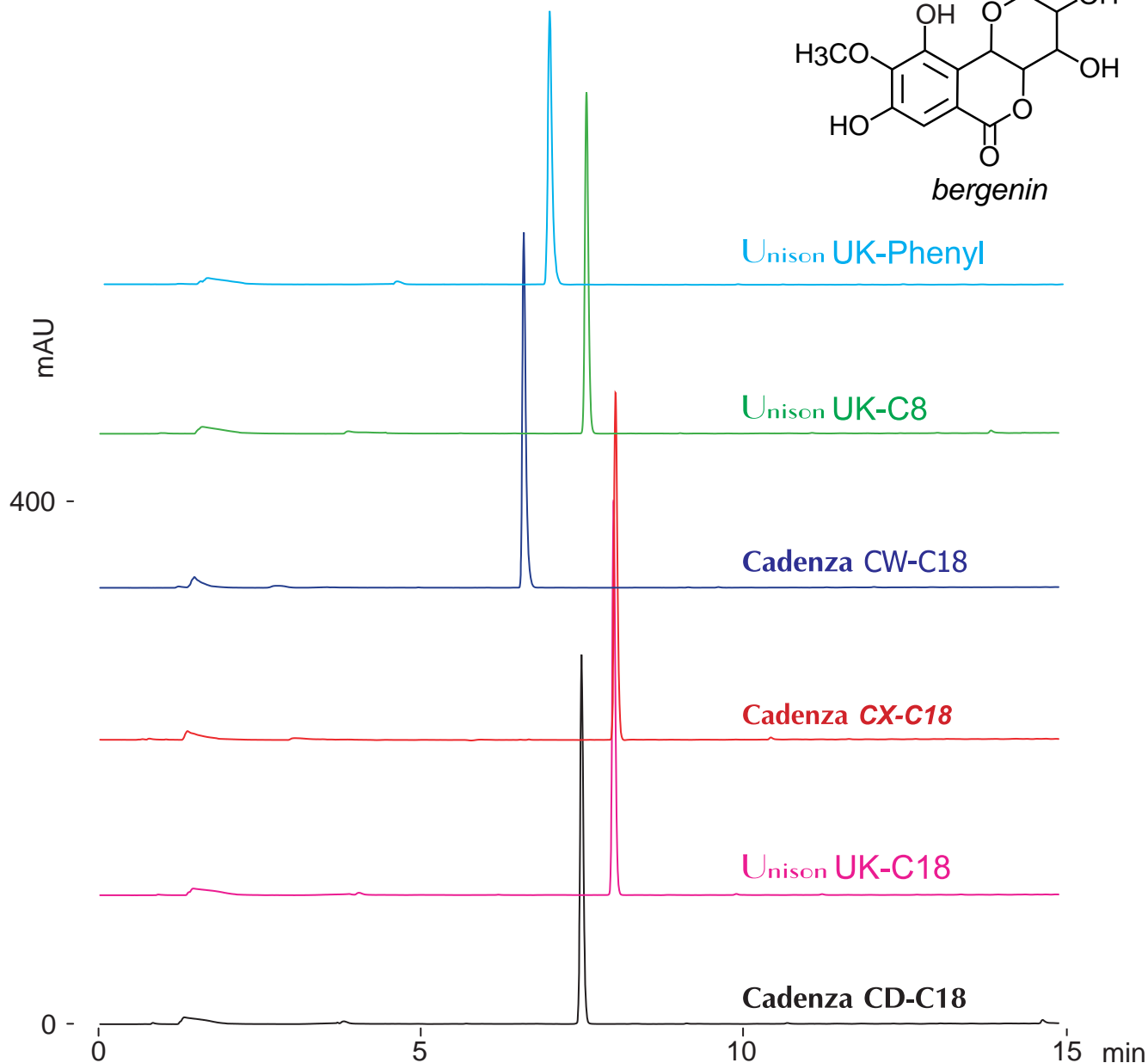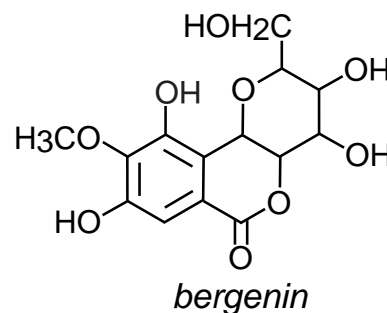

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

0%B (0-1min), 0-40%B (1-15min)

1.0mL/min (10MPa), 37°C, 270nm

10µL (0.15mg/mL, water / dimethyl sulfoxide / formic acid = 70 / 30 / 1)