


Scherzo SM-C18

150 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Aminoglycoside antibiotics

アミノグリコシド系抗生物質の LC-MS 分析

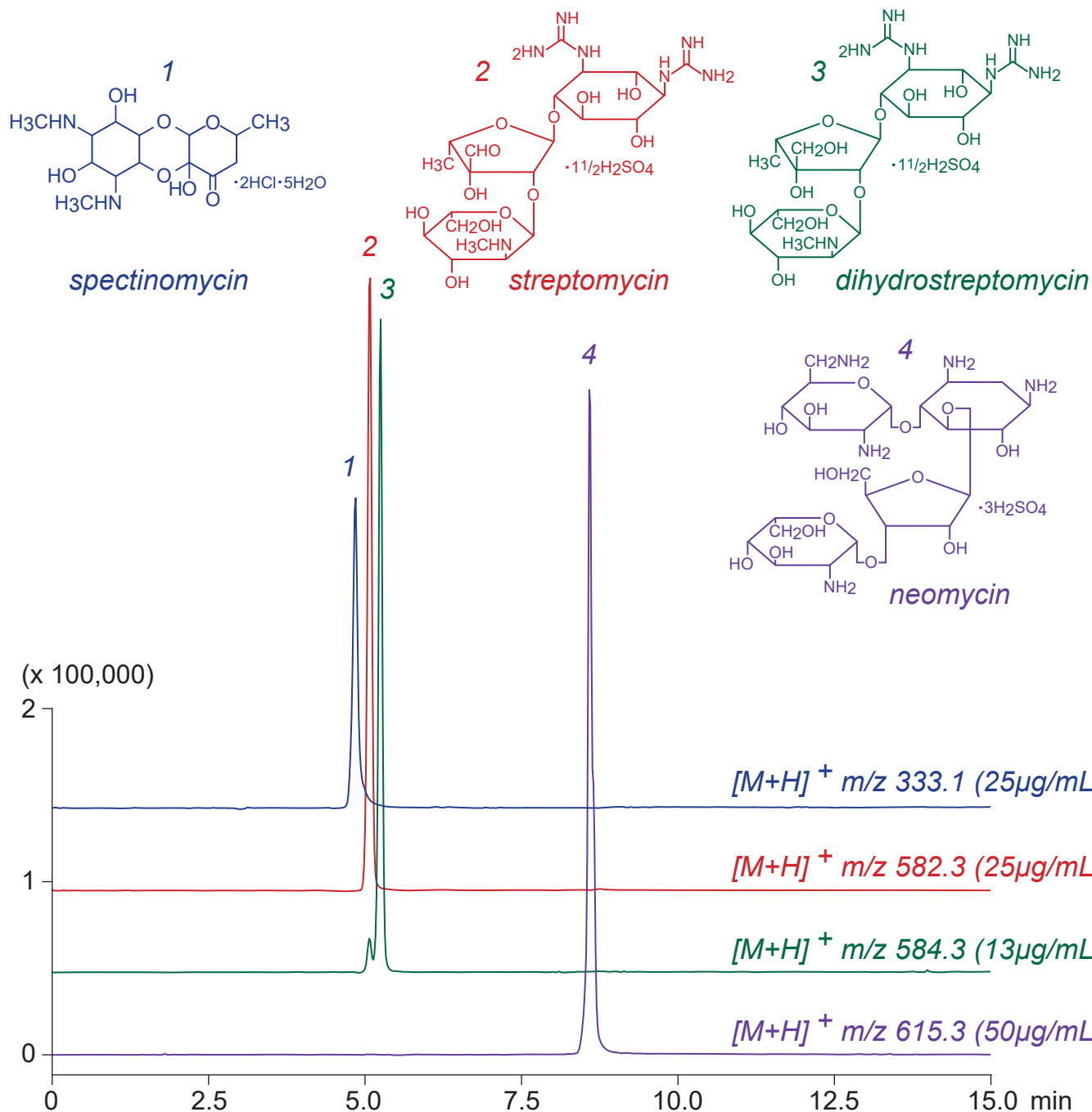

Scherzo SM-C18, 150 x 3 mm

A: 30mM ammonium acetate

B: acetonitrile / water / acetic acid = 20 / 80 / 0.5

5-50%B (0-6min), 90%B (6-15min)

0.4mL/min (9MPa), 37°C, Single Quad. MS (ESI)

10µL (13-50µg/mL, 100mM ammonium acetate)