






Nardis ND-RX

150 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Aminoglycoside antibiotics

アミノグリコシド系抗生物質の LC-MS 分析

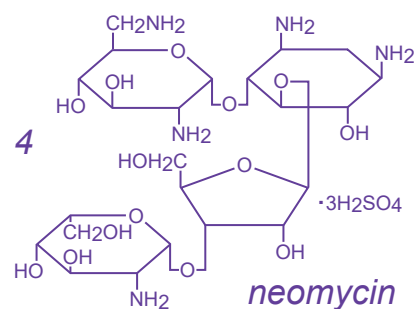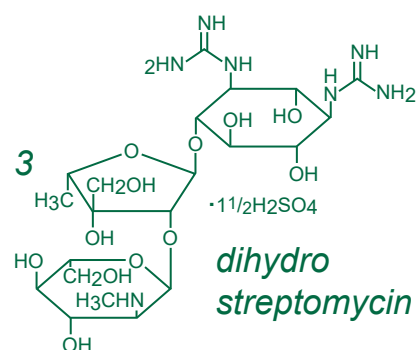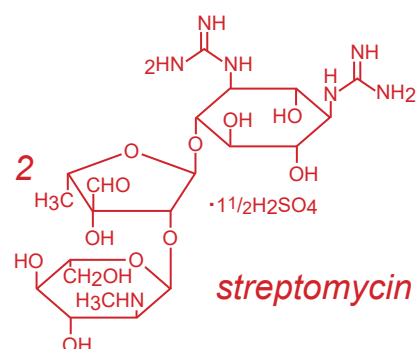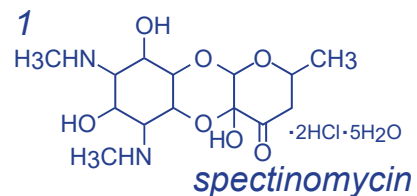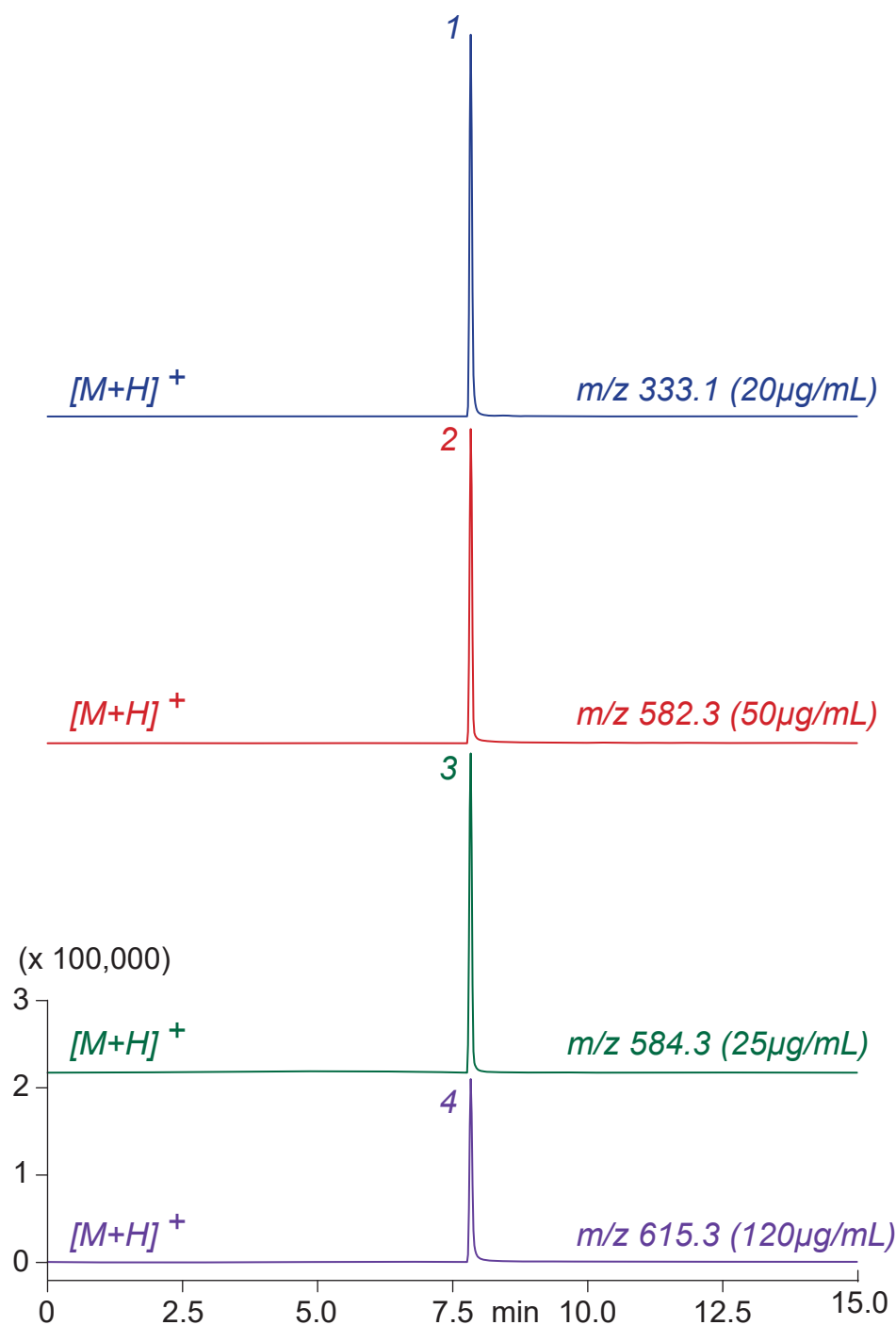

Nardis ND-RX, 150 x 3 mm
 A: 10mM ammonium formate
 B: acetonitrile / water / formic acid = 20 / 80 / 0.5
 0-50%B (0-15min)
 0.4mL/min (8MPa), 37°C, Single Quad. MS (ESI)
 10µL (20-120µg/mL, 100mM ammonium formate)