



Unison UK-Amino
 Nardis ND-NX
 Unison UK-Silica

150 x 4.6 mm

HPLC COLUMN APPLICATION

Arbutin (Phenolic glycosides)

アルブチン (フェノール性配糖体)

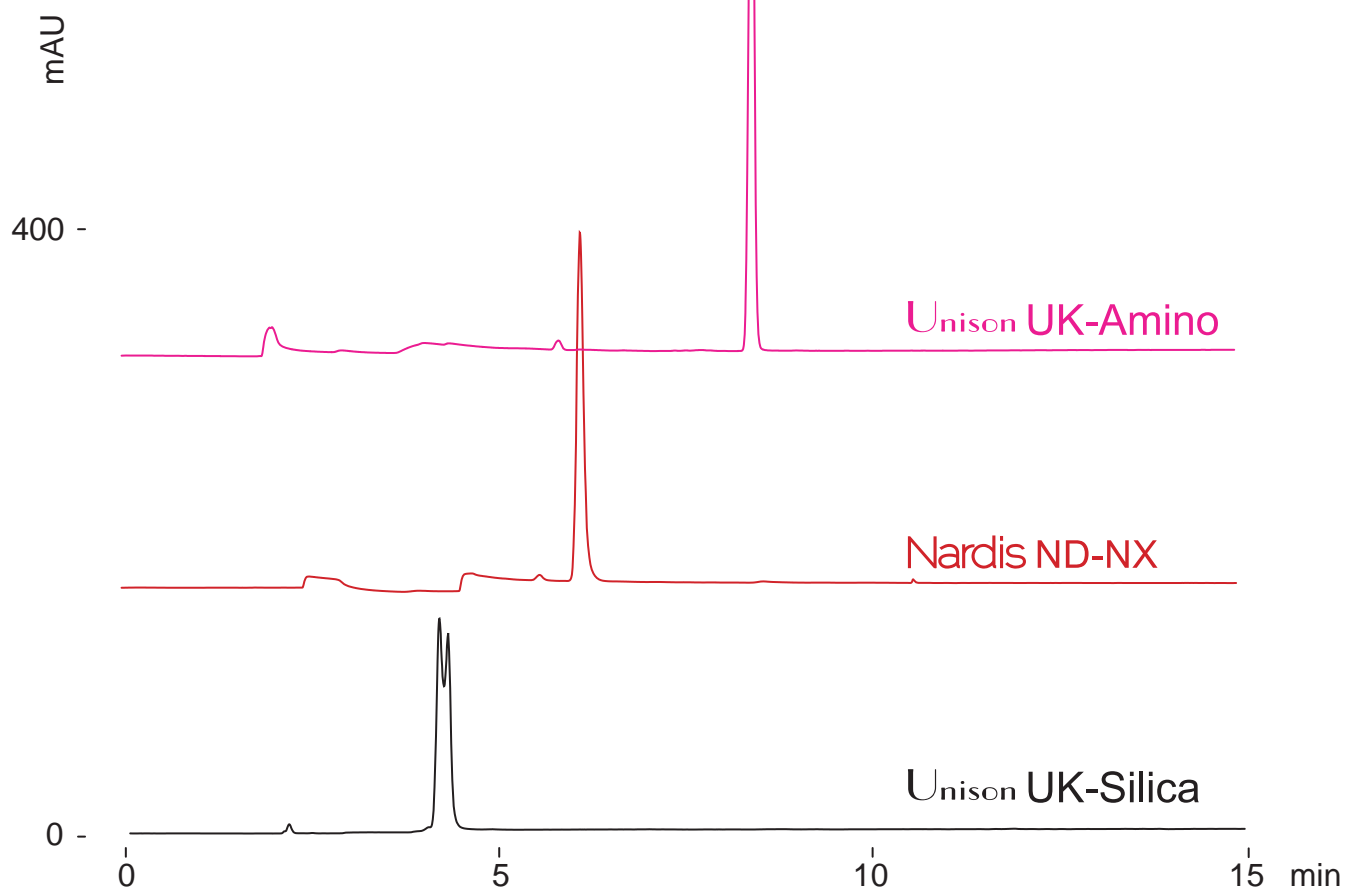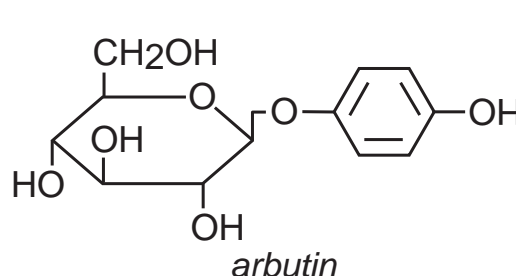

150 x 4.6 mm

A: acetonitrile / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 50 / 50 / 0.1

0-70%B (0-15min)

1.0mL/min (5MPa), 37°C, 280nm

10µL (1mg/mL, water / acetonitrile / formic acid = 30 / 70 / 1)