

Intrada Organic Acid



Nardis ND-NX MF

150 x 2 mm

HPLC COLUMN APPLICATION

LC-MS: Beta-nicotinamide mononucleotide

- ニコチンアミドモノヌクレオチドの LC-MS 分析

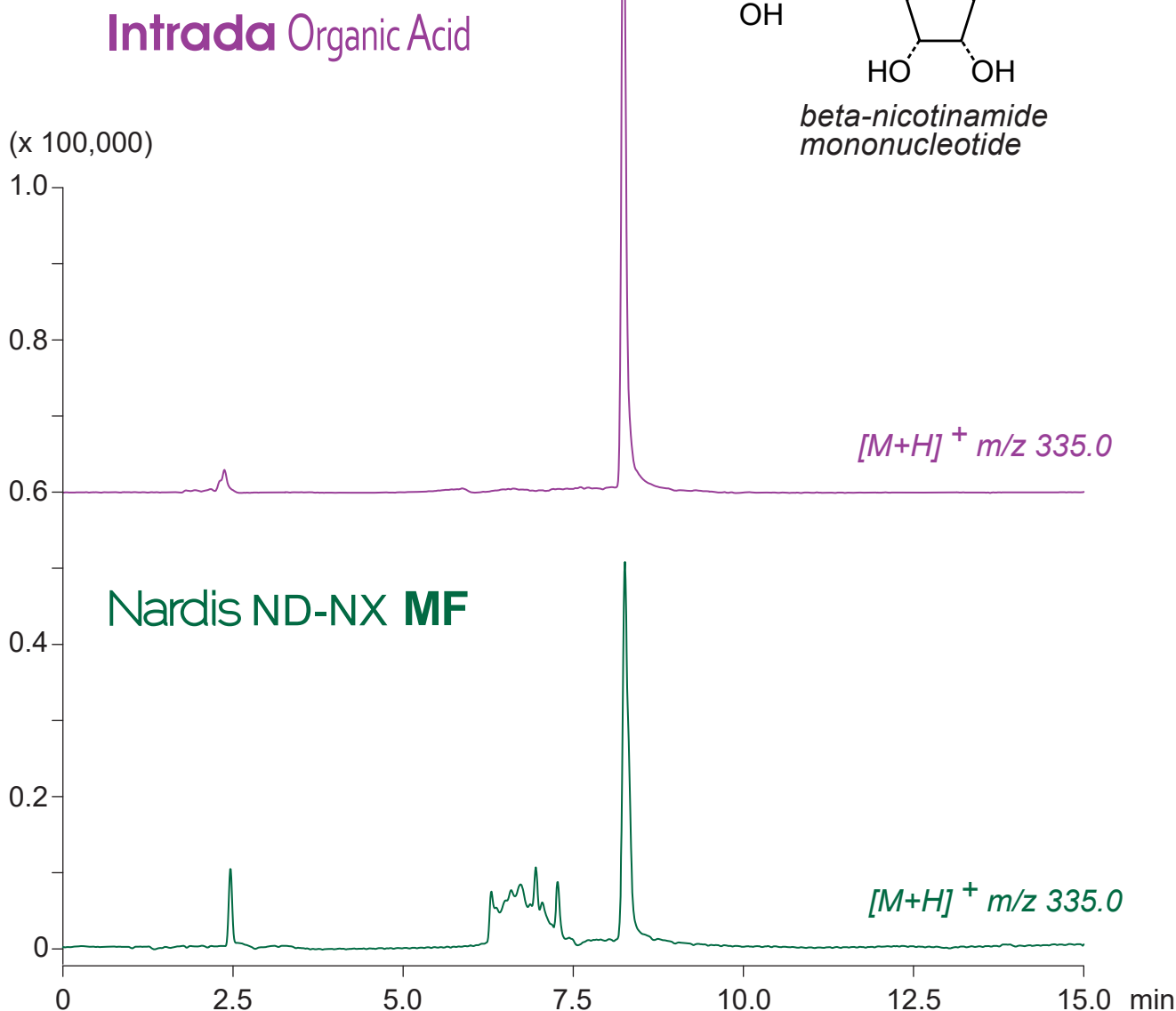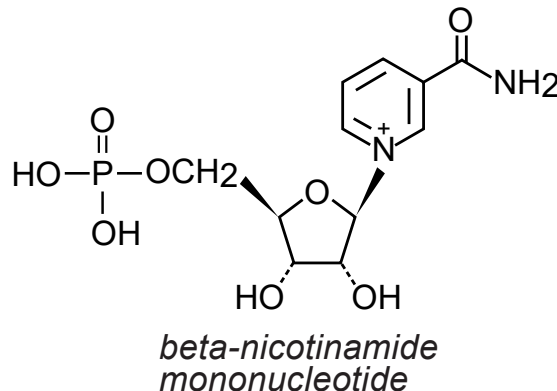

150 x 2 mm

A: acetonitrile / 100mM ammonium formate = 95 / 5

B: acetonitrile / 100mM ammonium formate = 10 / 90

0%B (0-1min), 0-100%B (1-7min), 100%B (7-15min)

0.2mL/min (4-5MPa), 37°C, Single Quad. MS (ESI)

5μL (10μg/mL, acetonitrile / water / formic acid = 50 / 50 / 1)