


**Nardis ND-NX MF**

Unison UK-Amino MF 150 x 2 mm

**HPLC COLUMN APPLICATION**

**Beta-nicotinamide mononucleotide(NMN)**

- ニコチンアミドモノヌクレオチド (NMN)の分析

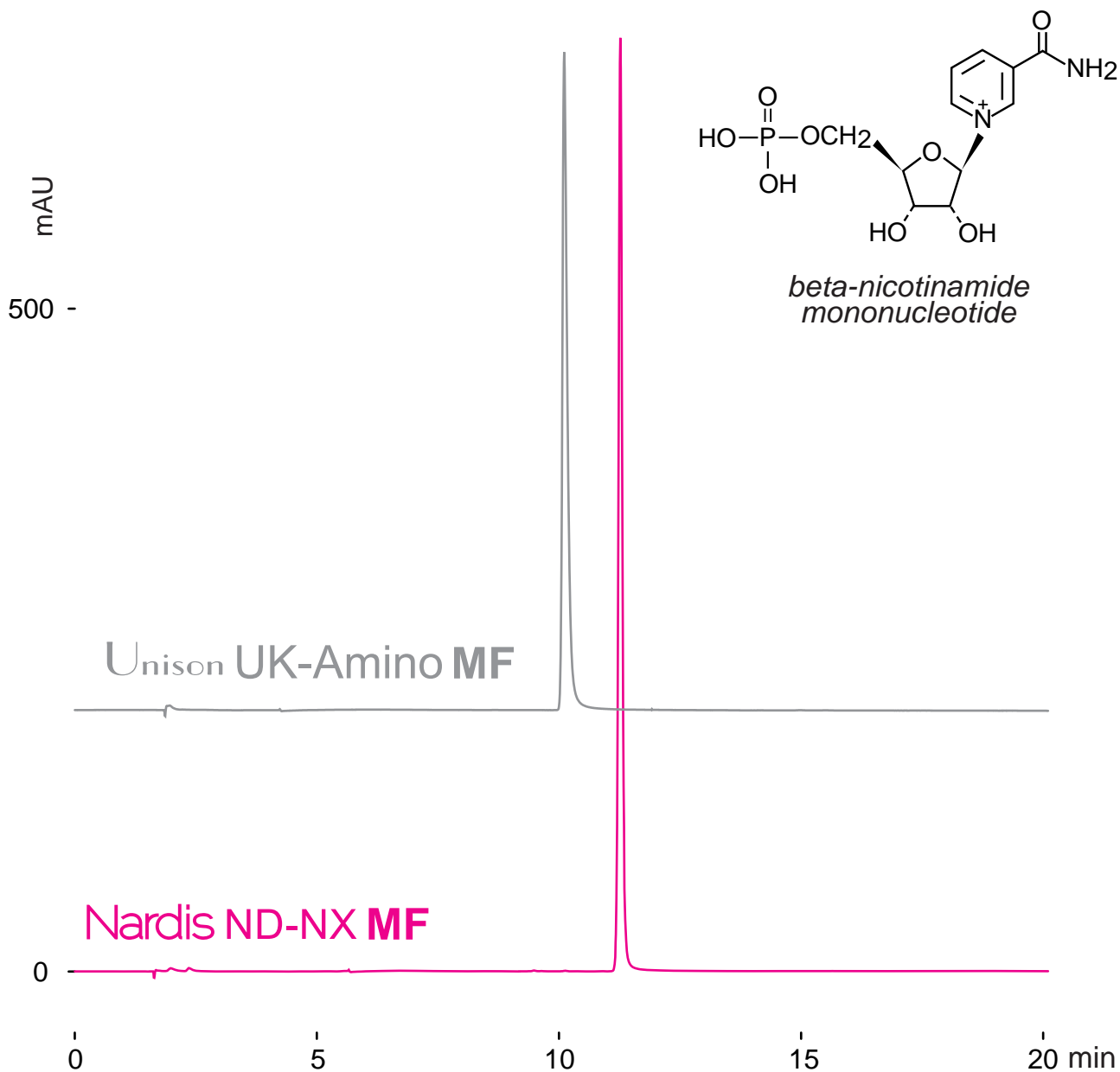

150 x 2 mm  
 A: acetonitrile  
 B: acetonitrile / 100mM ammonium acetate = 25 / 75  
 25 -100%B (0 -15 min), 100%B (15 -20 min)  
 0.2 mL/min (6-7MPa), 37°C, 260nm  
 3 µL (1.0mg/mL, water)