



Unison UK-C1
 Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Usnic Acid (Aromatic polyketide)

ウスニン酸 (芳香族ポリケチド)

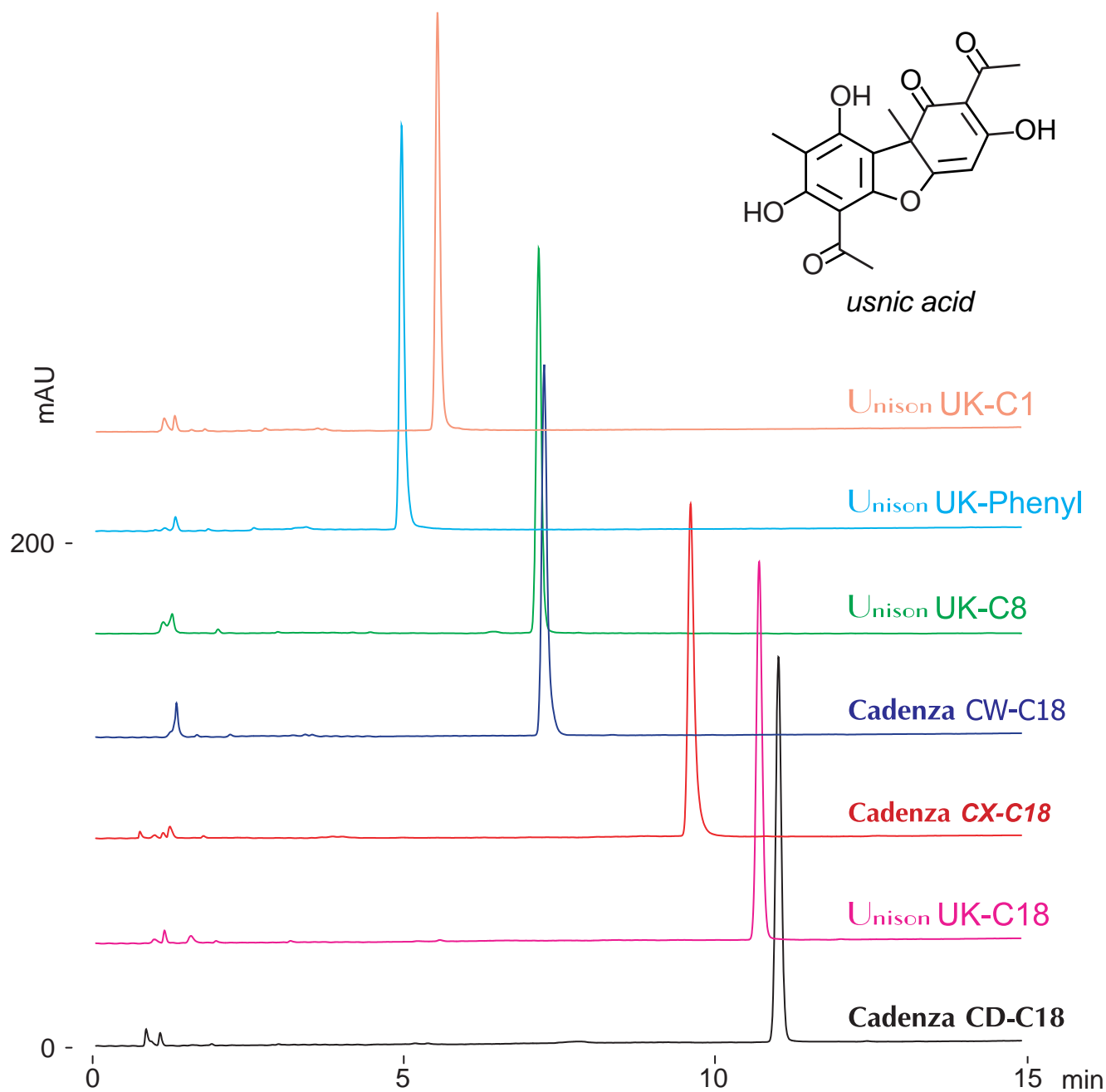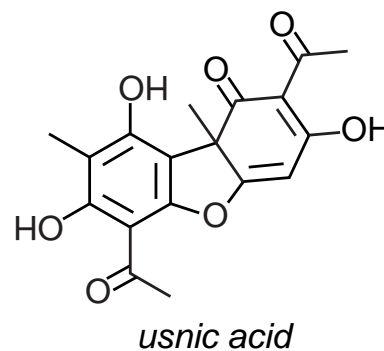

100 x 4.6 mm

A: water / acetonitrile / formic acid = 90 / 10 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

50-100%B (0-15min)

1.0mL/min (7MPa), 37°C, 260nm

10µL (0.1mg/mL, water / acetonitrile / formic acid = 25 / 75 / 1)