


Intrada Organic Acid

150 x 2 mm

HPLC COLUMN APPLICATION

LC-MS: Dioctyl sulfosuccinate (Double-tail anionic surfactant)

ジオクチルスルホコハク酸 (二鎖型アニオン界面活性剤) の LC-MS 分析

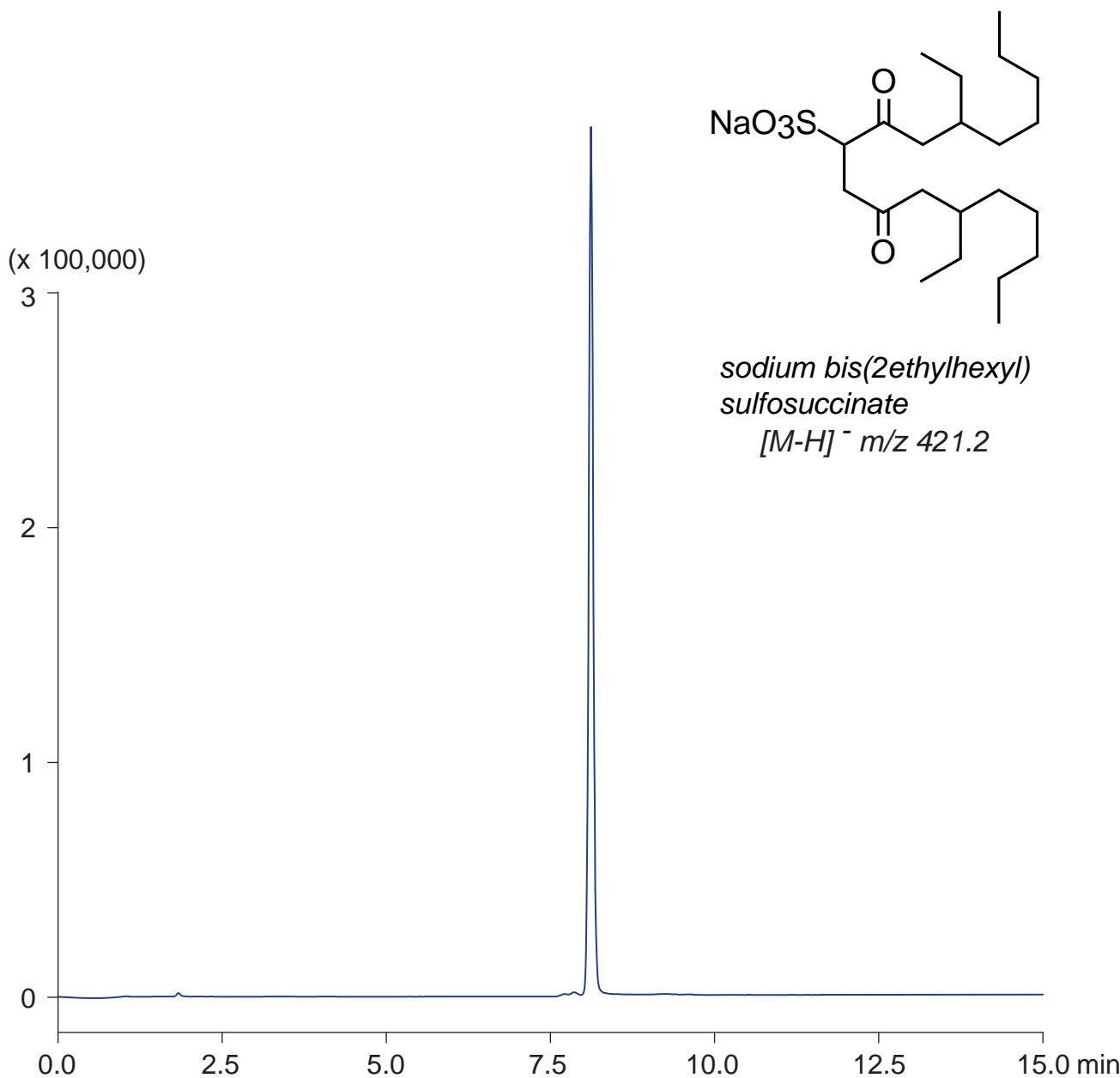

Intrada Organic Acid, 150 x 2 mm

A: acetonitrile / water / formic acid = 30 / 70 / 0.1

B: acetonitrile / 100mM ammonium formate = 60 / 40

0-100%B (0-10min), 100%B (10-15min)

0.2mL/min (13MPa), 37°C, Single Quad. MS (ESI)

5uL (50ng/mL, ACN / water / HCOOH = 50 / 50 / 1)