




Unison UK-C8

100 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Avilamycin

アビラマイシンの LC-MS/MS 分析

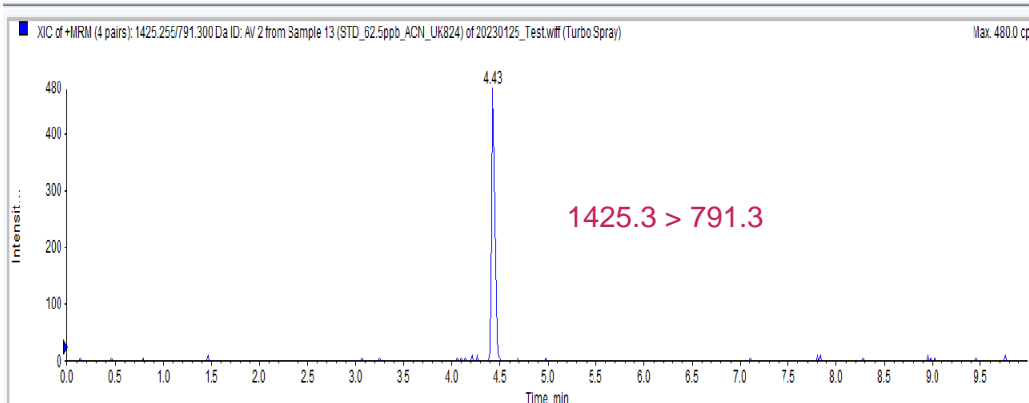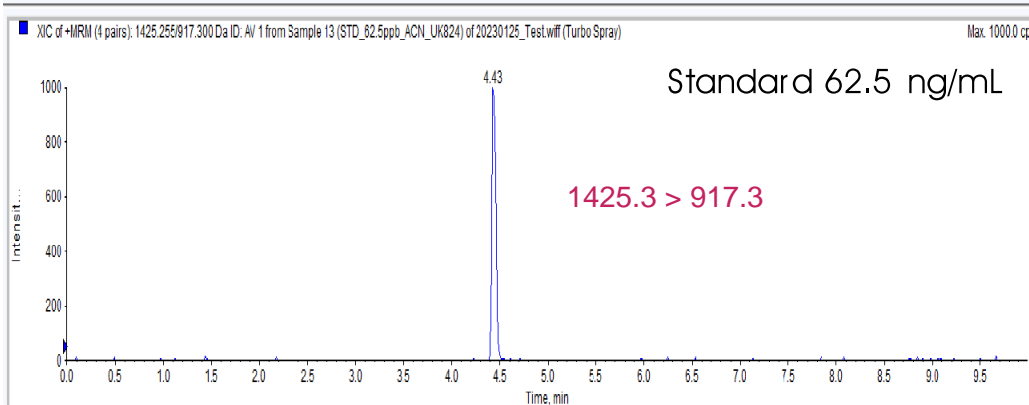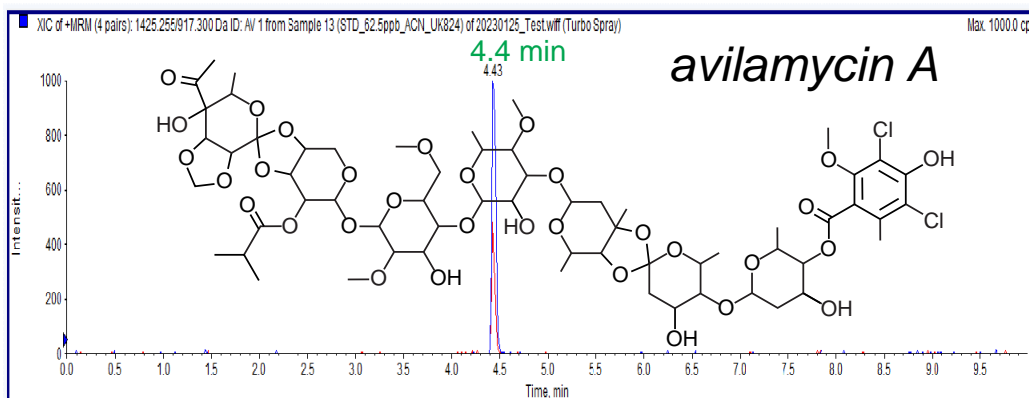

Unison UK-C8, 100 x 2 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

10%B (0-1min), 10-90%B (1-3min), 100%B (3-7min),

100-10%B (7-7.1min), 10%B (7.1-10min)

0.3mL/min, 40°C, Triple Quad. MS (ESI, Positive)

5µL (62.5 ng/mL)

Courtesy of PEAKMAN SP, KOREA