


Unison UK-C18

100 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Chlorogenic acid

クロロゲン酸の LC-MS/MS 分析

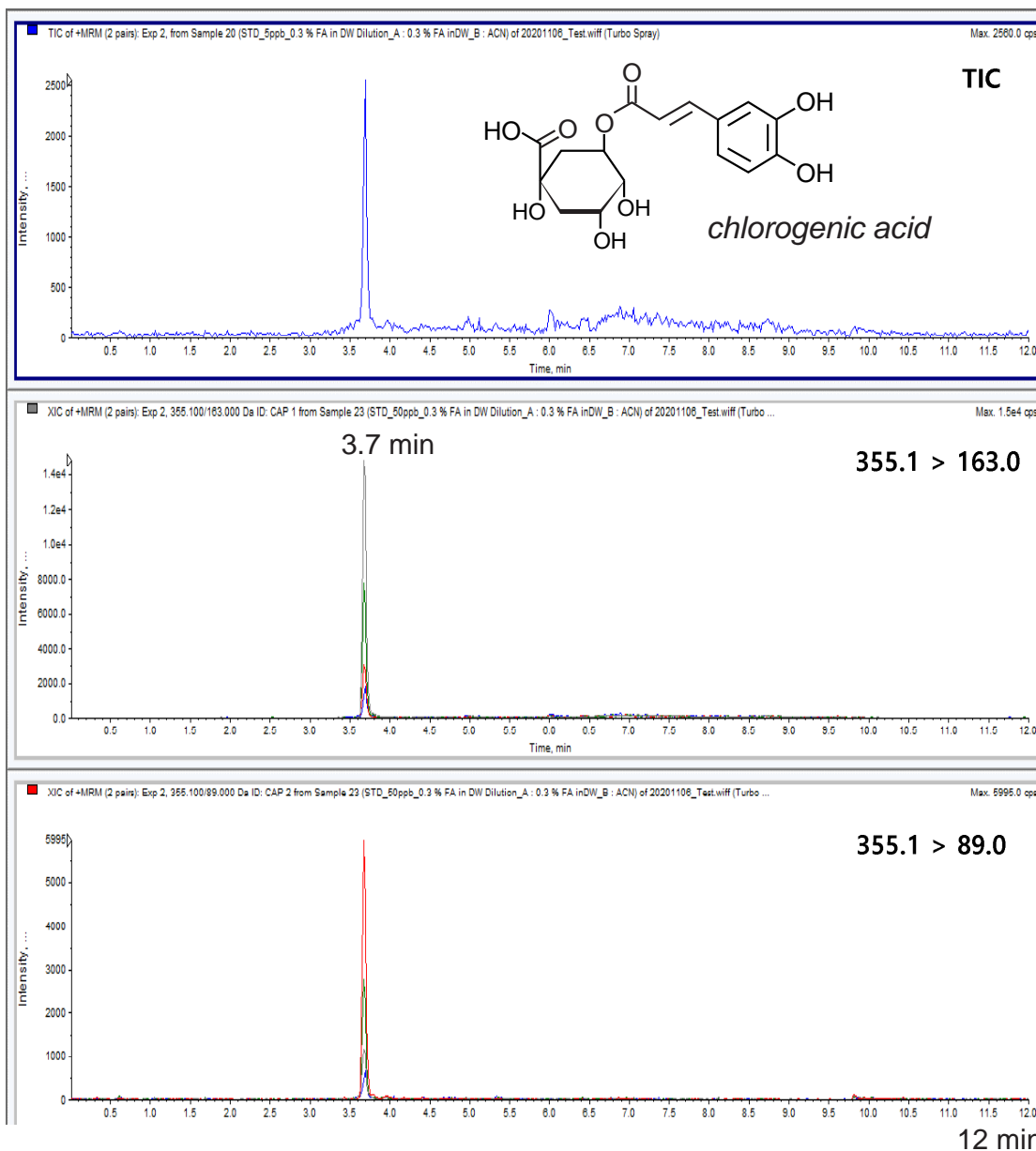

Unison UK-C18, 100 x 2 mm
 A: 0.3% formic acid in water
 B: acetonitrile
 5%B (0-1min), 5-100%B (1-5min), 100%B (5-9min),
 100-5%B (9-9.1min), 5%B (9.1-12min)
 0.3mL/min, 40°C, Triple Quad. MS (ESI, Positive)
 2µL (5, 10, 20, 50ppb, 0.3% formic acid in water)

Courtesy of PEAKMAN SP., KOREA