



Scherzo SM-C18 MF

100 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Glufosinate metabolites

グルホシネート代謝物の LC-MS/MS 分析

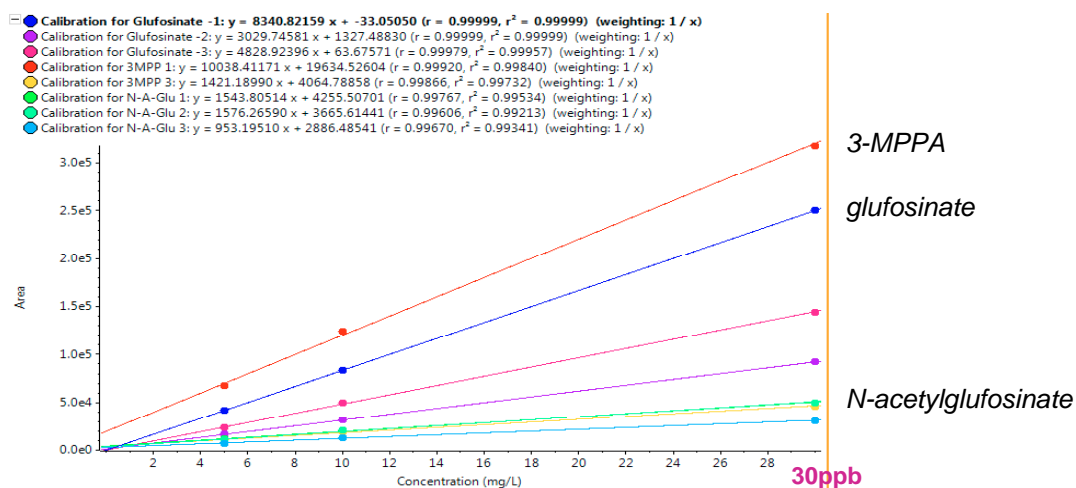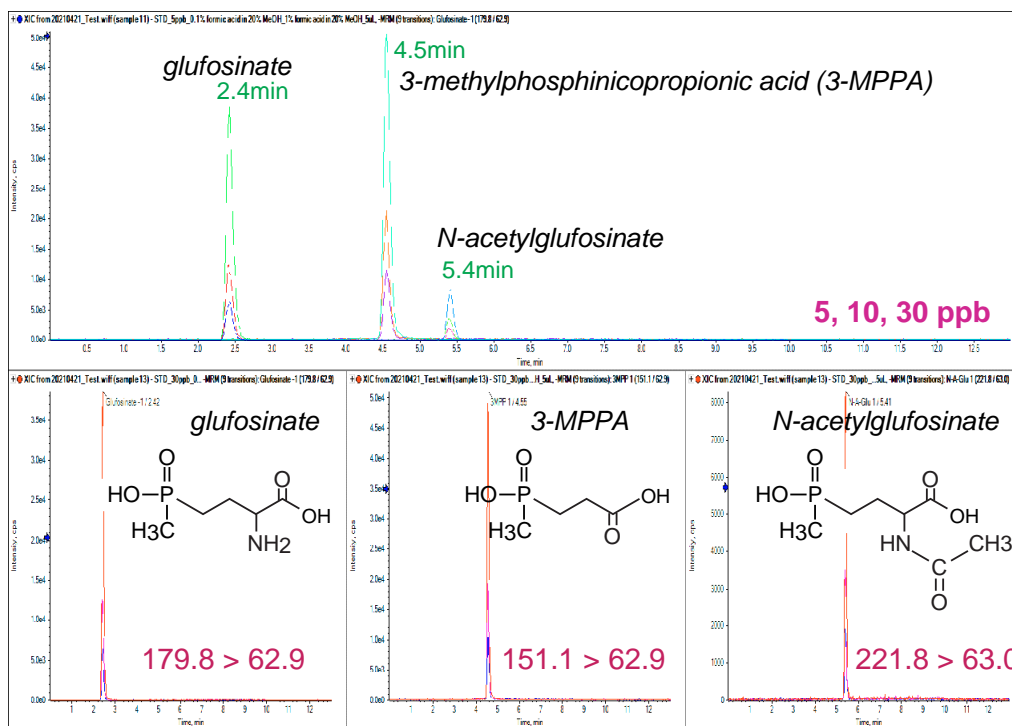

Scherzo SM-C18 MF, 100 x 2 mm
A: 0.1% formic acid in 20% methanol
B: 1% formic acid in 20% methanol
 0%B (0-1min), 0-100%B (1-4min), 100%B (4-7min),
 100-0%B (7-7.1min), 0%B (7.1-13min)
 0.2mL/min, 40°C, Triple Quad. MS (ESI, Negative)
 5μL (5, 10, 30 ppb)

Courtesy of PEAKMAN SP, KOREA