


Intrada Amino Acid 100 x 3 mm **HPLC COLUMN APPLICATION**

LC-MS/MS: Amino acids
アミノ酸の LC-MS/MS 分析

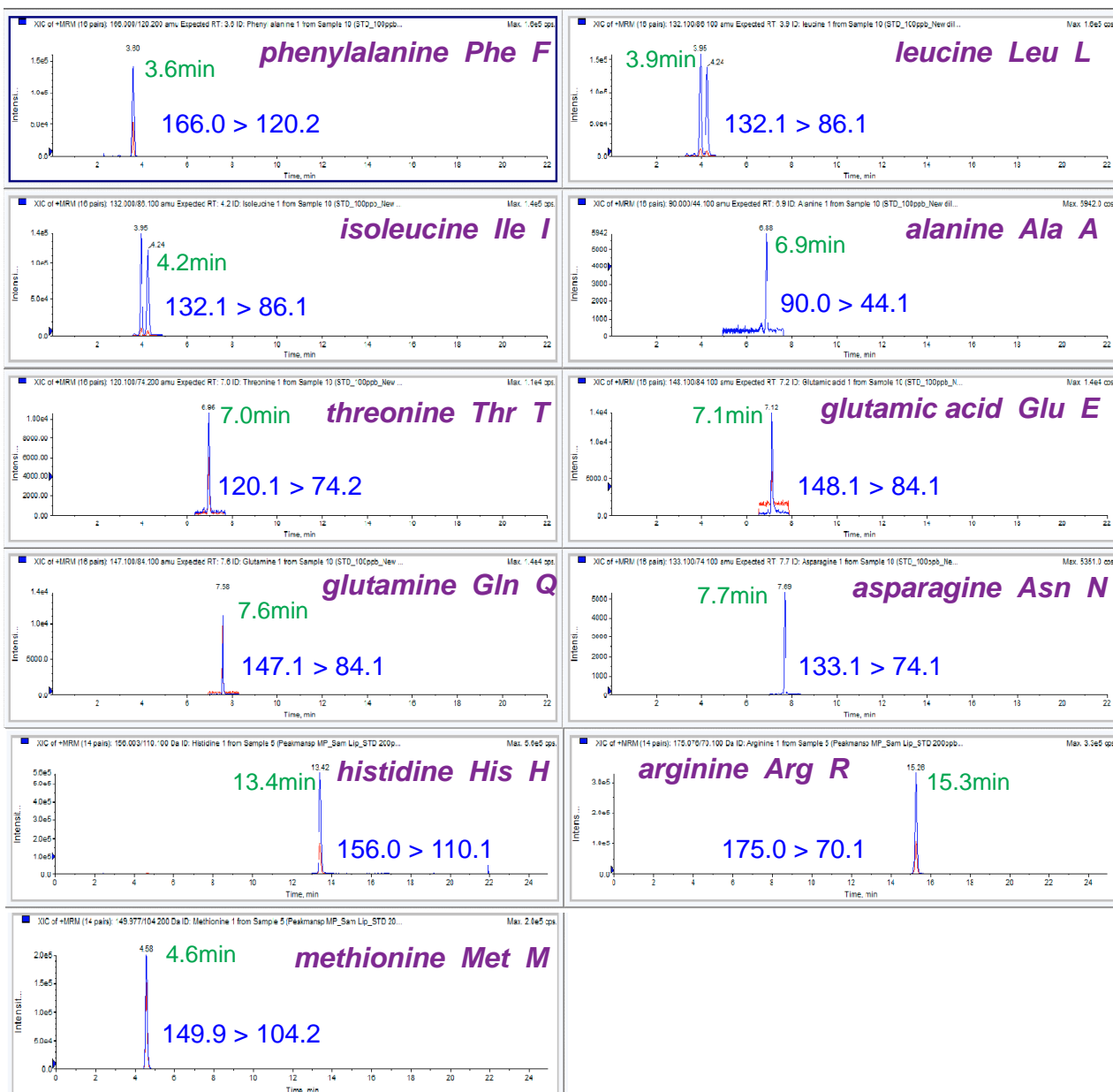

Intrada Amino Acid, 100 x 3 mm
A: 0.3% formic acid in acetonitrile
B: 100mM ammonium formate / acetonitrile = 80 / 20
20-95%B (0-4min), 95-100%B (4-14min), 100%B (14-18min),
100-20%B (18-18.1min), 20%B (18.1-25min)
0.6mL/min, 40°C, Triple Quad. MS (ESI, Positive)
5µL (100-200 ppb)

Courtesy of PEAKMAN SP, KOREA