

Unison UK-Amino

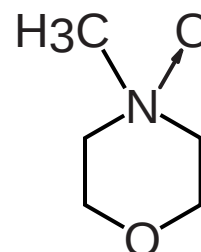
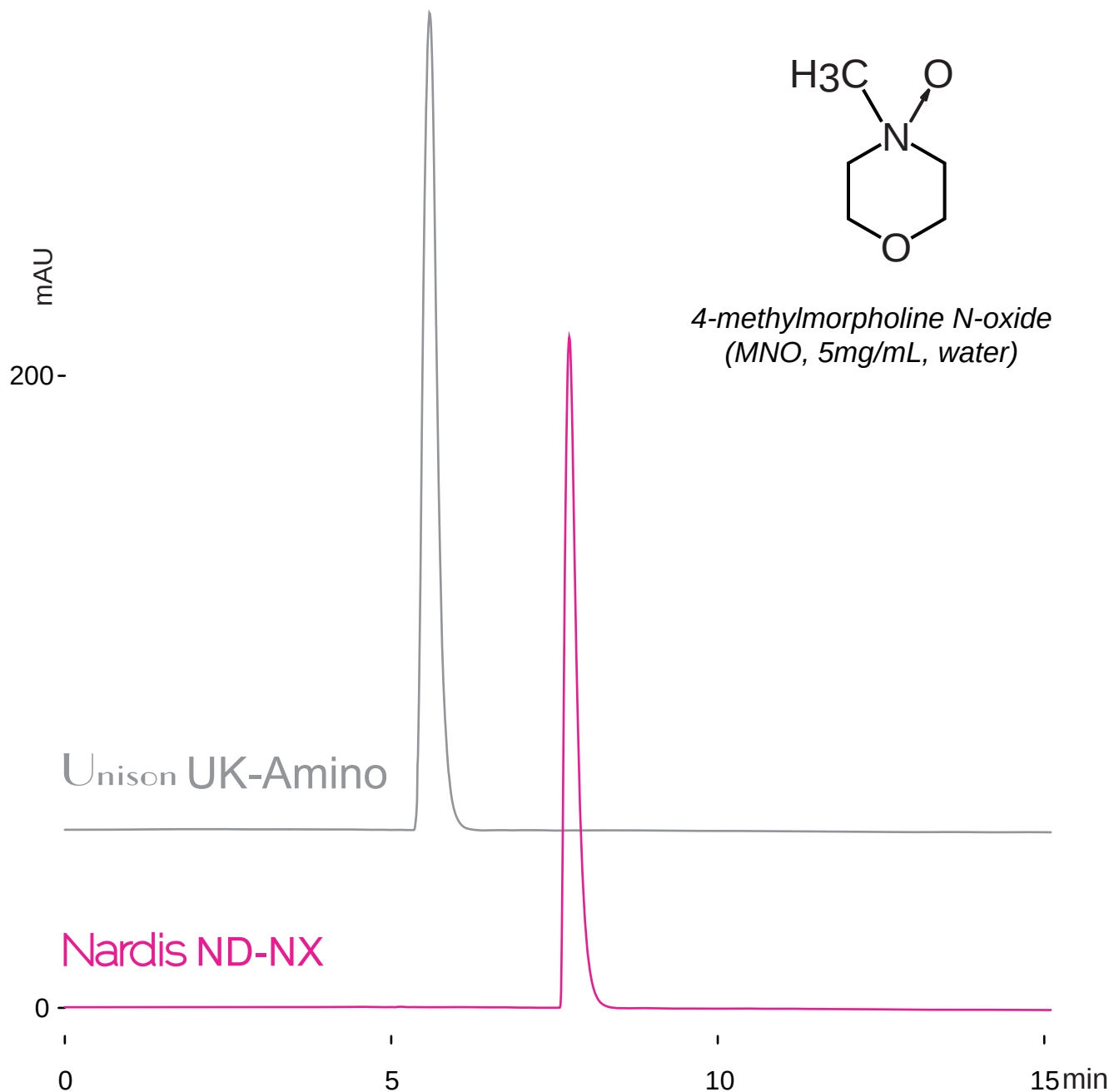
Nardis ND-NX

150 x 4.6 mm

HPLC COLUMN APPLICATION

4-Methylmorpholine N-Oxide(MNO)

4-メチルモルホリン N-オキシド (MNO) の分析


4-methylmorpholine N-oxide
(MNO, 5mg/mL, water)


150 x 4.6 mm

A: acetonitrile

B: 100mM ammonium acetate

0 -100%B (0 -15 min), 1.0 mL/min (6 MPa) , 37°C

ELSD (spray chamber 70°C, drift chamber 70°C)

10 µL (5.0mg/mL, water)